


Unison UK-Amino

150 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Guanidine hydrochloride salt

グアニジン塩酸塩の LC-MS 分析

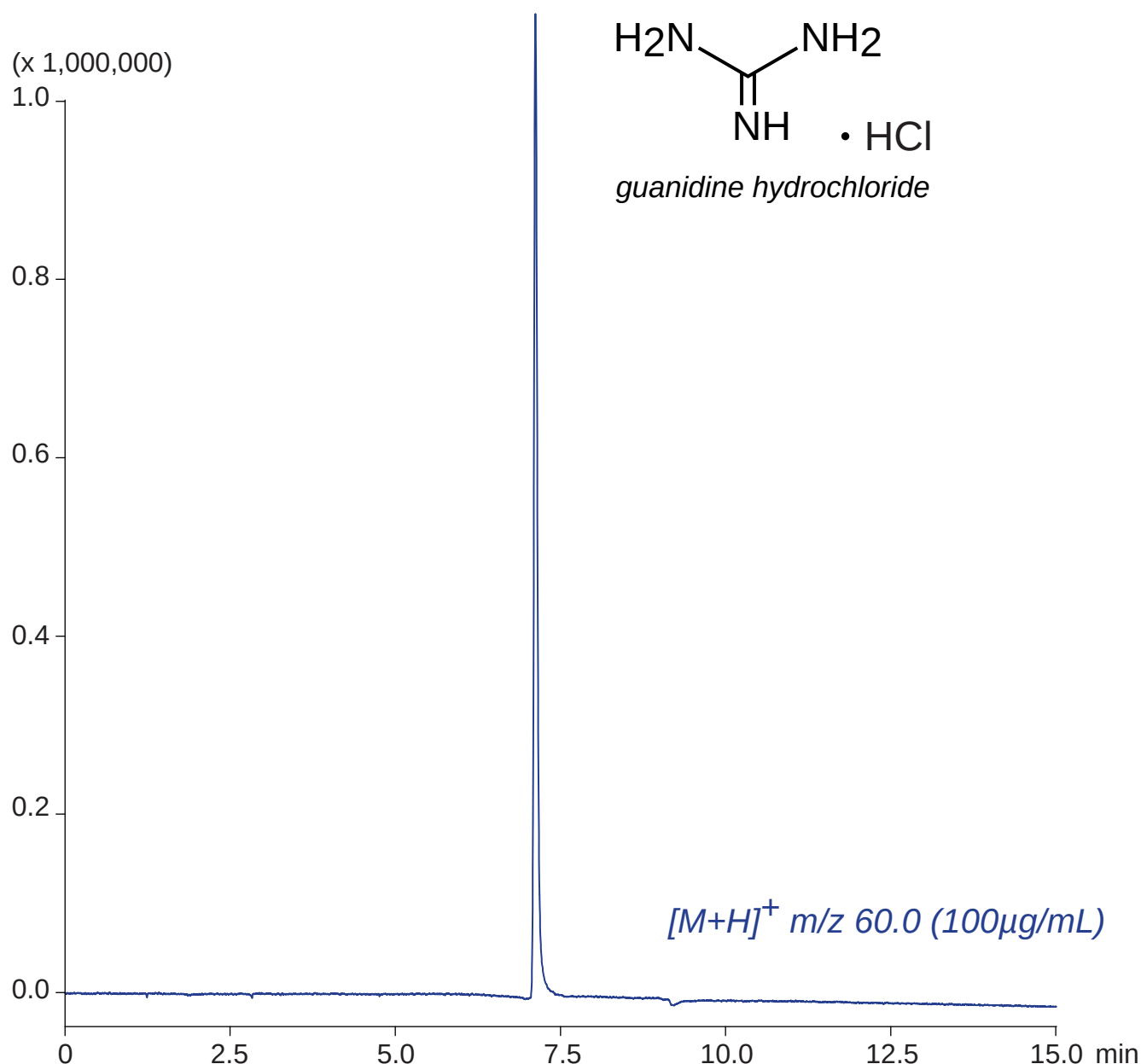

Unison UK-Amino, 150 x 3 mm

A: acetonitrile

B: acetonitrile / 100mM ammonium acetate = 10 / 90

5-100%B (0-15min)

0.4mL/min (4MPa), 37°C, Single Quad. MS (ESI)

5µL (50µg/mL, 100mM AcONH₄ / ACN = 30 / 70)