


Cadenza CD-C18

150 x 2 mm

Application

Peptides - BSA Tryptic Digest

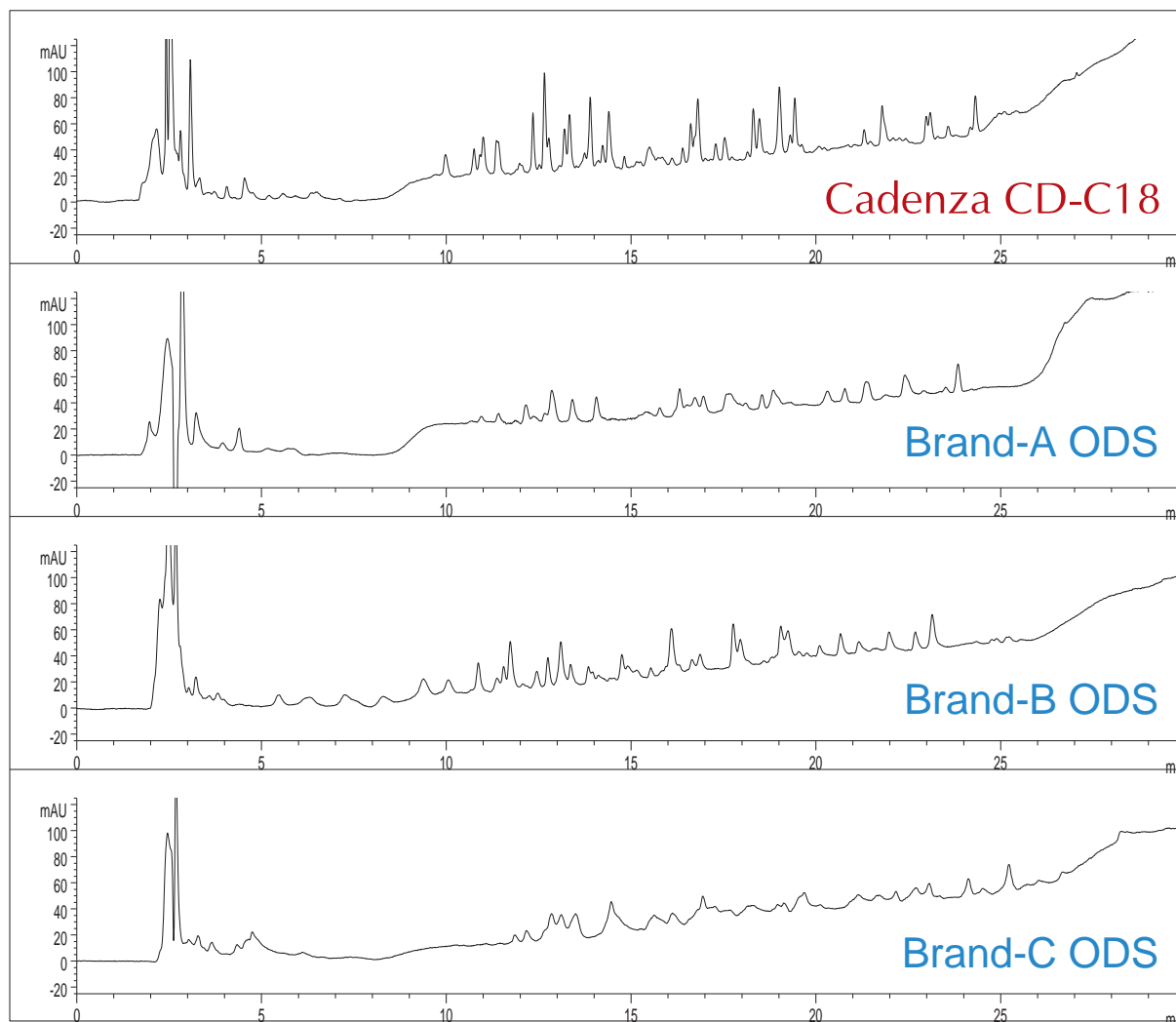

150 x 2 mm

A: acetonitrile / 1-propanol / water / acetic acid / TFA = 1 / 1 / 98 / 0.1 / 0.02

B: acetonitrile / 1-propanol / water / acetic acid / TFA = 70 / 20 / 10 / 0.09 / 0.02

5%B(0-2min), 5-40%B(2-20min), 40-90%B(20-25min), 90%B(25-30min), 90-5%B(30-33min)

0.2 mL/min

UV at 214 nm

40 °C

BSA Tryptic Digest, 18.5uL (ca.6.5ug)