


Cadenza CD-C18

150 x 2 mm

Application

LC-MS Application of Myoglobin Tryptic Digest

ミオグロビン酵素分解物のLC-MS

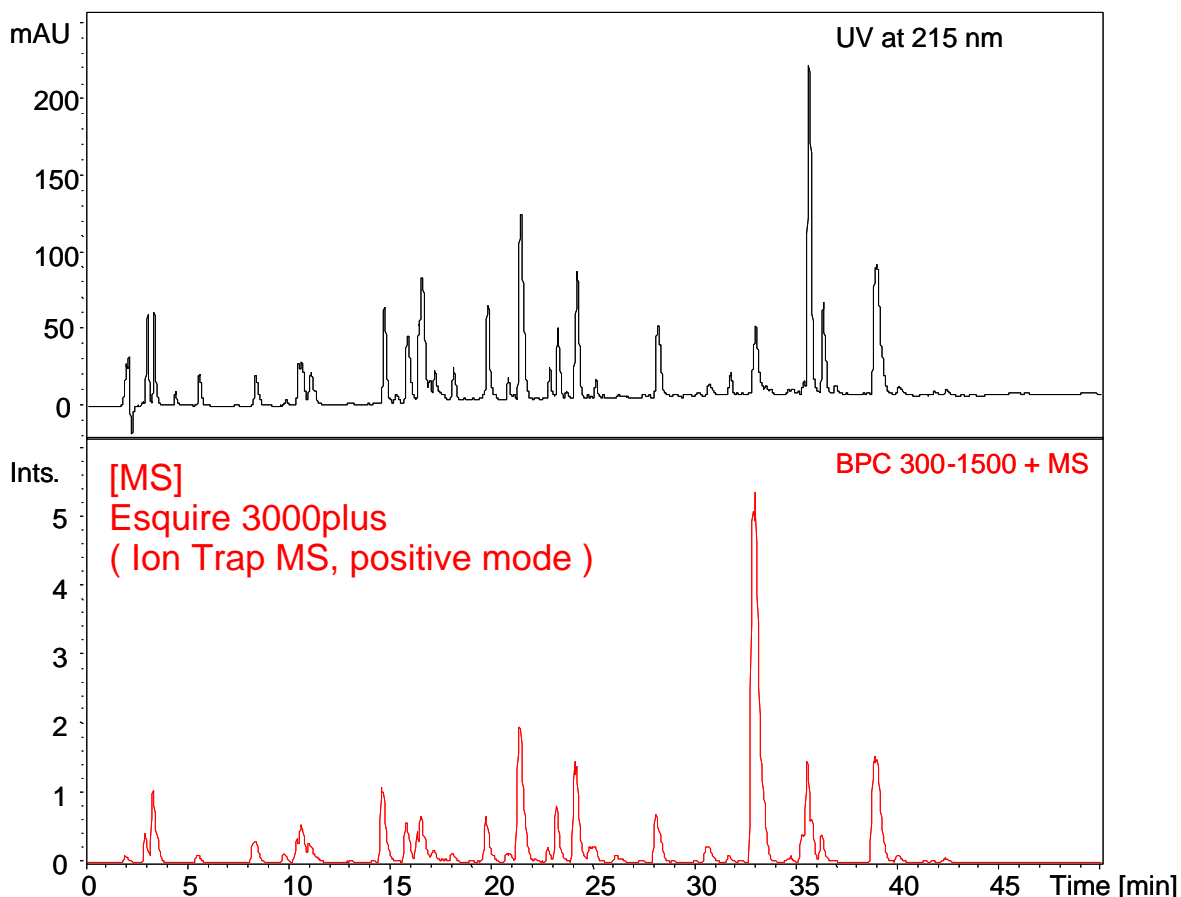

Cadenza CD-C18, 150 x 2 mm
 A: 0.01% TFA in water
 B: 0.01% TFA in acetonitrile
 2 - 50 %B (0 - 50 min)
 0.2 mL/min, 45 °C
 Myoglobin (Horse Heart), Tryptic digest

Courtesy of D.Higo, Bruker Daltonics K.K.