




Unison UK-Phenyl
 Unison UK-C8
 Cadenza CW-C18
 Cadenza HS-C18
 Cadenza CD-C18
 Cadenza CL-C18
 Unison UK-C18

150 x 4.6 mm

HPLC COLUMN APPLICATION

Minoxidil and Metabolite

ミノキシジルと代謝物

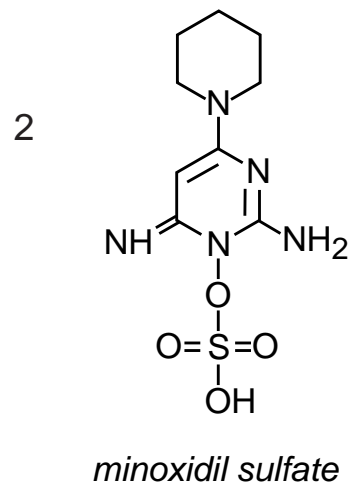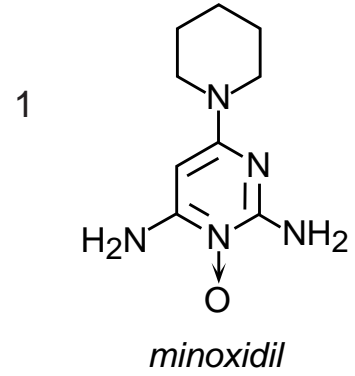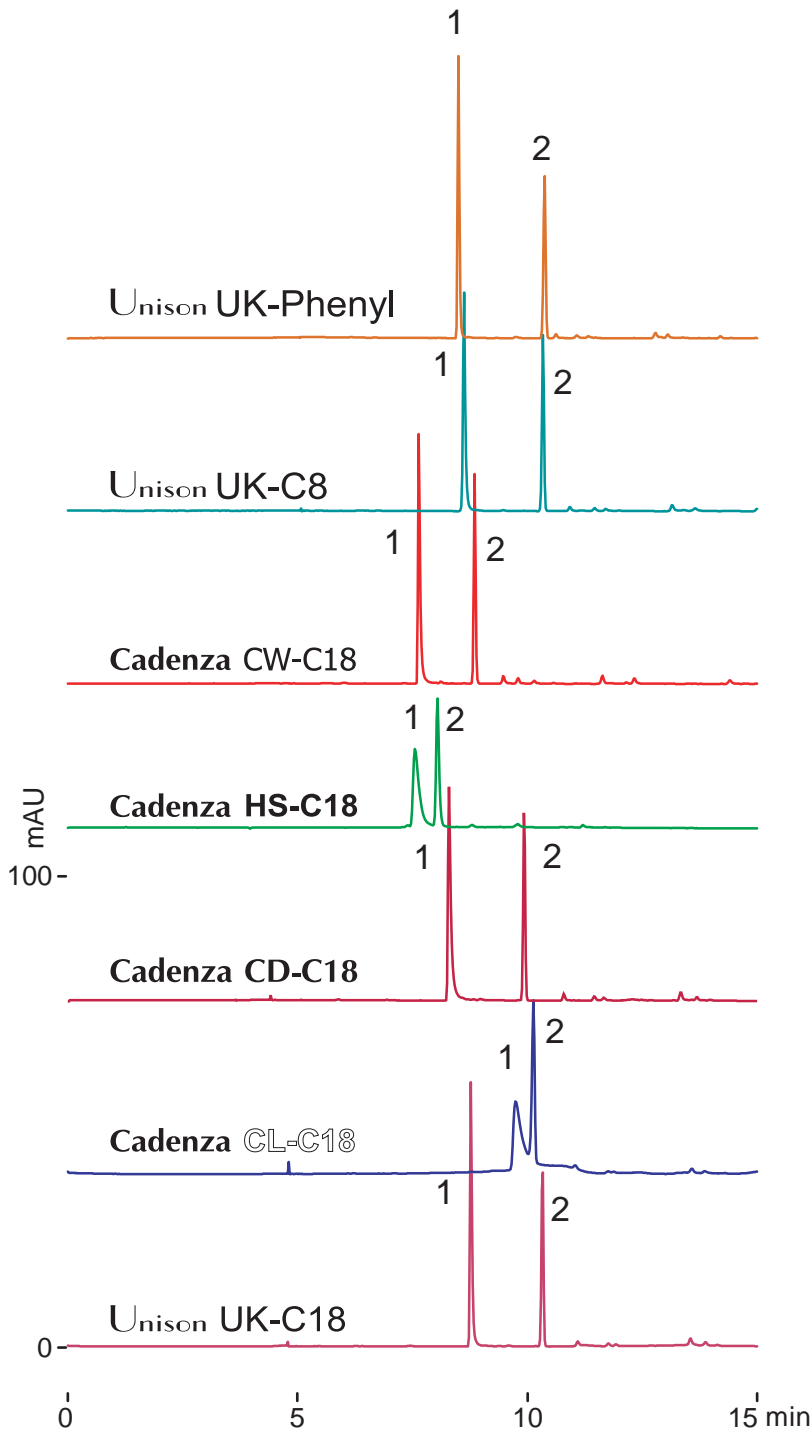

150 x 4.6 mm
 A: 50mM ammonium acetate
 B: acetonitrile
 0-70%B (0-14min)
 0.8mL/min, 37 °C
 260nm, 2µL (0.2-0.6µg)