

Cadenza CW-C18


150 x 3 mm

Application

Water-soluble vitamins


水溶性ビタミン

1

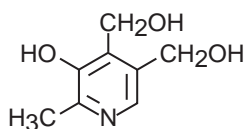


pyridoxine

3

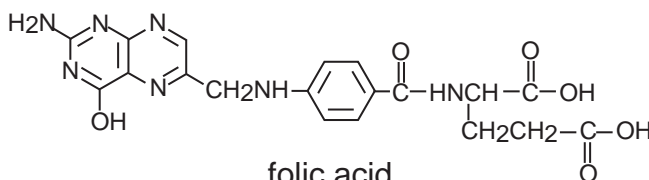


folic acid

2

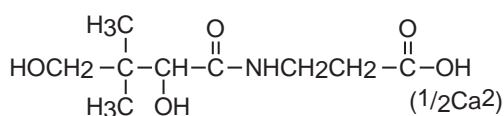


calcium pantothenate

4

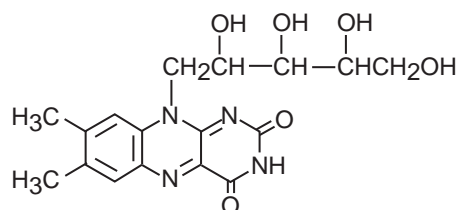

riboflavin

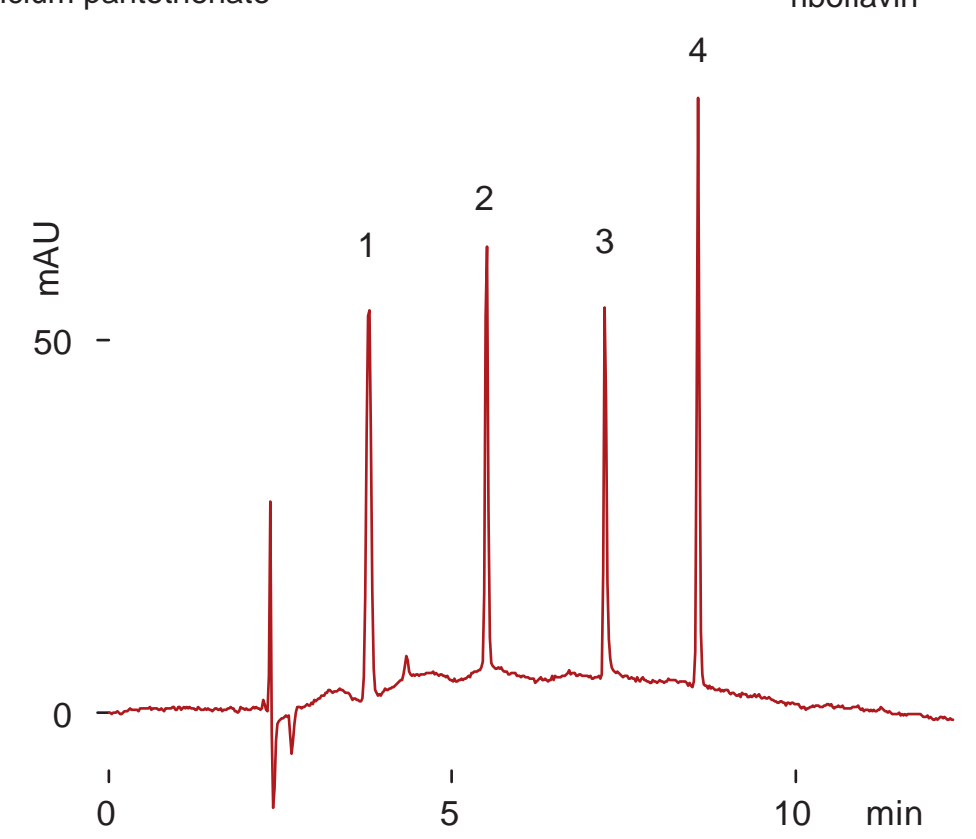

Cadenza CW-C18, 150 x 3 mm

A: water / trifluoroacetic acid = 100 / 0.1

B: acetonitrile / trifluoroacetic acid = 100 / 0.07

3-30 %B (0-10min), 0.4 mL/min (9MPa), 37 deg.C

210 nm, 1 uL (0.1-1.2ug, ethanol, water)