


Unison UK-Amino

250 x 3 mm

Application

## Saccharide analysis for beverages using Corona CAD

Corona CAD による飲料中の糖類分析

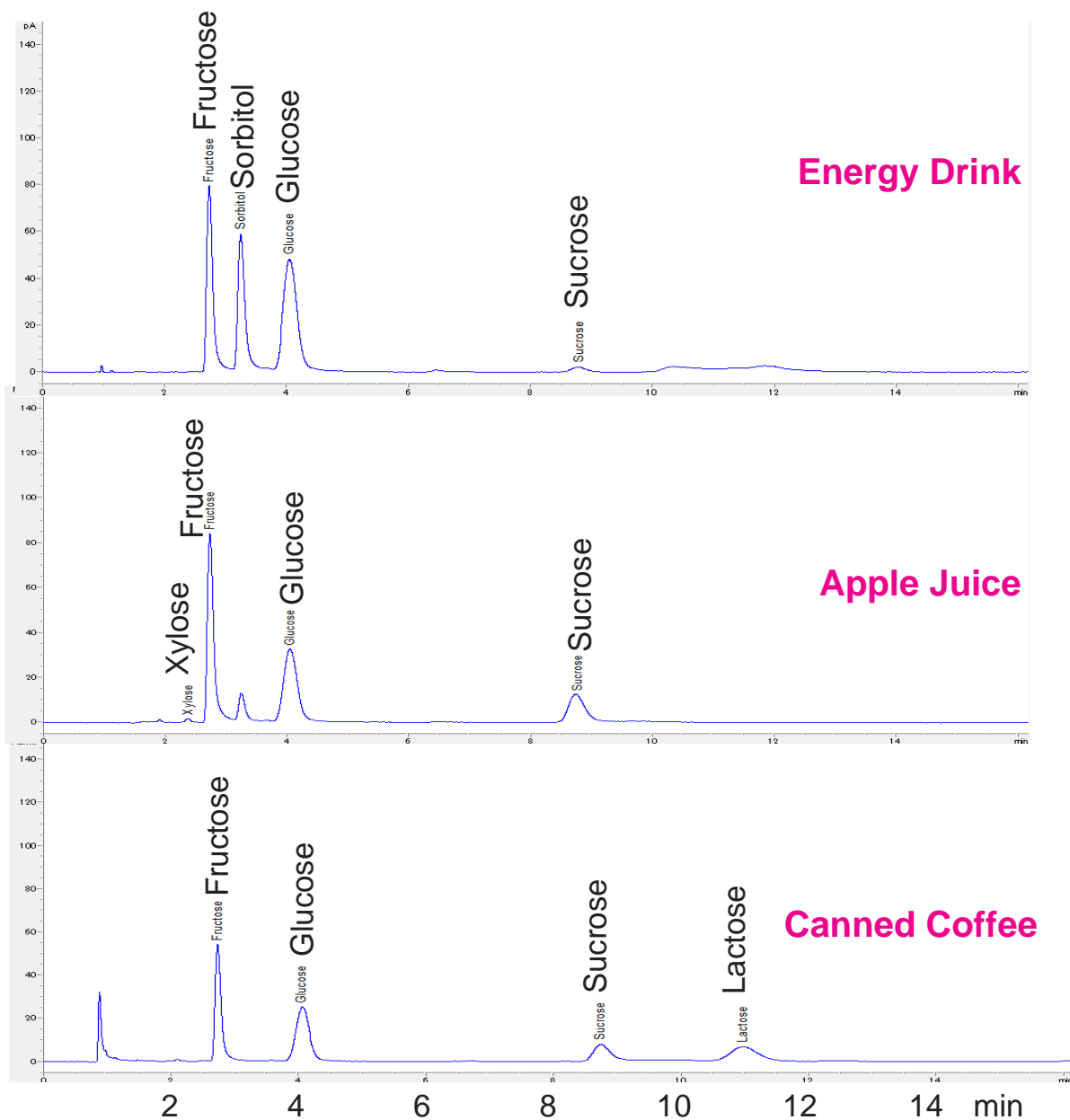

Unison UK-Amino, 250 x 3 mm  
 acetonitrile / water = 92 / 8  
 1.5 mL/min  
 Pre-column temp.: 65 °C, Post-column temp.: 30 °C  
 Corona ultra ( Filter: High, Nebulizer Temp.: OFF )  
 1.5 µL

Courtesy of Kyoko Fukushima, LMS CO., LTD. / JAPAN