

**Scherzo SM-C18**



Unison UK-C18

150 x 3 mm


Application

**Coenzymes (Nicotinamide adenine dinucleotides)**


補酵素 (ニコチンアミド アデニン ジヌクレオチド)

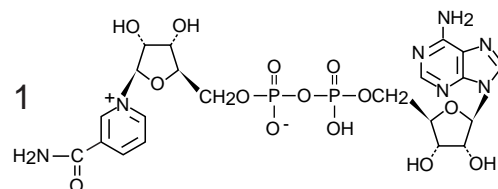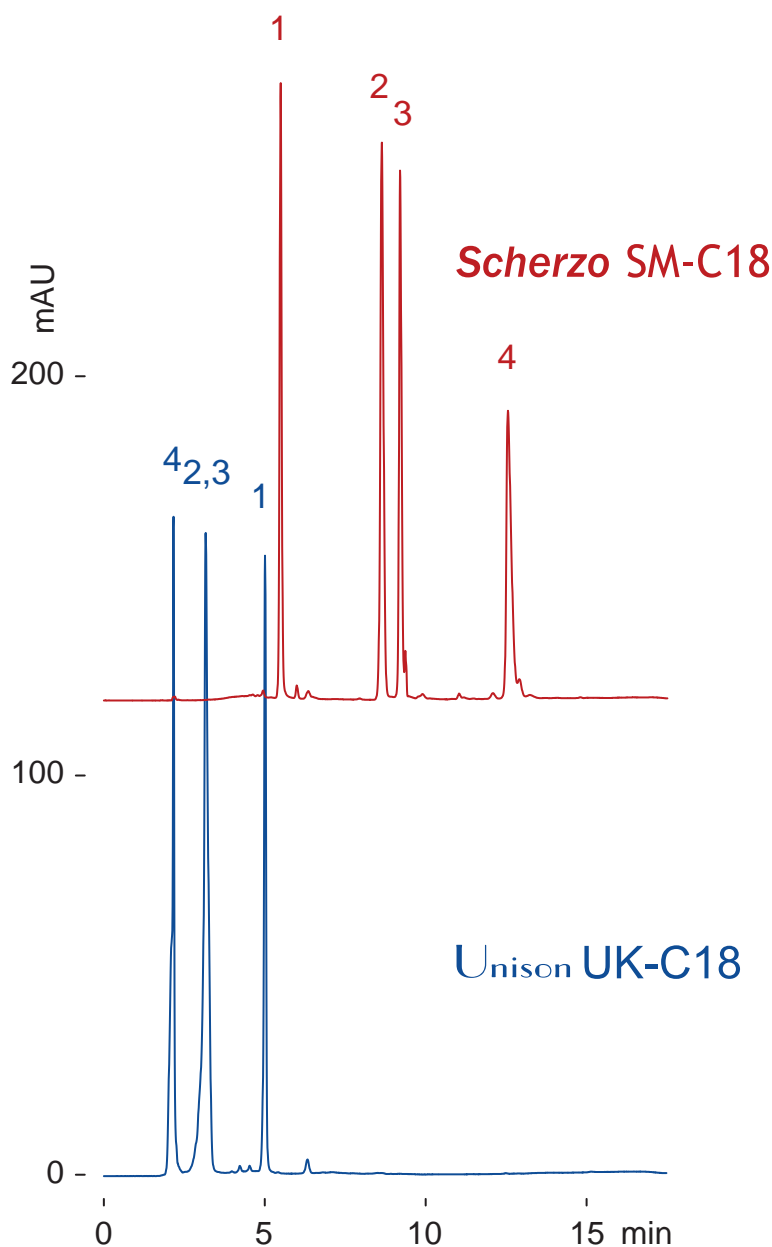


NAD

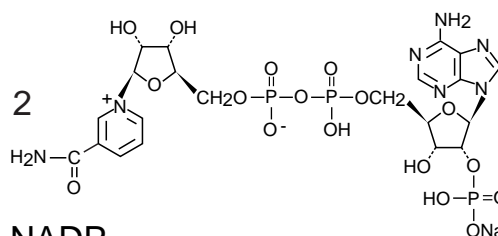

NADP

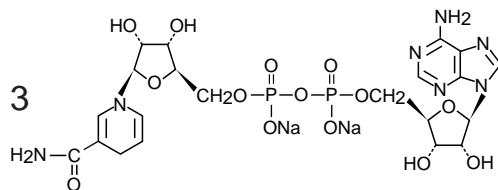

NADH

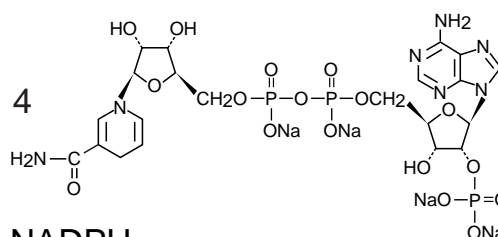

NADPH

150 x 3 mm

A: 20 mM NH<sub>4</sub>H<sub>2</sub>PO<sub>4</sub>

B: (50 mM NH<sub>4</sub>H<sub>2</sub>PO<sub>4</sub> : 50 mM (NH<sub>4</sub>)<sub>2</sub>HPO<sub>4</sub> = 1 : 1) / ACN = 70 / 30

0-100 %B (0-15 min)

0.4 mL/min (9 MPa), 37 deg.C, 260 nm, 1.6 uL (0.3 - 0.6 ug, 1 % formic acid)