



Unison UK-Phenyl

50 x 3 mm


Application

Zinc pyriothione, copper pyriothione

ピリチオン亜鉛, ピリチオン銅

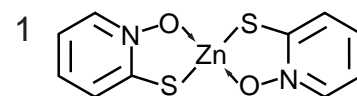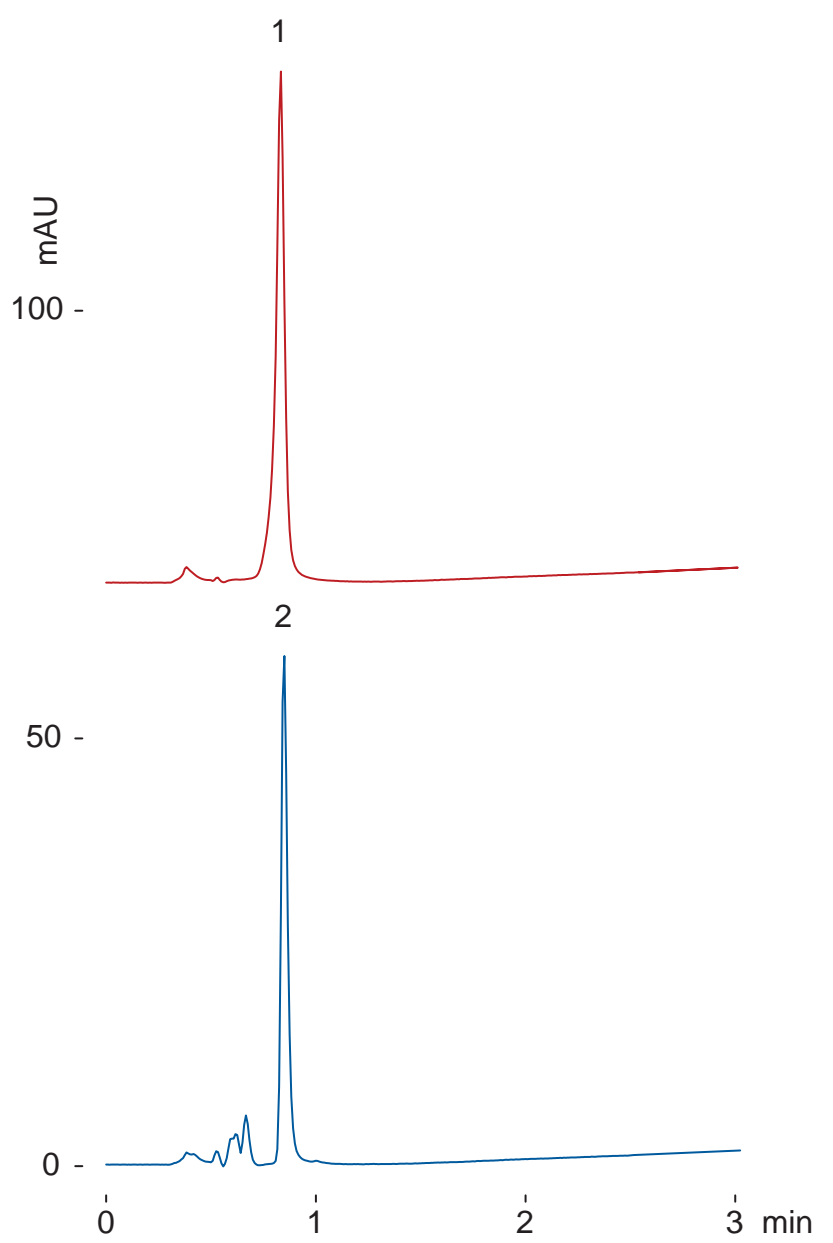

zinc pyriothione

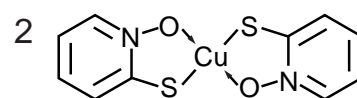

copper pyriothione

Unison UK-Phenyl, 50 x 3 mm
 A: 5mM ammonium hydrogencarbonate
 B: tetrahydrofuran
 60-95 %B (0-3 min)
 0.4 mL/min (7 MPa), 37 deg.C, 330 nm
 0.2 uL (0.1 ug, tetrahydrofuran)