


Scherzo SW-C18
 Scherzo SM-C18
 Scherzo SS-C18
 Unison UK-C18

150 x 3 mm

Application

Shikimic acid and malic acid

シキミ酸とリンゴ酸

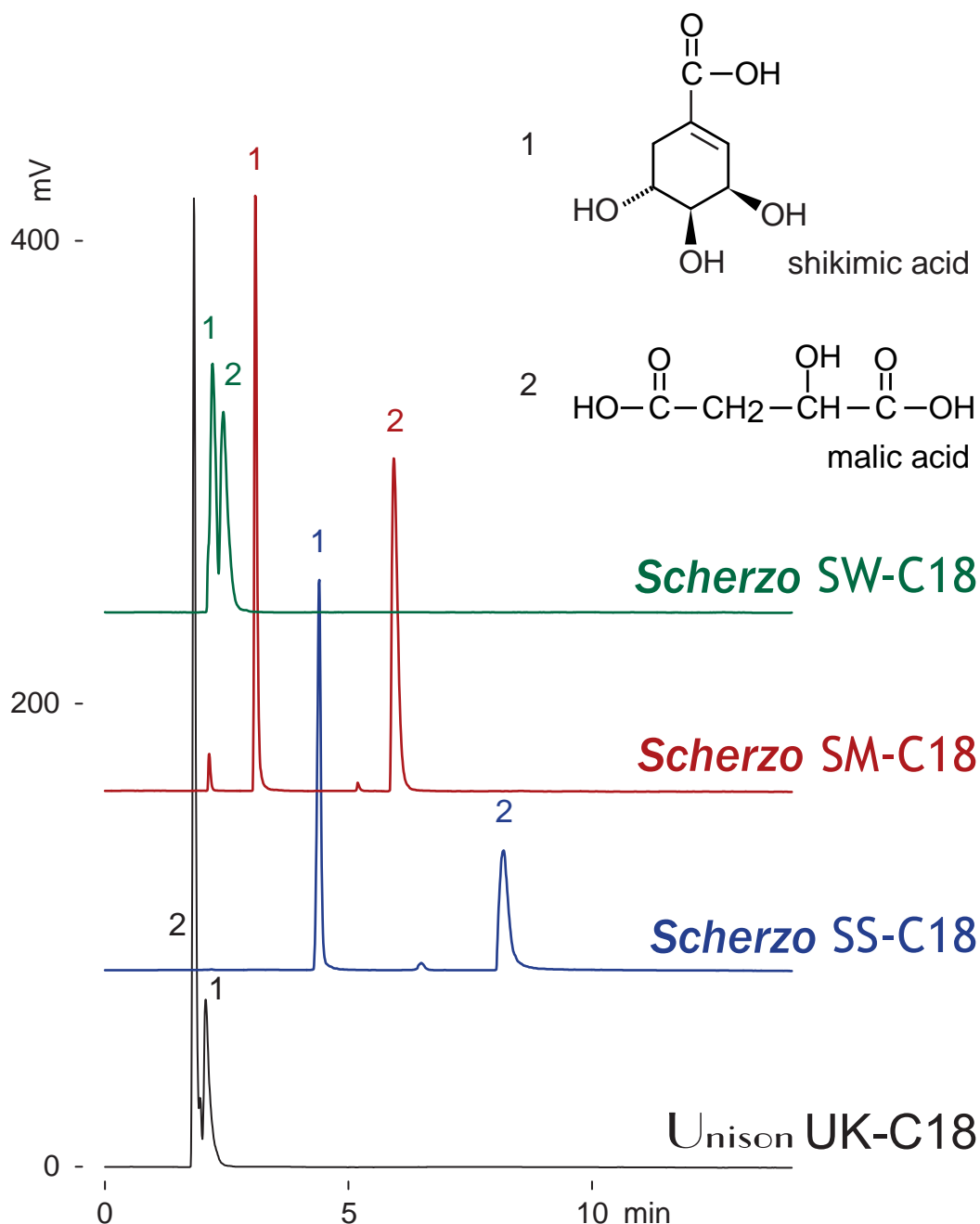

150 x 3 mm

A: 10 mM ammonium formate

B: water / acetonitrile / formic acid = 75 / 25 / 0.5

55-65 %B (0-10 min)

0.4 mL/min (8-9 MPa), 37 deg.C, ELSD(spray chamber 19 deg.C

drift tube 43 deg.C), 1 uL (0.2-0.3 ug, 50% acetonitrile)