

Intrada Amino Acid


150 x 2 mm

Application

LC-MS: 20 standard amino acids

20種アミノ酸(標準試料)のLC-MS分析

150 x 2 mm

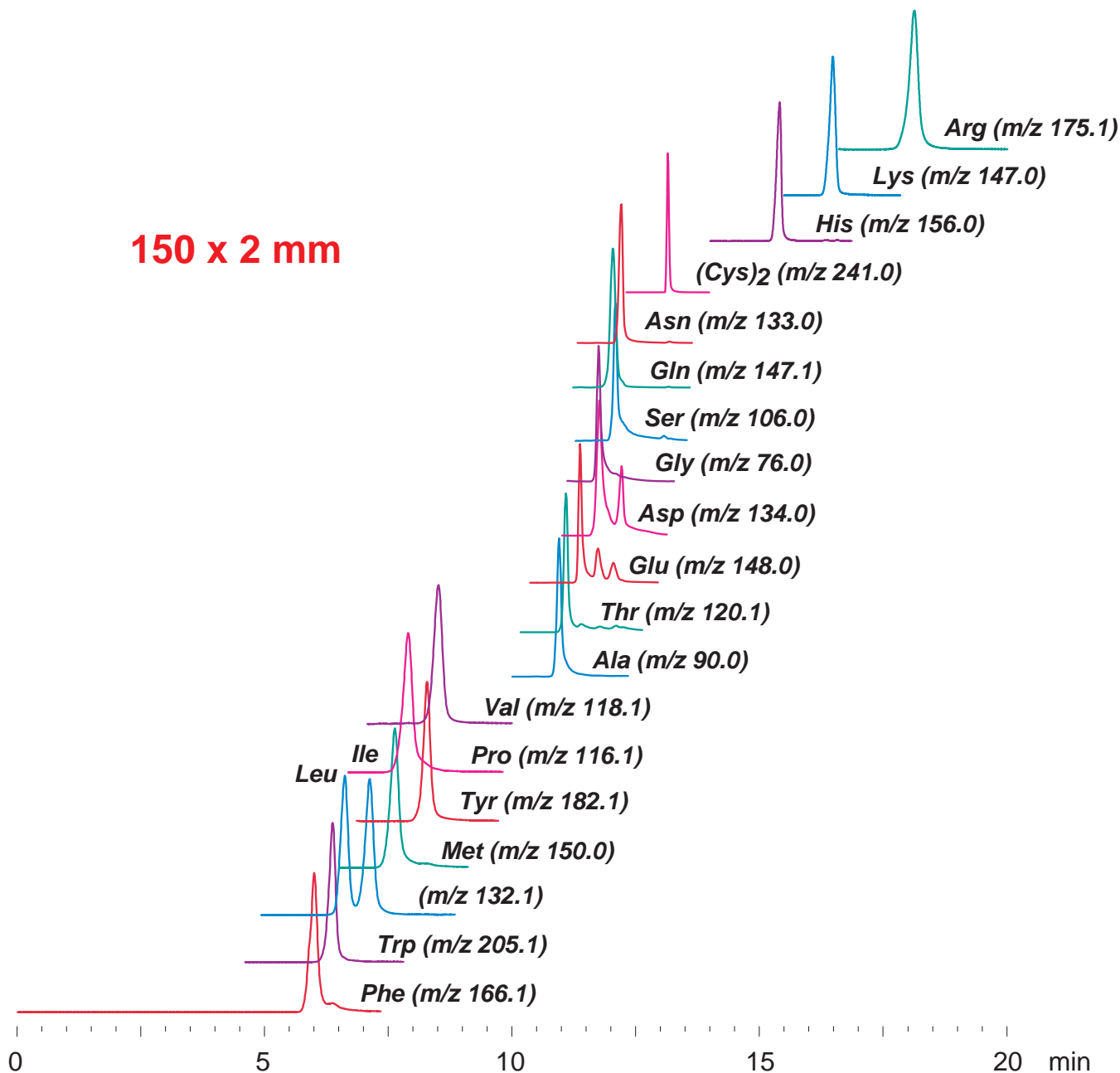

Intrada Amino Acid, 150 x 2 mm

A: acetonitrile / formic acid = 100 / 0.3

B: acetonitrile / 100mM ammonium formate = 20 / 80

18 %B (0-6 min), 18-45 %B (6-11 min), 100 %B (11-20 min)

0.3 mL/min (7MPa), 37°C, 5 uL (1000nmol/mL, 0.1N HCl)

Single Quad. MS: ESI (SIM, positive)