


Scherzo SW-C18

150 x 3 mm

Application

## LC-MS: Nicotine related compounds (2)

ニコチン関連化合物 (2)

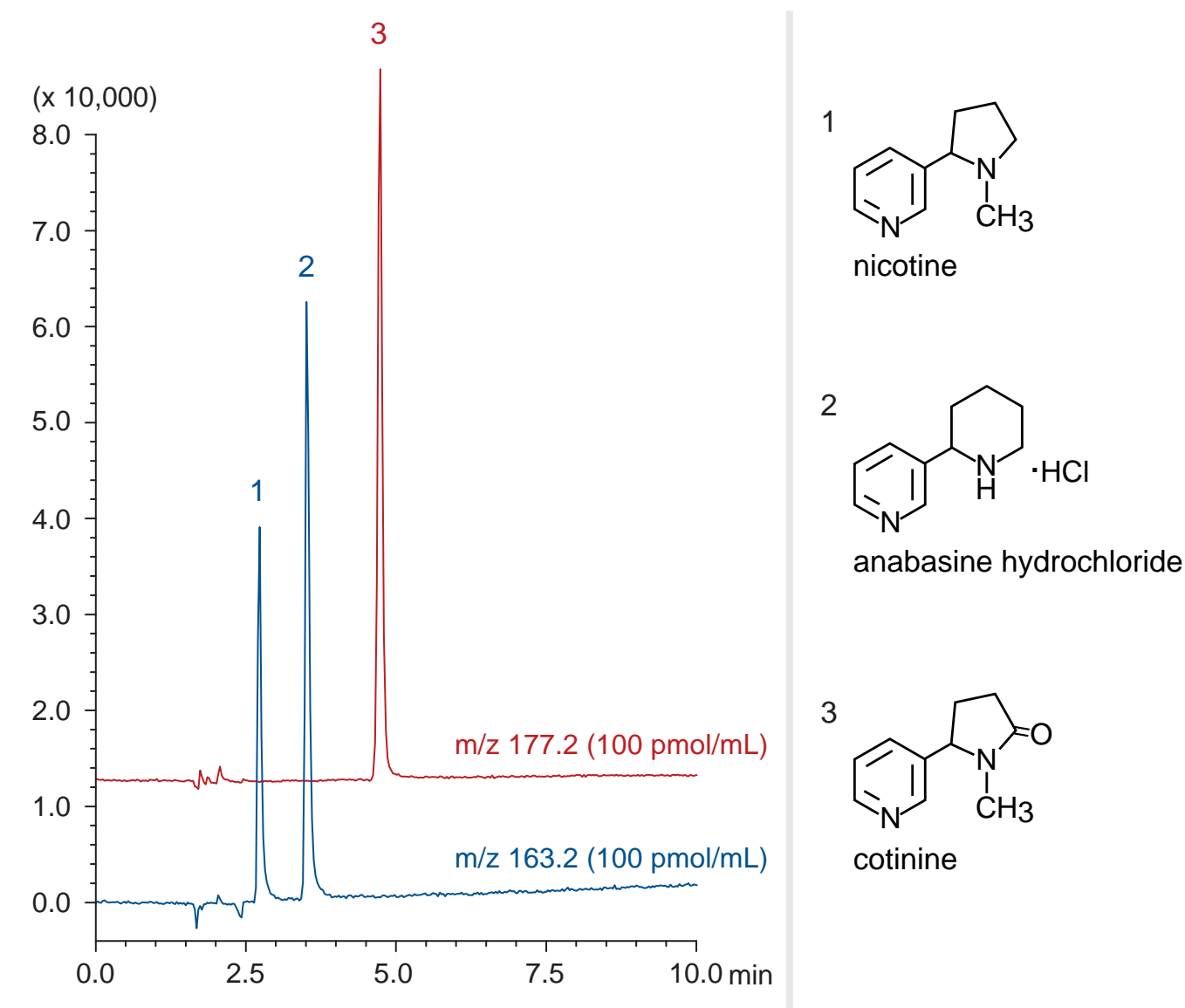

Scherzo SW-C18, 150 x 3 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / 100mM ammonium formate = 20 / 80

50-100 %B (0-10 min)

0.4 mL/min (10 MPa), 37°C, 10  $\mu$ L (100 pmol / mL, 1% HCOOH)

Single Quad. MS (ESI, positive)