


Scherzo SW-C18

150 x 3 mm

Application

LC-MS: Aspartame and diastereomer

アスパルテームとジアステレオマーのLC-MS分析

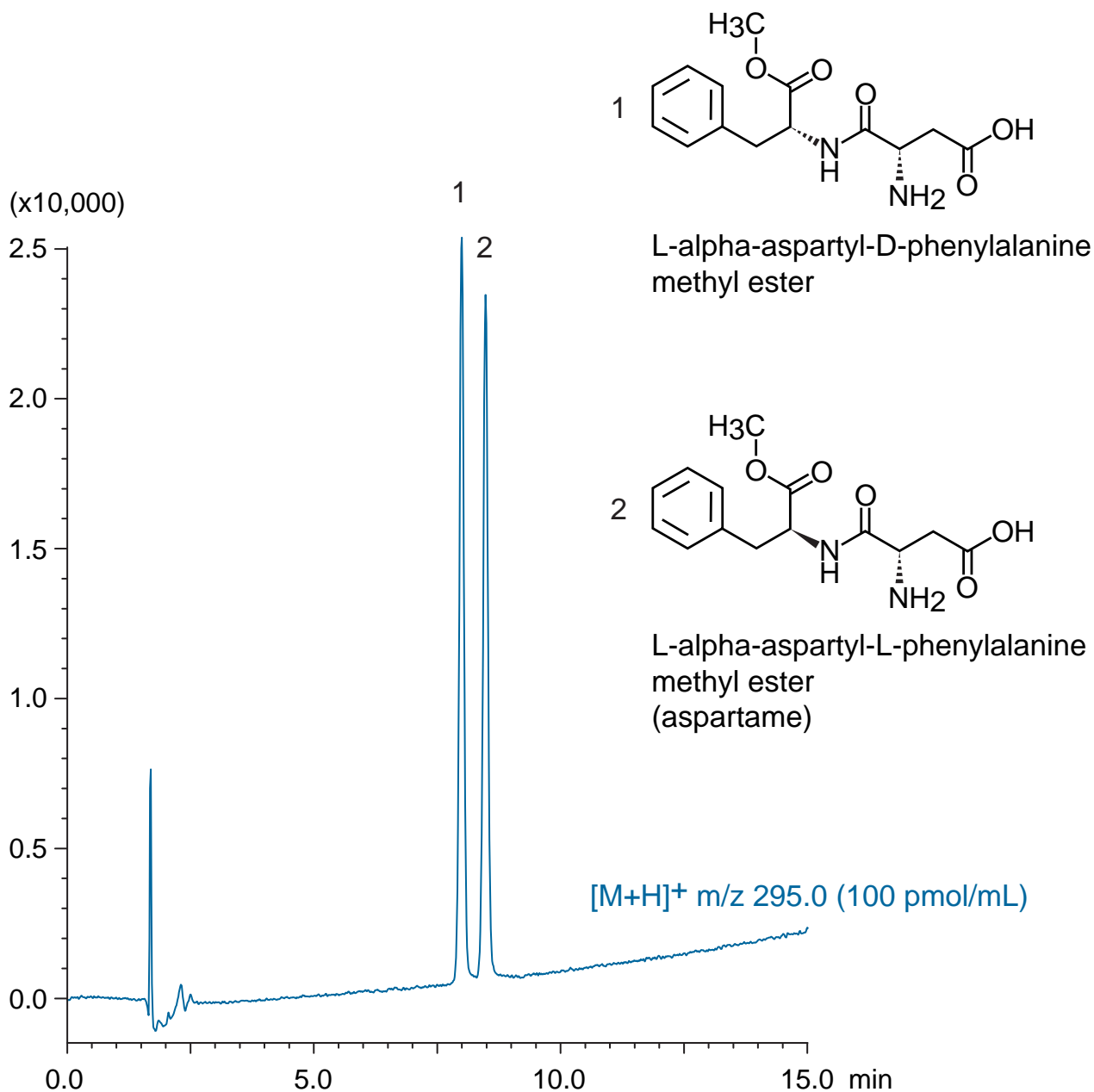

Scherzo SW-C18, 150 x 3 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / 100mM ammonium formate = 20 / 80

60-100 %B (0-15 min)

0.4 mL/min (10 MPa), 37°C, 10 uL (100pmol/mL, 1% HCOOH)

Single Quad. MS (ESI, positive)