


Intrada Amino Acid

100 x 3 mm

Application

LC-MS: choline metabolites and related compounds

コリン代謝物質と関連化合物のLC-MS分析

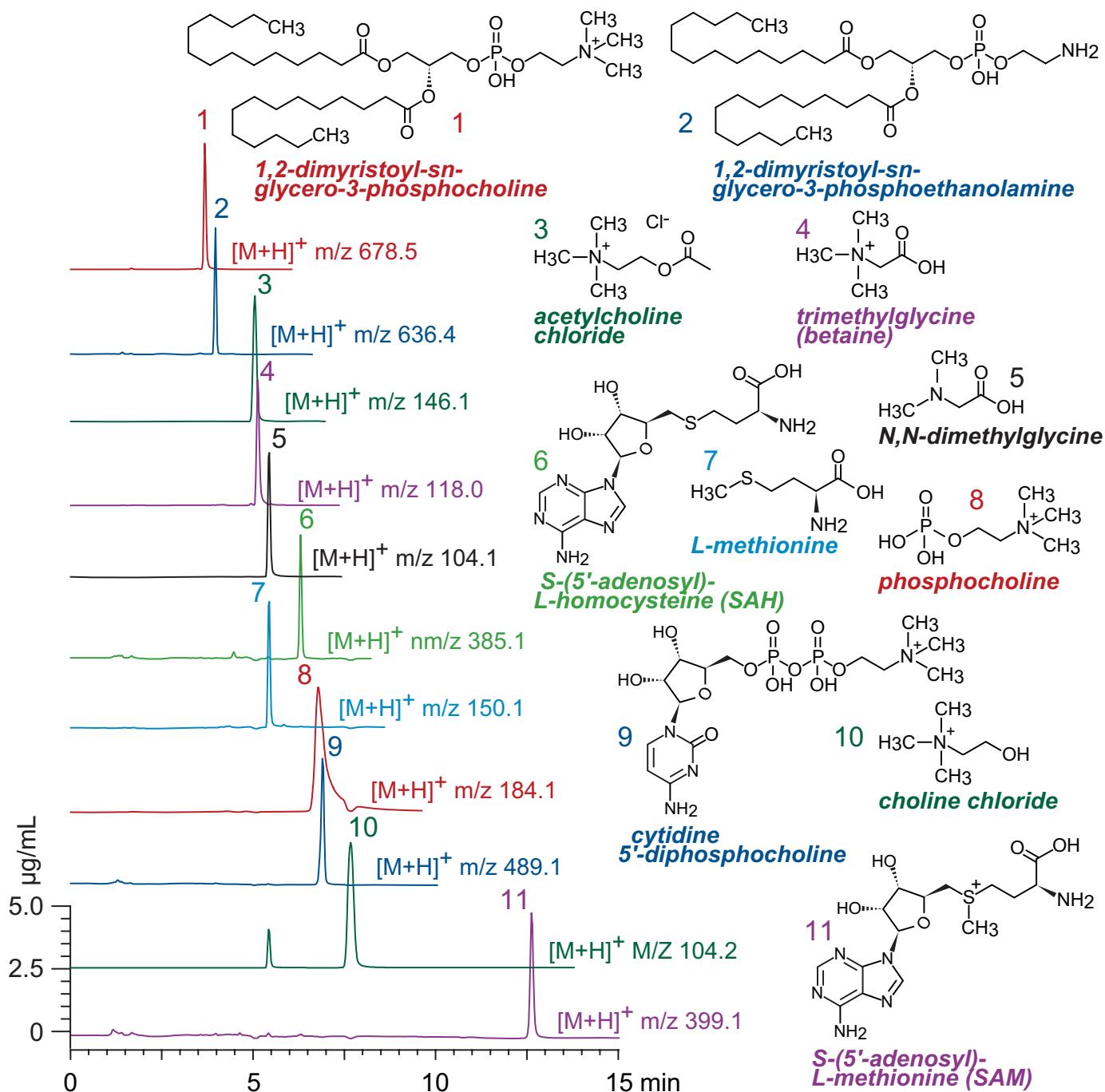

Intrada Amino Acid , 100 x 3 mm

A: acetonitrile / formic acid = 100 / 0.3

B: acetonitrile / 100mM ammonium formate = 20 / 80

10-100 %B (0-10min), 100 %B (10-15min)

0.4 mL/min (3 MPa), 37°C, 10 µL (5µg/mL, 2% HCOOH)

Single Quad. MS (ESI, positive)