



Scherzo SW-C18

50 x 3 mm


Application

LC-MS: SAH and SAM

SAHおよびSAMのLC-MS分析

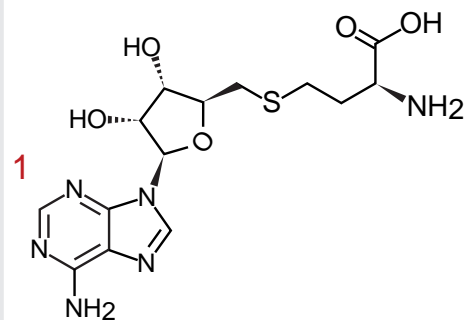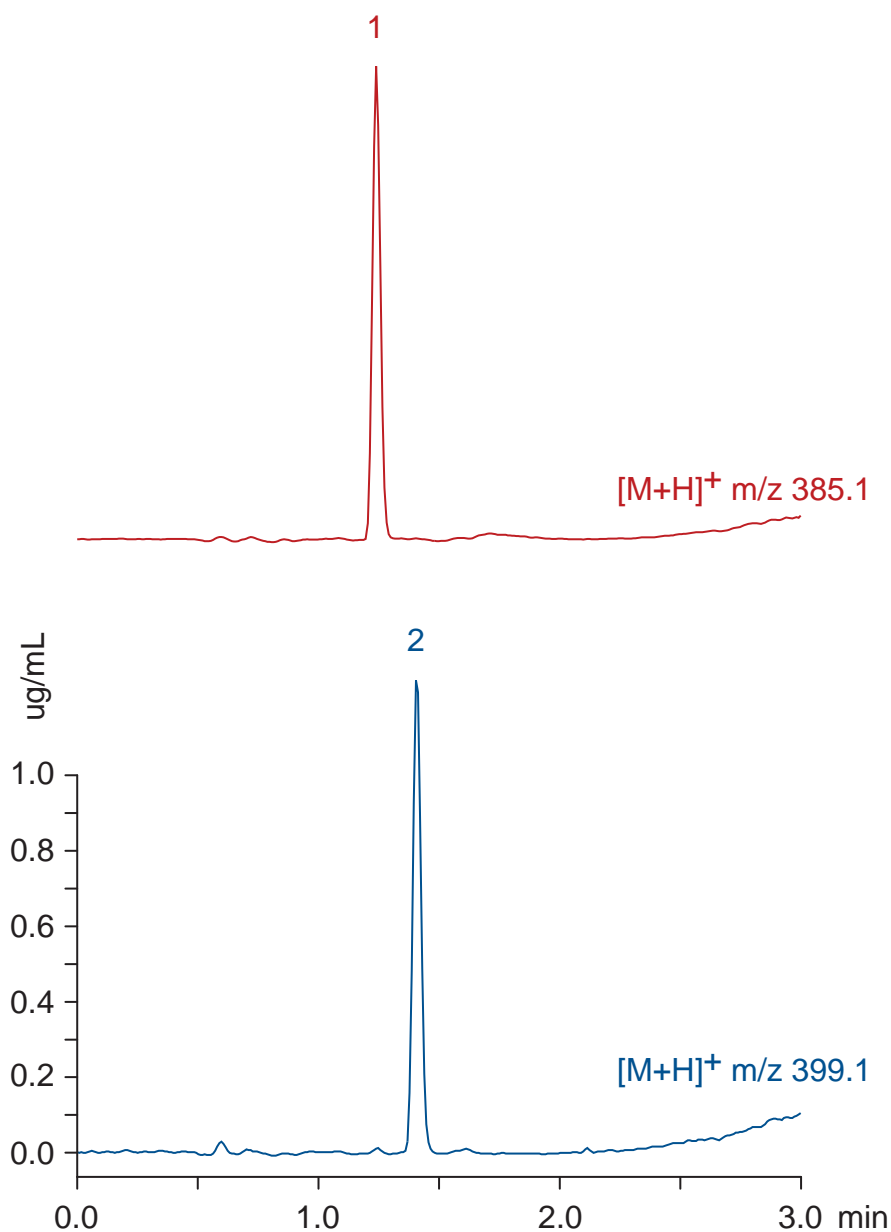

**S-(5'-adenosyl)
-L-homocysteine (SAH)**

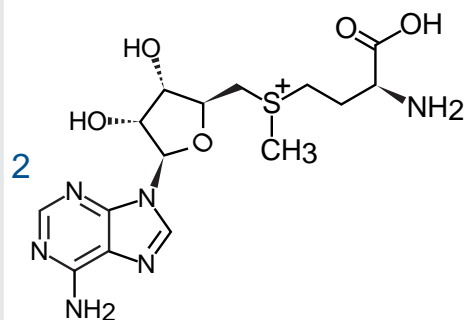

**S-(5'-adenosyl)
-L-methionine (SAM)**

Scherzo SW-C18, 50 x 3 mm

A: acetonitrile / water / formic acid = 10 / 90 / 0.1

B: acetonitrile / 50mM ammonium formate = 30 / 70

0-100 %B (0-3 min)

0.4 mL/min (4 MPa), 37°C, 10 µL (1µg/mL, 1% HCOOH)

Single Quad. MS (ESI)