







Unison UK-Amino

150 x 3 mm

Application

LC-MS: Nucleosides and bases

核酸塩基・ヌクレオシドのLC-MS分析

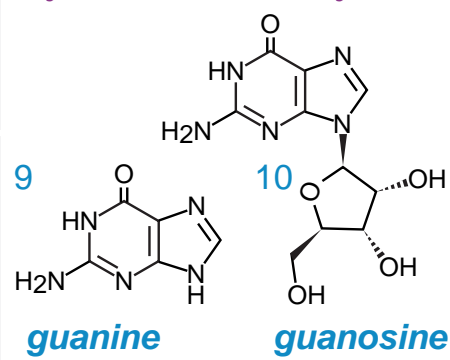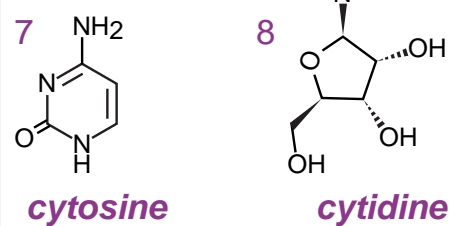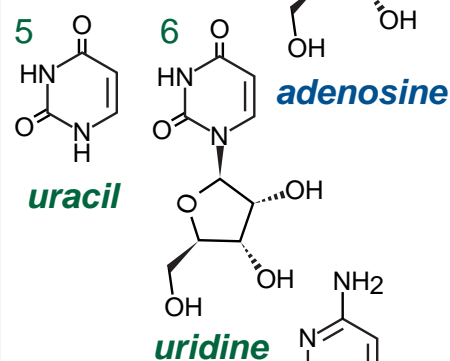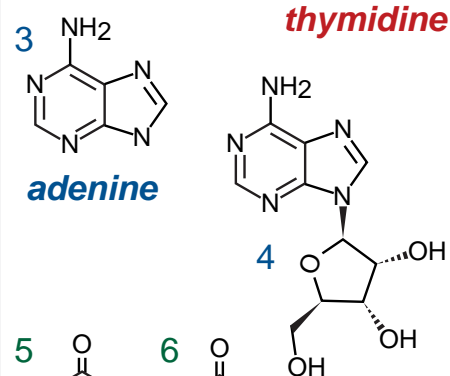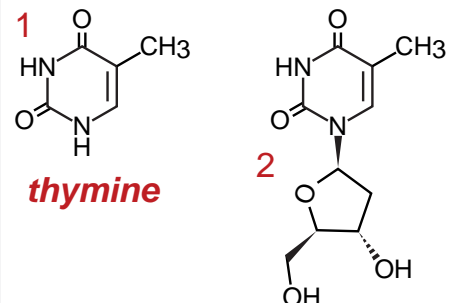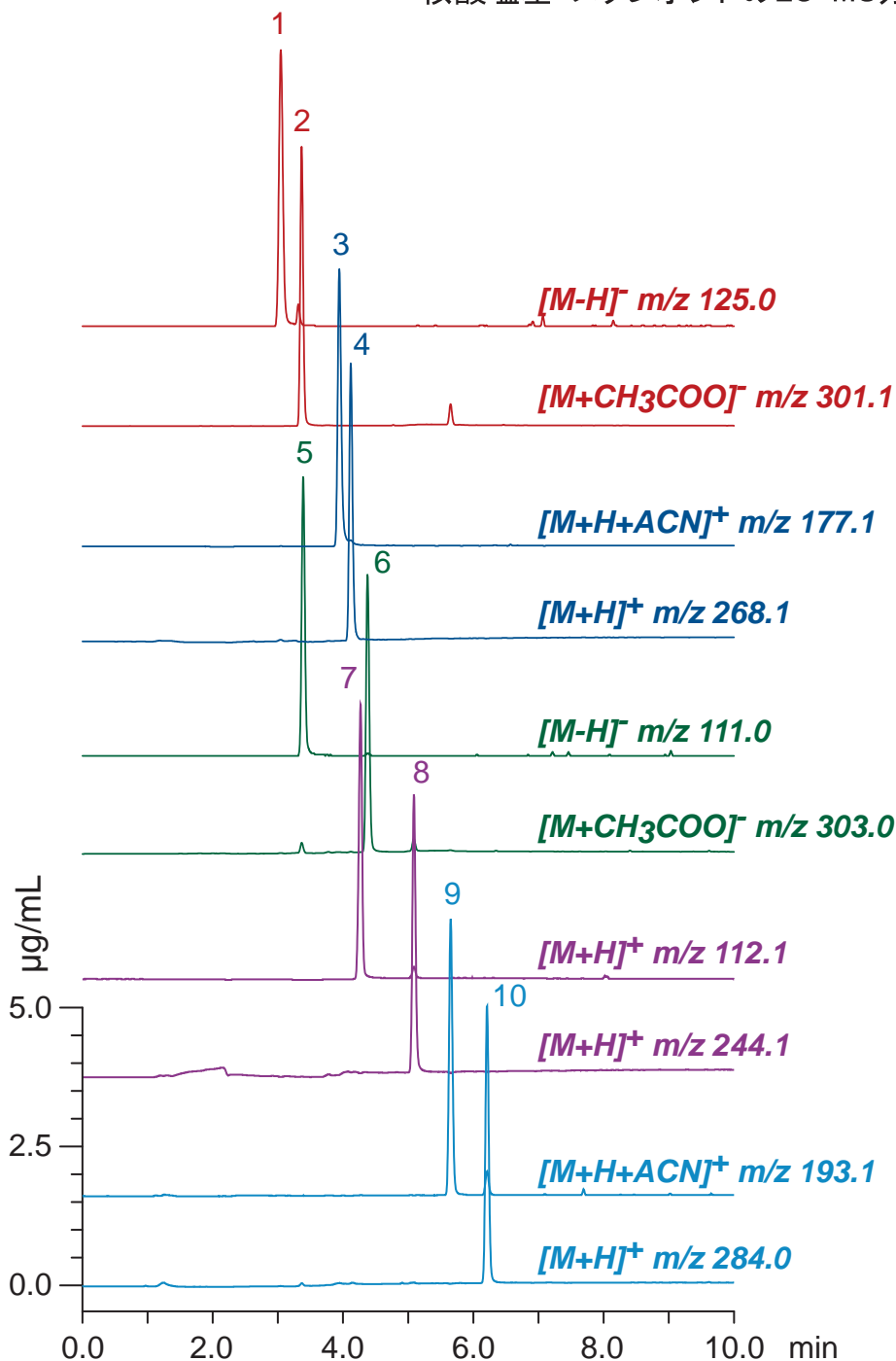

Unison UK-Amino, 150 x 3 mm  
 A: acetonitrile / water / acetic acid = 95 / 5 / 0.3  
 B: water / acetic acid = 100 / 0.3  
 0-60 %B (0-10 min)  
 0.4 mL/min (3 MPa), 60°C  
 2 μL (5μg/mL, diluted with ACN / AcOH = 100 / 2)  
 Single Quad. MS (ESI)