

Intrada SEC


250 x 4.6 mm

Application



SEC: Glycosaminoglycans (GAGs)

グリコサミノグリカン(GAGs)のSEC分離



Glycosaminoglycans (GAGs) Intrada SEC

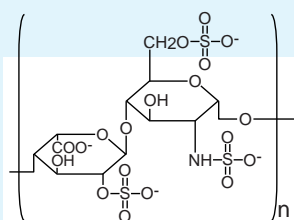


heparin

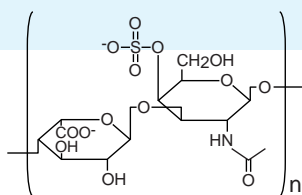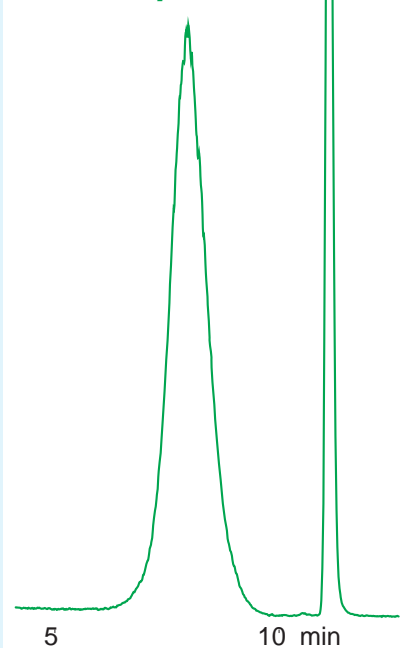

dermatan sulfate

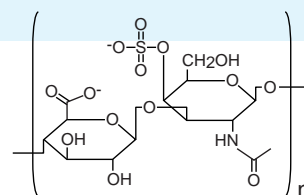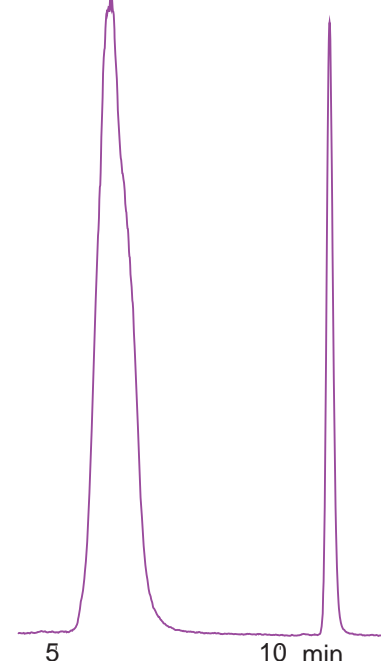

chondroitin sulfate

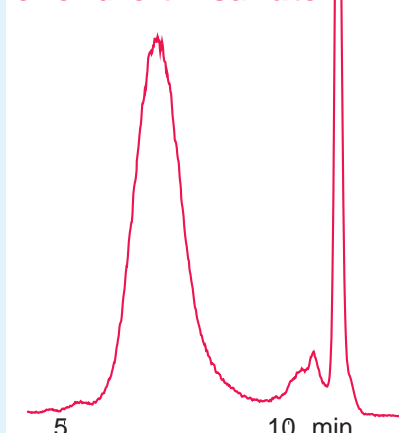

*Intrada SEC, 250x4.6mm
 100mM ammonium formate
 0.3 mL/min(5MPa), 37°C
 ELSD (spray chamber 50°C, drift tube 100°C)
 10 μL (5mg/mL)*