



Cadenza CD-C18 MF

50 x 2 mm

Application

**Metal Free Column: Hydrocortisone phosphate**

メタルフリーカラム：ヒドロコルチゾンリン酸エステル

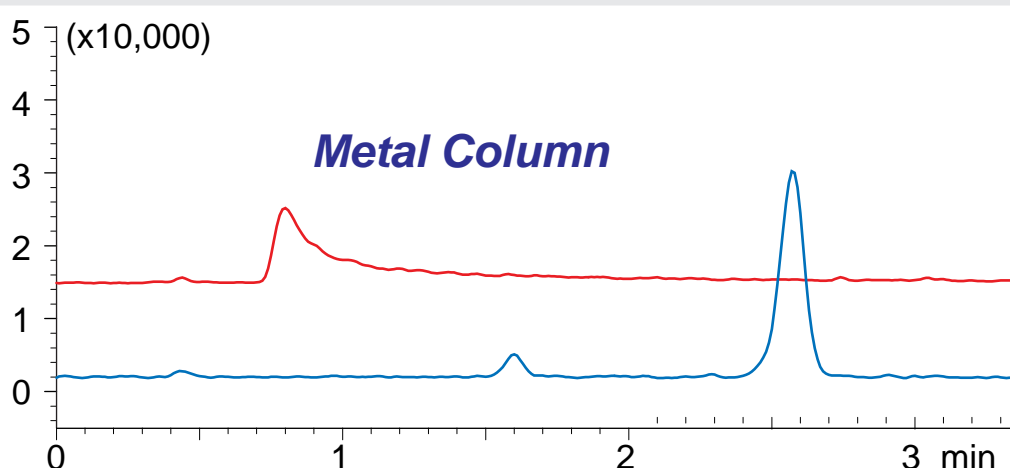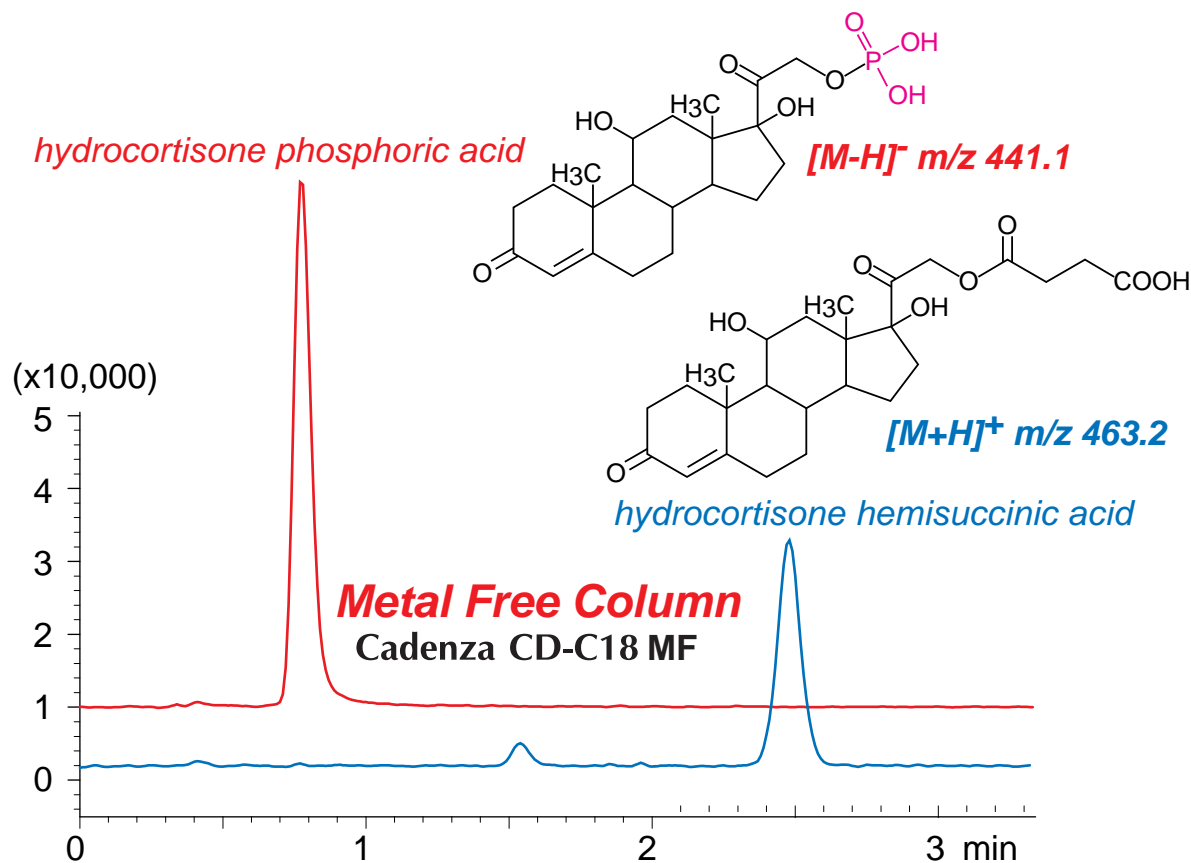

Cadenza CD-C18 MF, 50 x 2 mm  
 acetonitrile / water / formic acid = 35 / 65 / 0.5  
 0.2 mL/min (4MPa), 37°C  
 4  $\mu$ L (1 $\mu$ g/mL in 2% HCOOH)  
 Single Quad. MS (ESI)