

Cadenza CW-C18

Unison UK-C8


Unison UK-Phenyl

50 x 3 mm

Application


LC-MS: Cereulide (depsipeptide, cytotoxin)

サイトキシン, 嘔吐毒, セレウリド

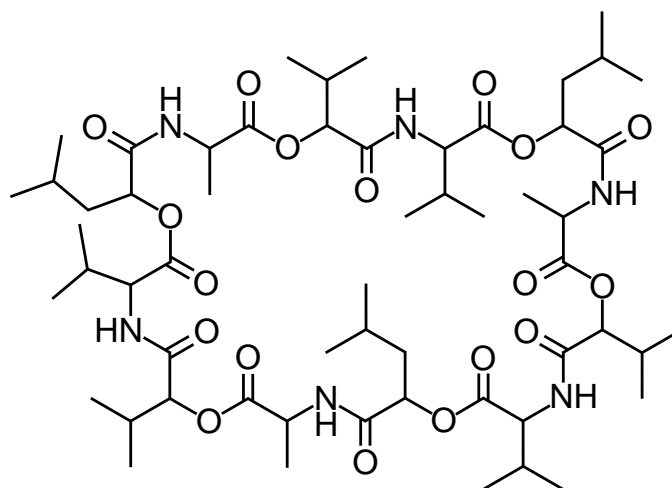

cereulide

[M+NH₄]⁺ m/z 1170.7

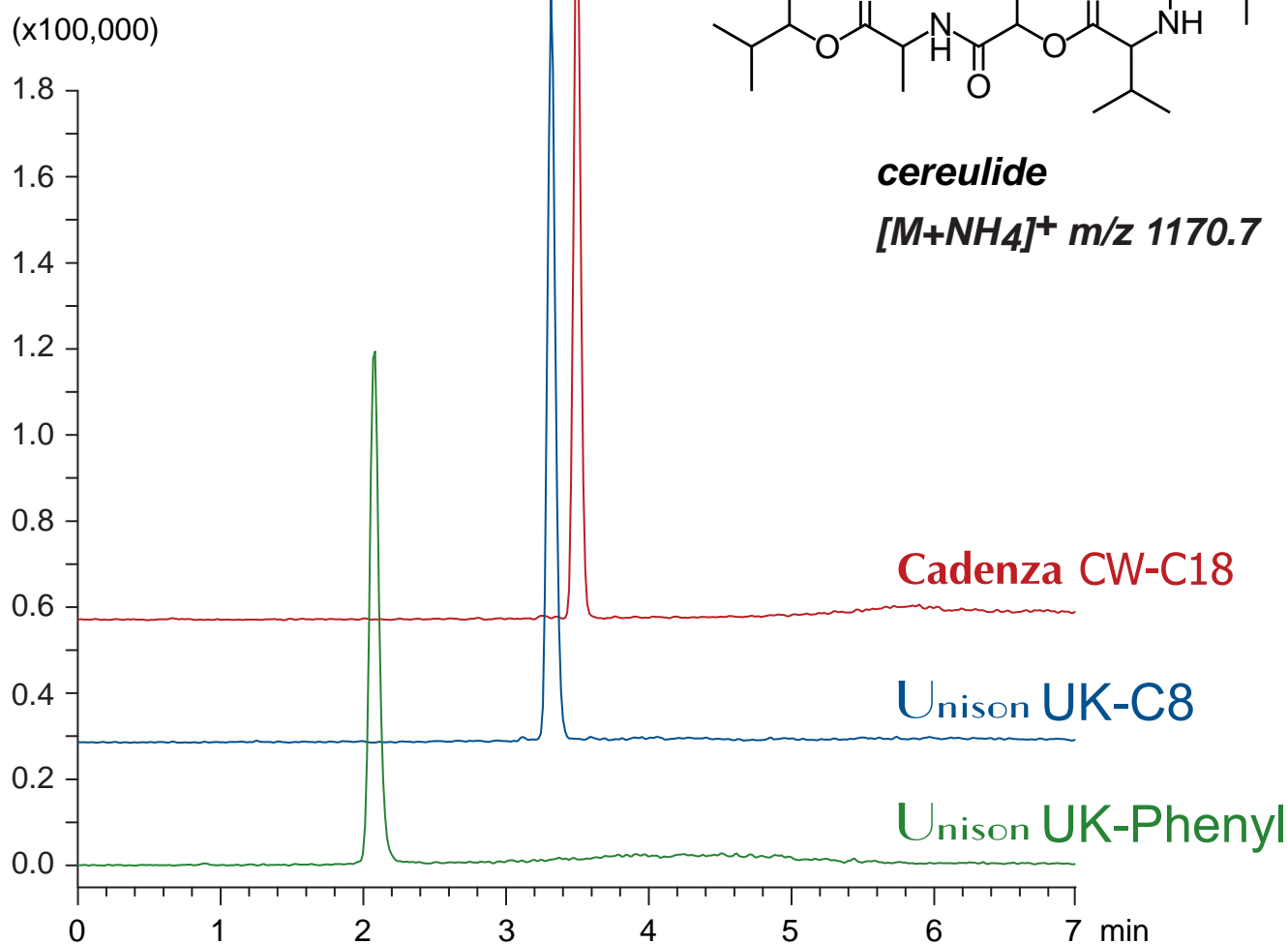

50 x 3 mm

A: 10mM ammonium acetate

B: methanol

80-100%B (0-5 min), 100%B (5-7min)

0.4 mL/min (4-5 MPa), 37°C

2 μL (500ng/mL, ACN / 10mM AcONH₄ = 90 / 10), Single Quad. MS (ESI)