


Scherzo SM-C18

100 x 3 mm

Application

Benzethonium chloride, Benzalkonium chloride


塩化ベンゼトニウム, 塩化ベンザルコニウム

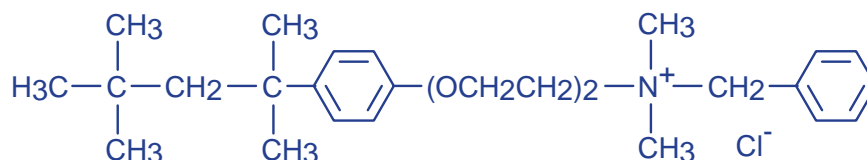

benzethonium chloride
(2mg/mL, water)
2.5 μ L

mAU

60

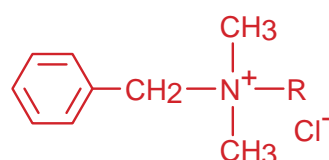

benzalkonium chloride
(2mg/mL, water)
15 μ L

40

20

0

0

5

10

15

20 min

Scherzo SM-C18, 100 x 3 mm
A: 10mM ammonium formate
B: methanol / formic acid = 100 / 0.5
25-100%B (0-15min), 100%B(15-20min)
0.4 mL/min (9MPa), 37°C, 260nm