

Unison UK-Amino MF


Unison UK-Amino

150 x 2 mm

Application

LC-MS: Ascorbic acid and related compounds

アスコルビン酸関連化合物のLC-MS分析

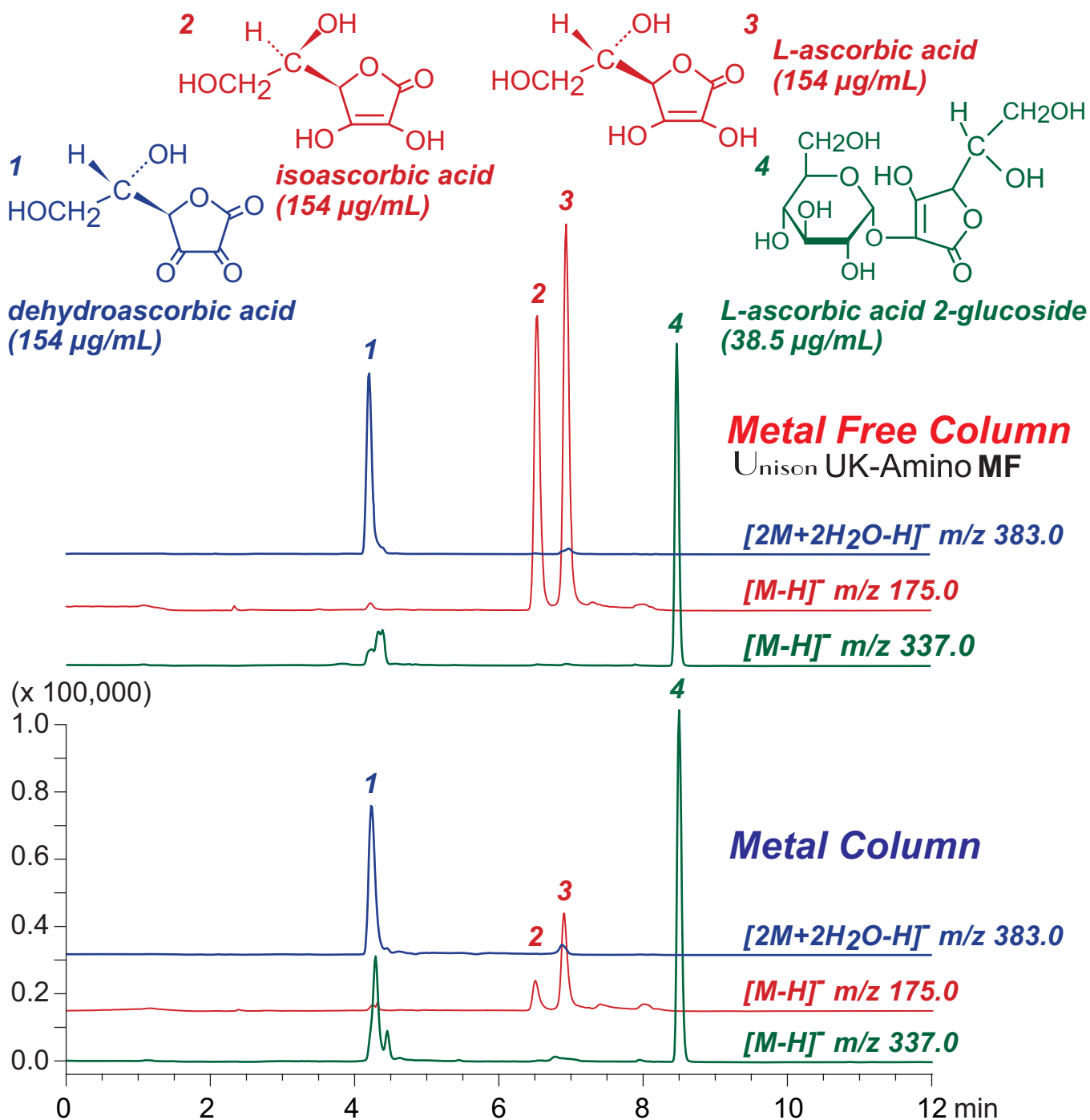

Unison UK-Amino MF, 150 x 2 mm
 A: 100mM ammonium formate
 B: acetonitrile / water / formic acid = 90 / 10 / 0.3
 100-70%B (0-8min), 0%B (8-12min)
 0.2 mL/min (4MPa), 37°C, 1 μL (38.5-154 $\mu\text{g/mL}$, water)
 Single Quad. MS (ESI)