


Intrada Organic Acid

150 x 2 mm

Application

LC-MS: Polyvalent / polycarboxylic acids (Normal-Phase Mode)

多価酸, ポリカルボン酸のLC-MS分析 (順相モード)

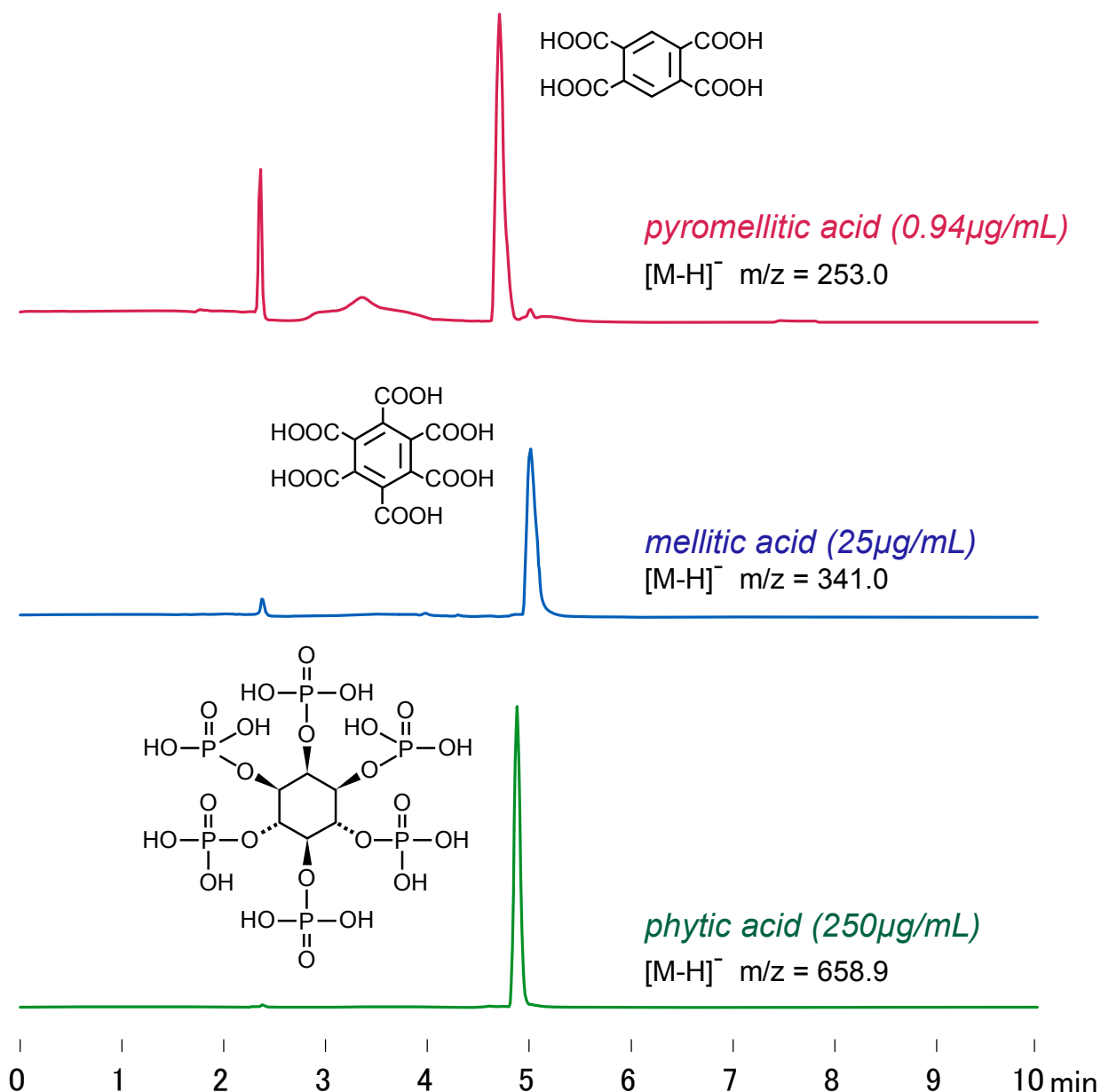

Intrada Organic Acid, 150 x 2 mm

A: acetonitrile / 100mM ammonium acetate = 95 / 5

B: acetonitrile / 100mM ammonium acetate = 10 / 90

0%B (0-1 min), 0-100%B (1-4 min), 100%B (4-10 min)

0.2 mL/min (5-6MPa), 37 °C, MS(ESI, SIM)

5 μL (50mM ammonium acetate solution)