


Cadenza CD-C18 Cadenza CW-C18
 Unison UK-C18 Cadenza CX-C18

50 x 2 mm

Application

LC-MS: Omega-3 and Omega-6 fatty acids

オメガ-3およびオメガ-6脂肪酸のLC-MS分析

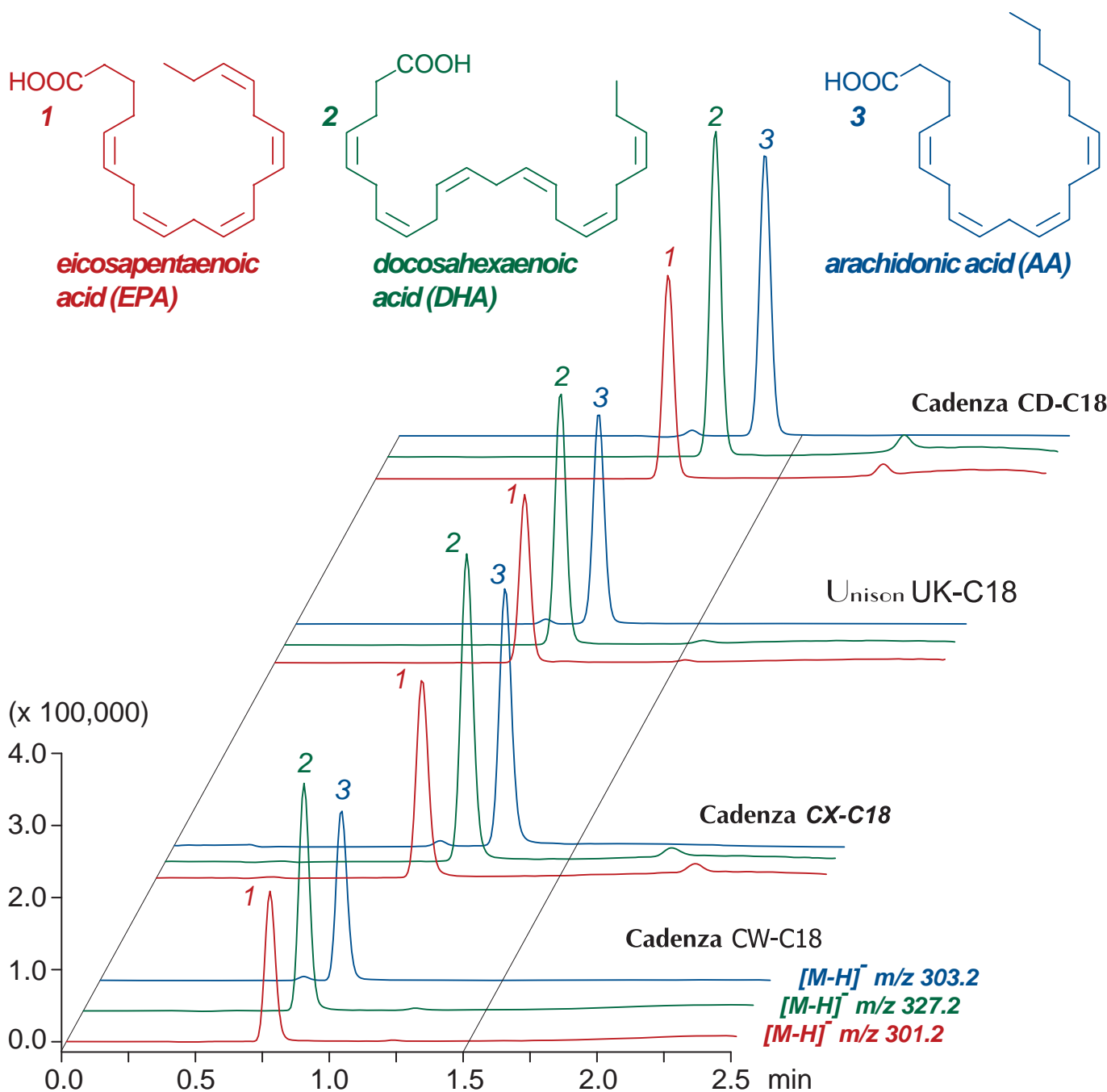

50 x 2 mm

A: 50mM ammonium formate

B: acetonitrile

85-100%B (0-1.7min), 100%B (1.7-2.5min)

0.3mL/min (3-4MPa), 37°C, MS (ESI, SIM, Negative)

5μL (1μg/mL, ACN/20mM HCOONH₄ = 50/50)