



Unison UK-C1  
Unison UK-C8  
Unison UK-C18

75 x 4.6 mm

Application

**Levofloxacin (new quinolone, antibiotics)**

レボフロキサシン (ニューキノロン系抗生物質)

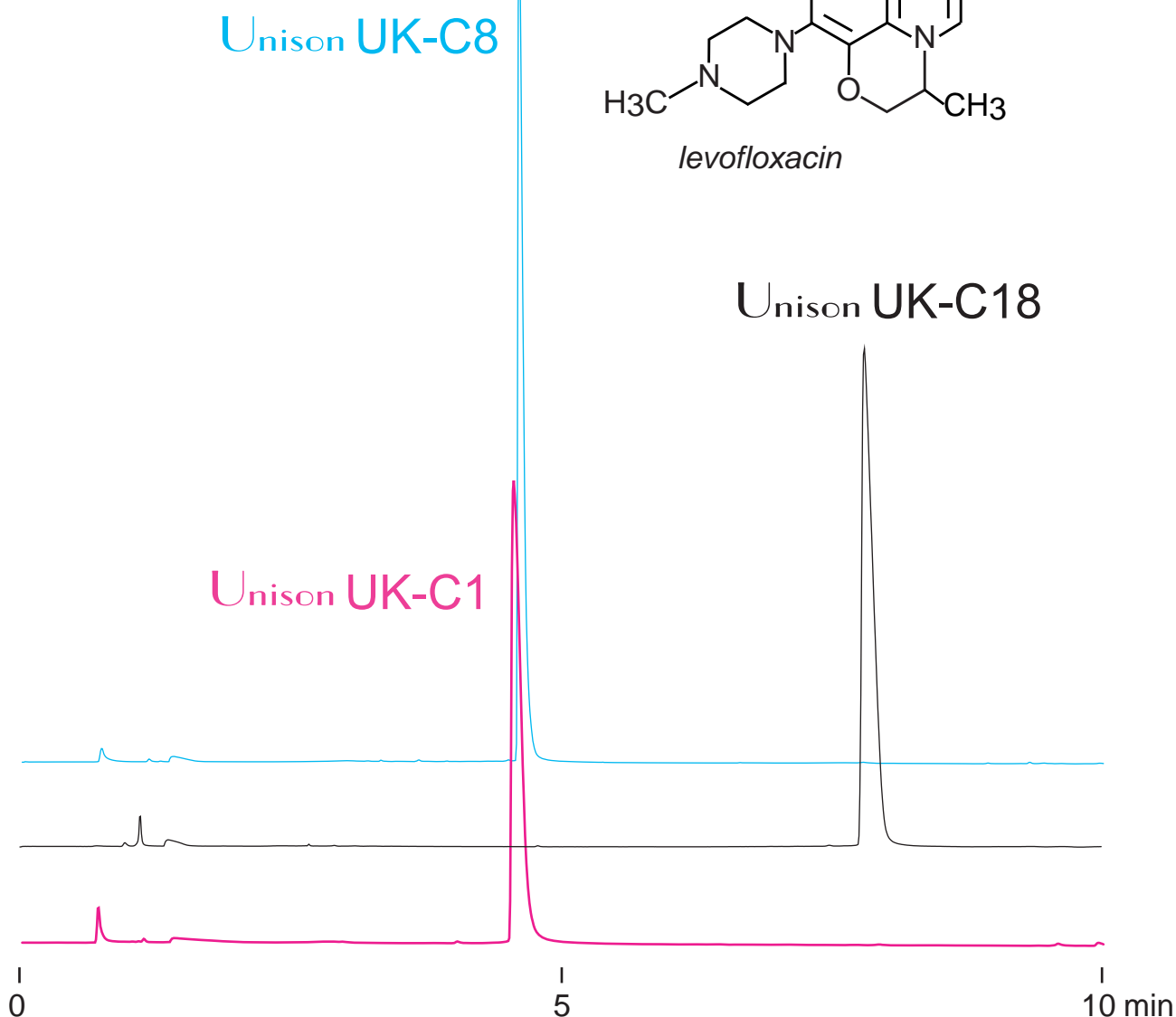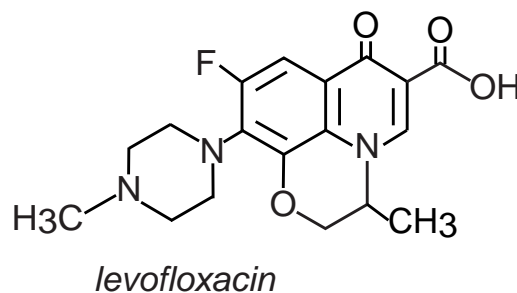

75 x 4.6 mm

A: water / formic acid= 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-50%B (0-10min)

1.0 mL/min (7MPa), 37°C

2  $\mu$ L (0.1mg/mL, 1%COOH), 260nm