


Unison UK-C1  
 Unison UK-C8  
 Unison UK-C18

50 x 3 mm

Application

Rizatriptan benzoate (treatment of migraine headaches)

リザトリプタン 安息香酸塩 (片頭痛薬)

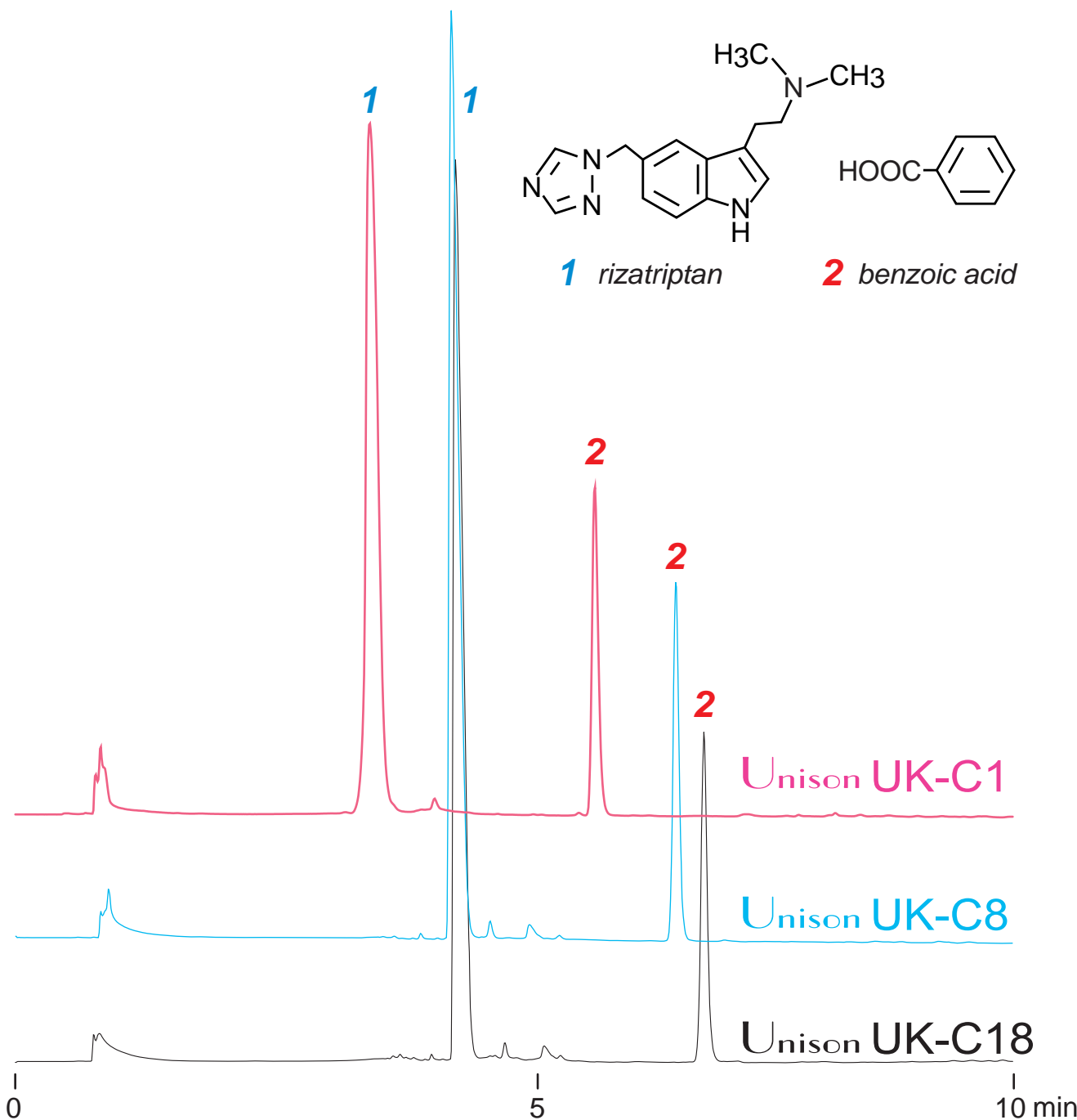

50 x 3 mm

A: water / formic acid= 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-50%B (0-10min)

0.4 mL/min (5MPa), 37°C

2 µL (1mg/mL, 1%HCOOH), 260nm