

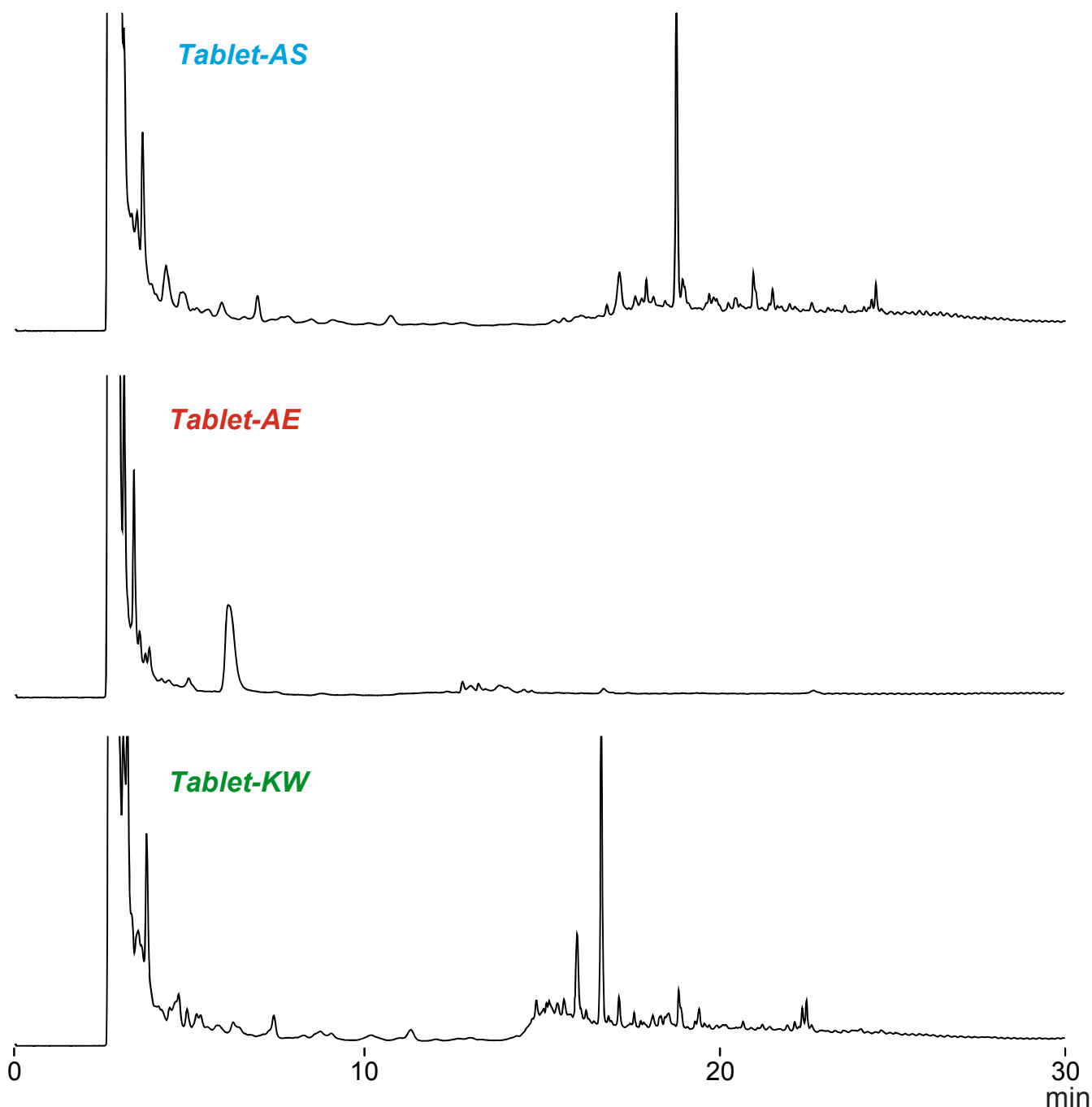
Presto FF-C18

150 x 4.6 mm

Application

Non-porous ODS : Beer yeast tablet extracts

ノンポーラスODS:ビール酵母錠剤成分の逆相分離

*Presto FF-C18, 150 x 4.6 mm**A: water / acetic acid = 100 / 0.1**B: acetonitrile / acetic acid = 100 / 0.1**0-100%B (0-30min), 0.4mL/min (12MPa)**37°C, 280nm, 10µL (1 tablet / 1mL of 100mM ammonium formate)*