


Intrada Organic Acid

150 x 2 mm

Application

LC-MS: Nelfinavir mesylate (antiretroviral drug)

ネルフィナビル メシル酸塩(抗HIV薬)のLC-MS分析

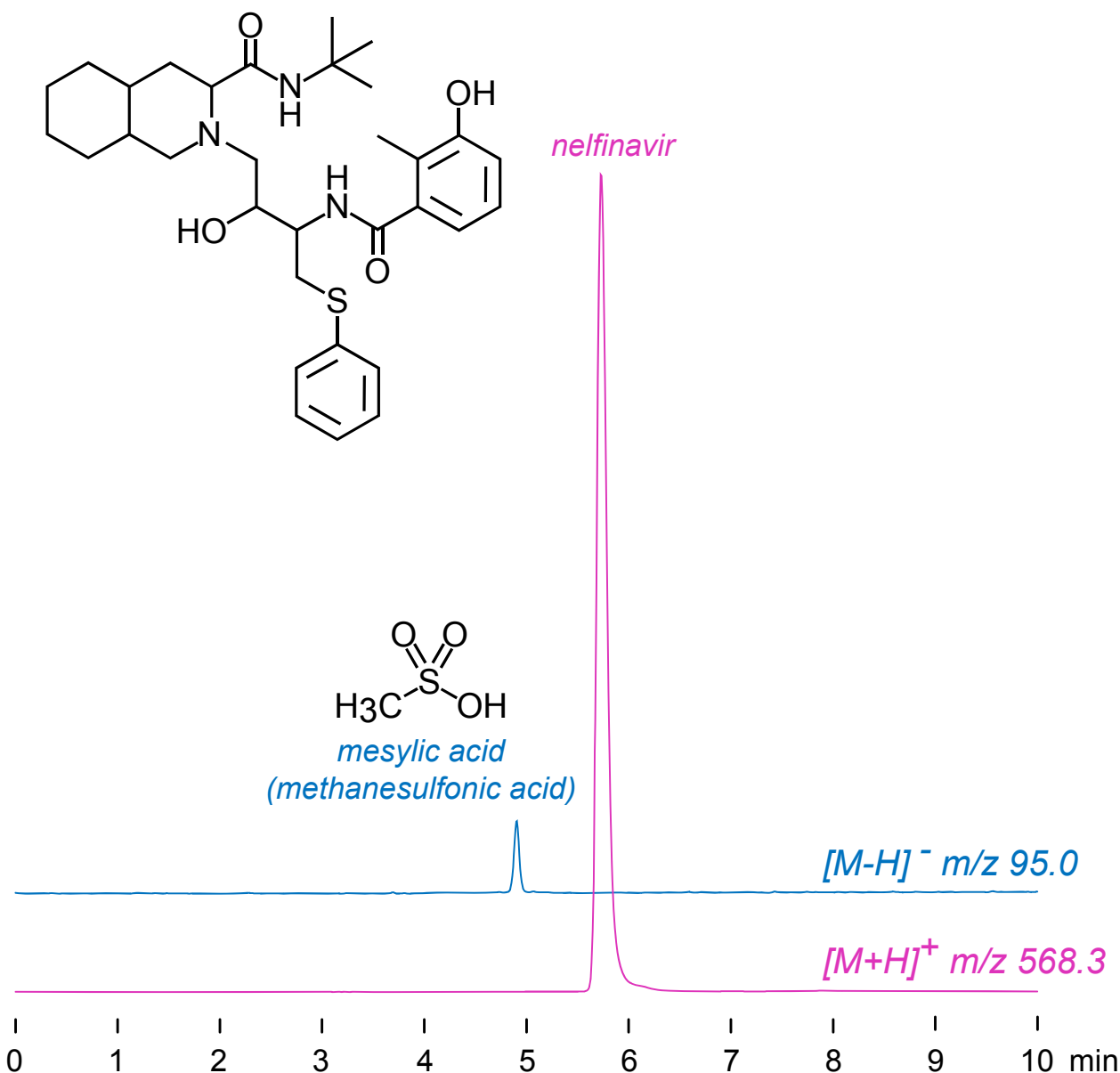

Intrada Organic Acid, 150 x 2 mm

A: acetonitrile / water / HCOOH = 10 / 90 / 0.1

B: acetonitrile / 100mM HCOONH₄ = 60 / 40

0%B (0-1 min), 0-100%B (1-7 min), 100%B (7-10 min)

0.2mL/min (12MPa), 37°C, MS (ESI, SIM)

5µL (5µg/mL, ACN/ water/ HCOOH = 10/ 90/ 0.1)