


Scherzo SS-C18
 Scherzo SM-C18
 Scherzo SW-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Amlodipine besilate (Antidepressant)


アムロジピン ベシル酸塩 (血圧降下薬)

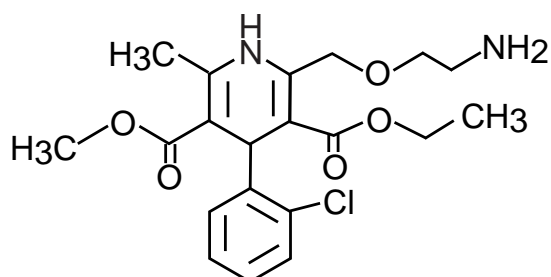

1. amlodipine

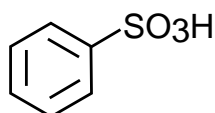

2. benzenesulfonic acid

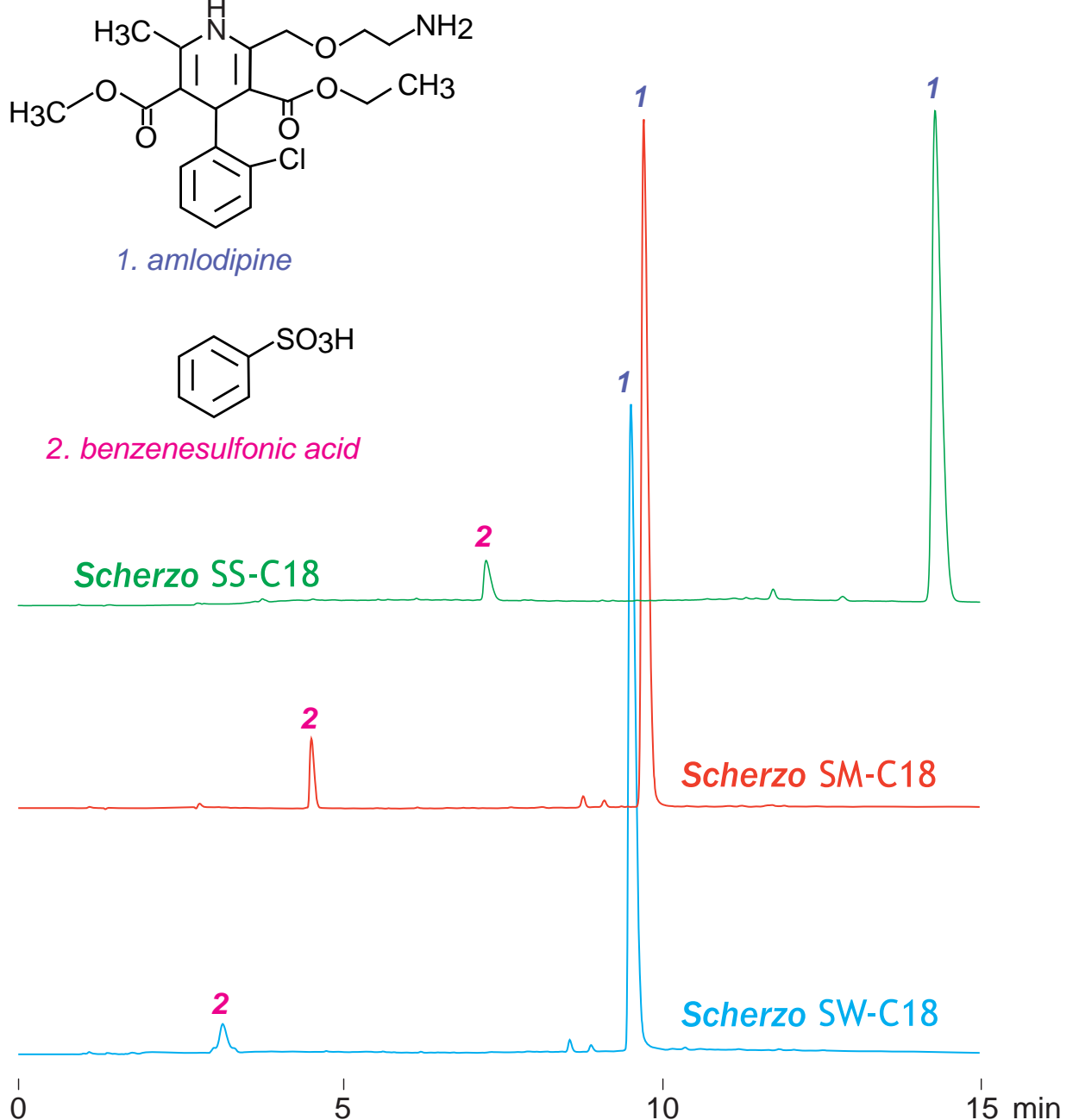

100 x 4.6 mm
 A: 10mM ammonium formate
 B: acetonitrile / 100mM ammonium formate= 70 / 30
 0-100%B (0-10min), 100%B (10-15min)
 1.0mL/min (8MPa), 37°C, 260nm
 10µL (0.5mg/mL, methanol / water / formic acid = 50 / 50 / 1)