



Scherzo SS-C18  
 Scherzo SM-C18  
 Scherzo SW-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Glycyrrhizic acid ( Anti-inflammatory )

グリチルリチン酸 ( 抗炎症薬 )

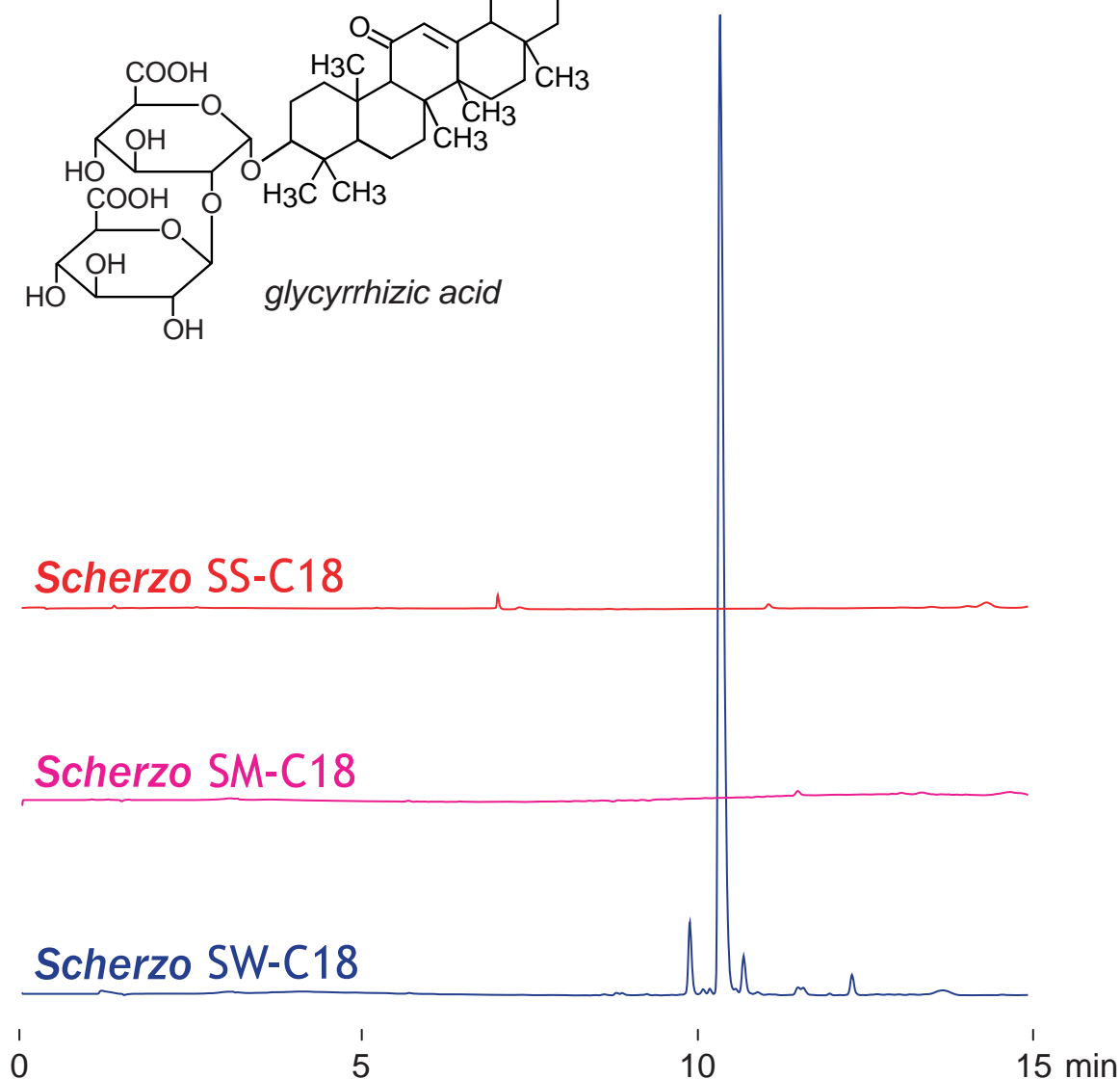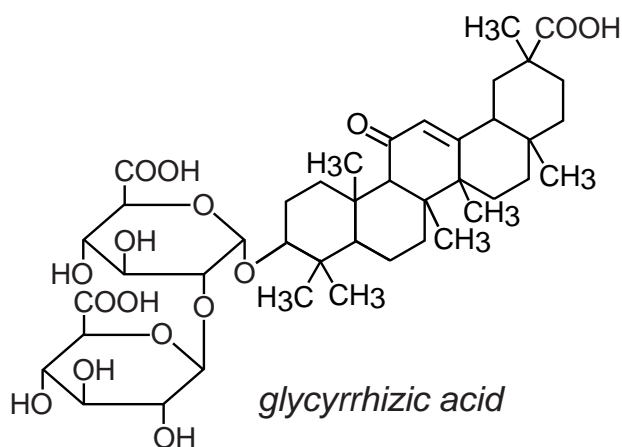

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / 100mM ammonium formate = 70 / 30

0-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (1mg/mL, 100mM ammonium formate)