


Scherzo SS-C18  
Scherzo SM-C18  
Scherzo SW-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

## Tranexamic acid ( Anti-inflammatory )

トラネキサム酸 ( 抗炎症薬 )

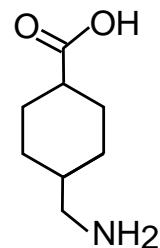

tranexamic acid

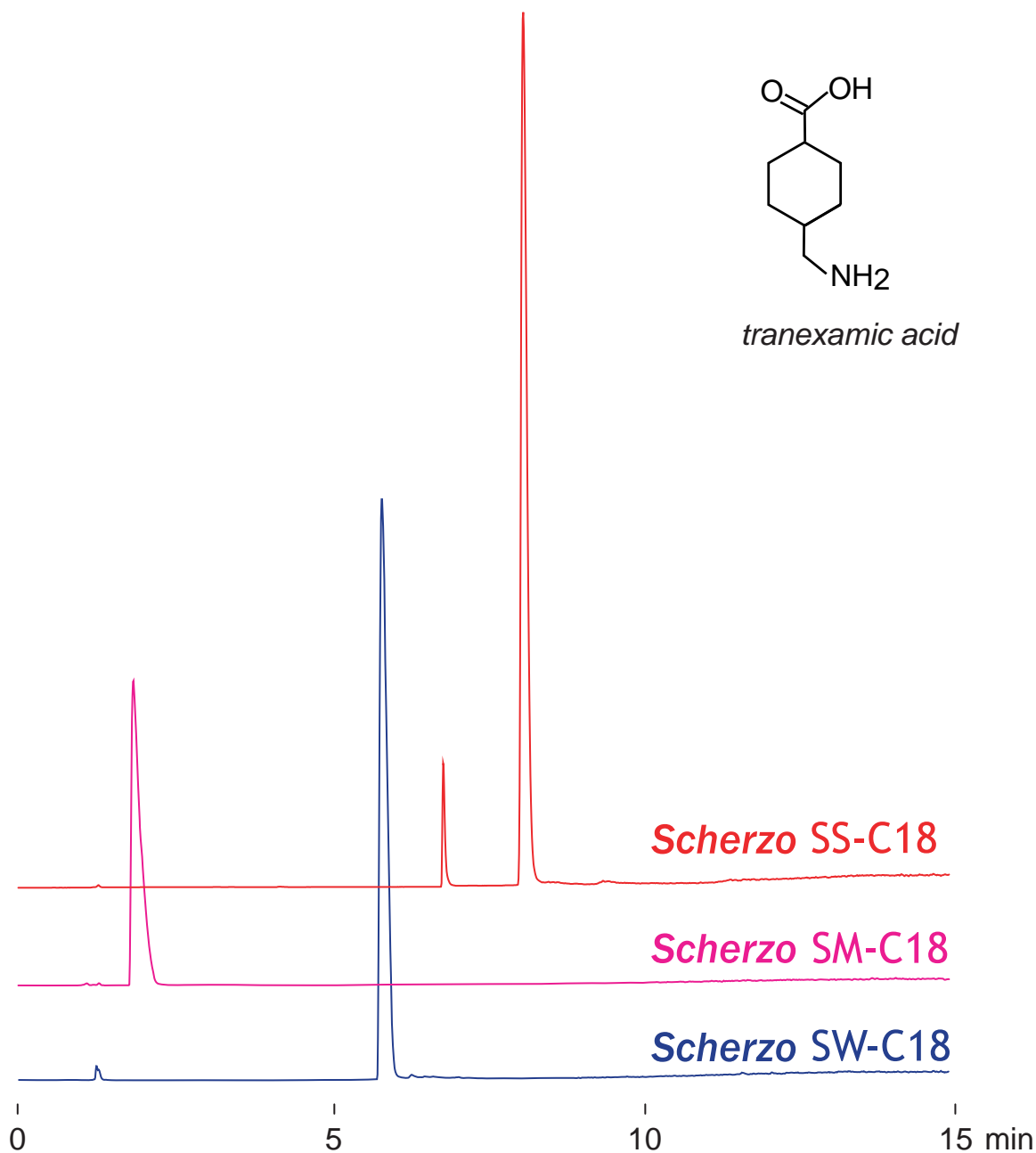

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / 100mM ammonium formate = 50 / 50

0%B (0-2min), 0-100%B (2-12min), 100%B (12-15min)

1.0mL/min (9MPa), 37°C,

ELSD (spray chamber 30°C, drift chamber 70°C)

10 $\mu$ L (1mg/mL, 100mM ammonium formate)