



Unison UK-C1

100 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Distearoylphosphatidylcholine (DPSC)

ジステアロイルホスファチジルコリン (脂質二重膜成分) の LC-MS 分析

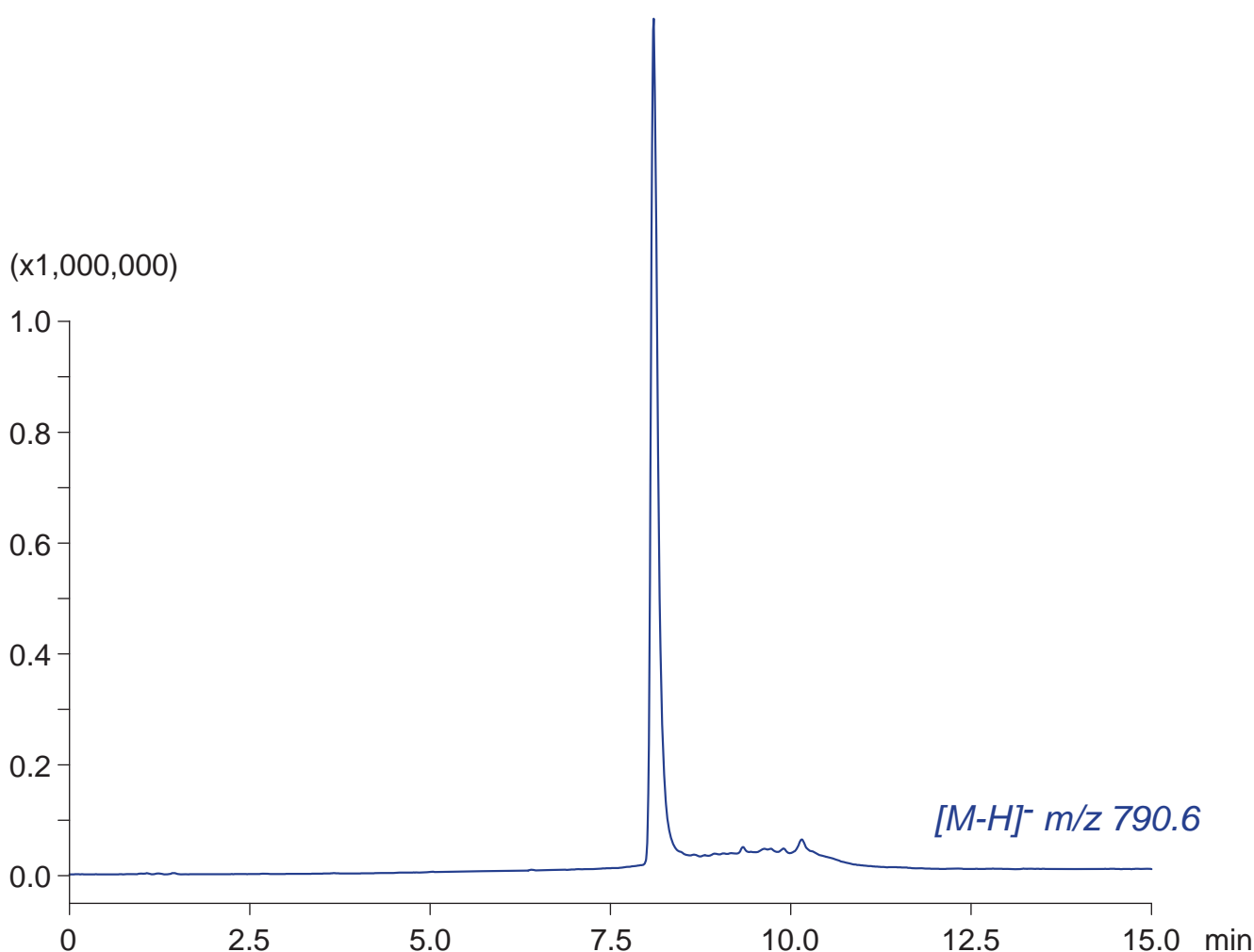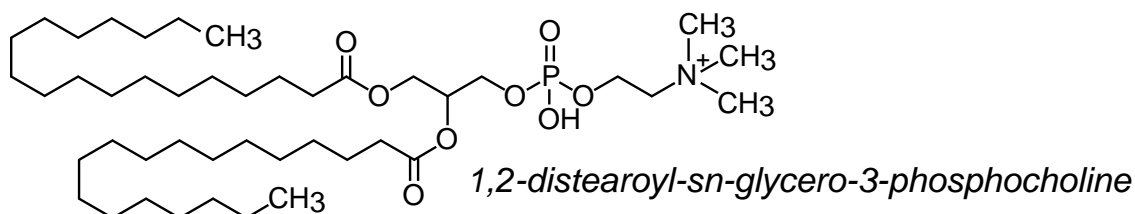

Unison UK-C1, 100 x 3 mm

A: water / formic acid = 100 / 0.3

B: methanol / formic acid = 100 / 0.3

60-95%B (0-10min), 95%B (10-15min)

0.4mL/min (7MPa), 60°C, MS (ESI, SIM, Positive)

10μL (1μg/mL, EtOH / water / HCOOH = 50 / 50 / 0.1)