



Intrada Organic Acid

150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS: Sucralfate (Antiulcer drug)

スクラルファート (胃炎・消化性潰瘍治療薬) の LC-MS 分析

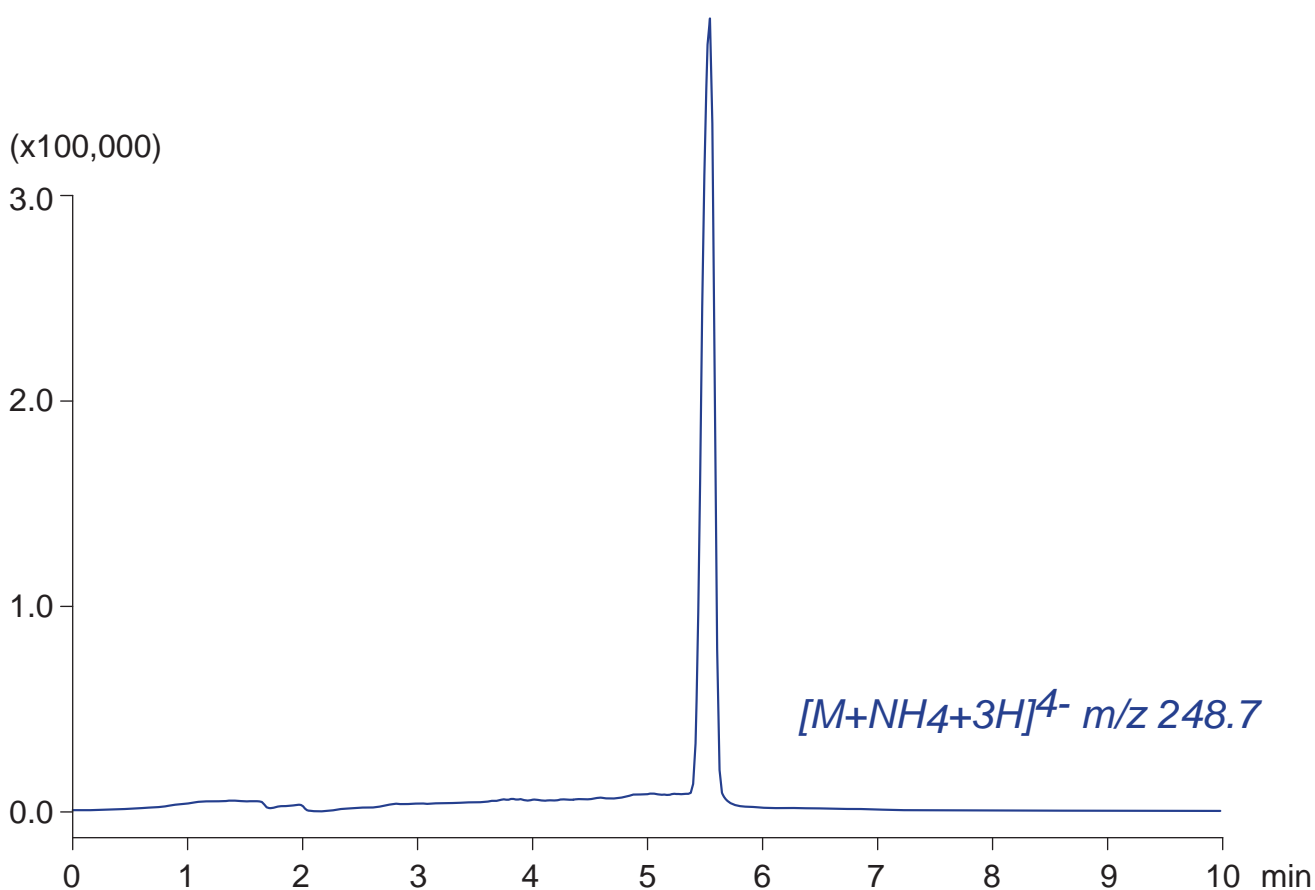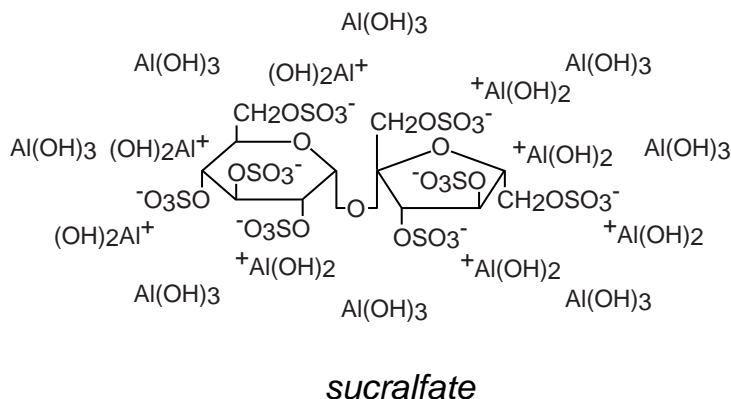

Intrada Organic Acid, 150 x 2 mm

A: methanol / 100mM ammonium formate = 95 / 5

B: methanol / 100mM ammonium formate = 10 / 90

0%B (0-1min), 0-100%B (1-7min), 100%B (7-10min)

0.2mL/min (9MPa), 37°C, MS (ESI, SIM, Negative)

5µL (1mg/mL, 0.02M NaOH)