



Scherzo SS-C18
 Scherzo SM-C18
 Scherzo SW-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Picosulfate (Laxative)

ピコスルファート (緩下薬)

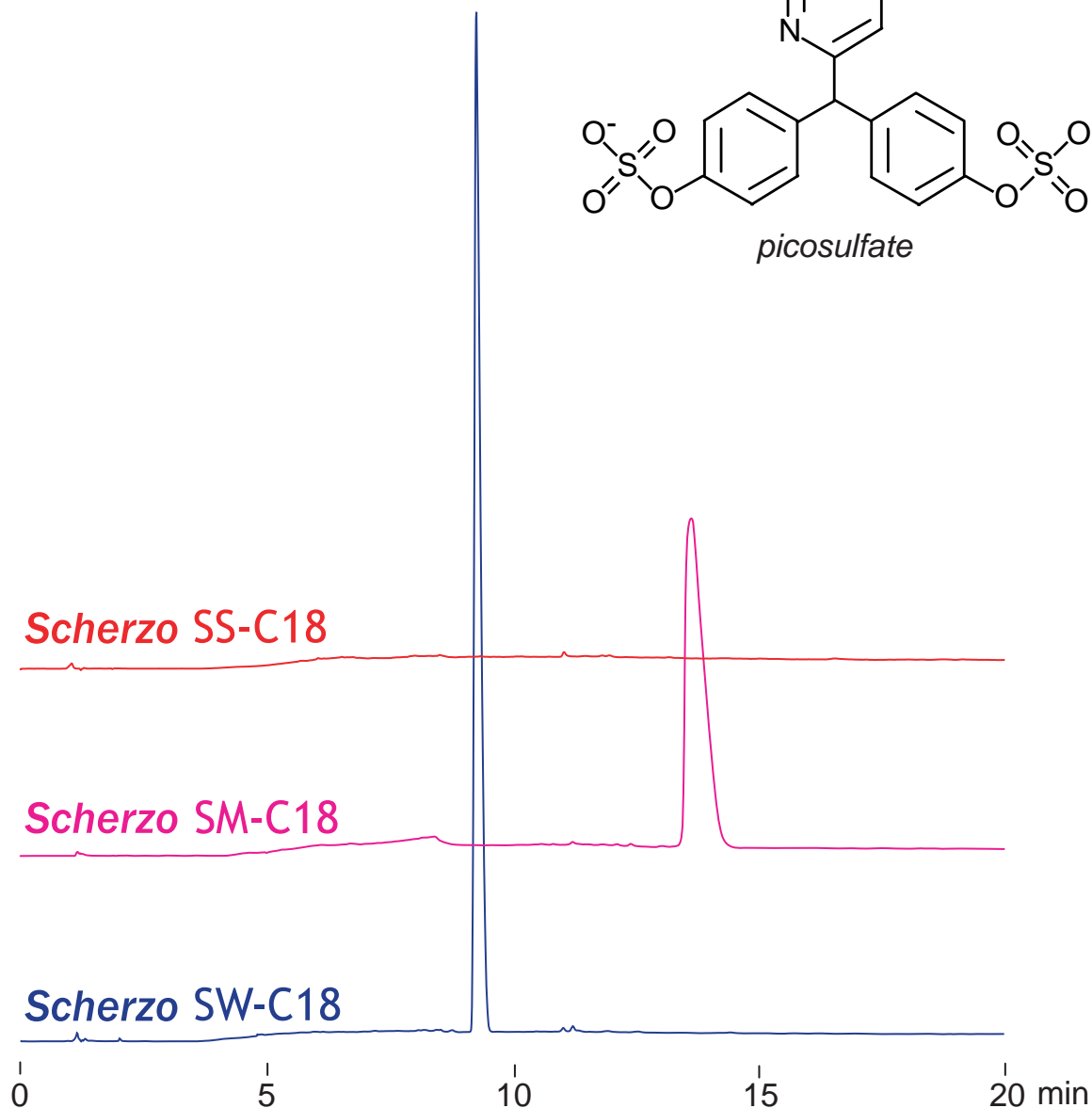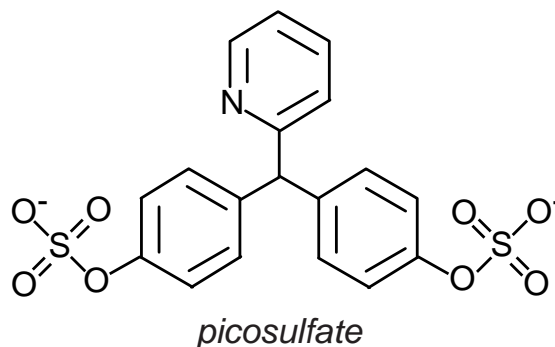

A: 5mM ammonium formate
 B: acetonitrile / 100mM ammonium formate = 30 / 70
 0%B (0-2min), 0-100%B (2-10min), 100%B (10-20min)
 1.0mL/min (9MPa), 37°C, 260nm
 10µL (1mg/mL, 100mM ammonium formate)