


Presto FF-C18

150 x 4.6 mm


HPLC COLUMN APPLICATION

Oligonucleotides ( Oligothymidylic acid )

オリゴヌクレオチド (オリゴチミジル酸)

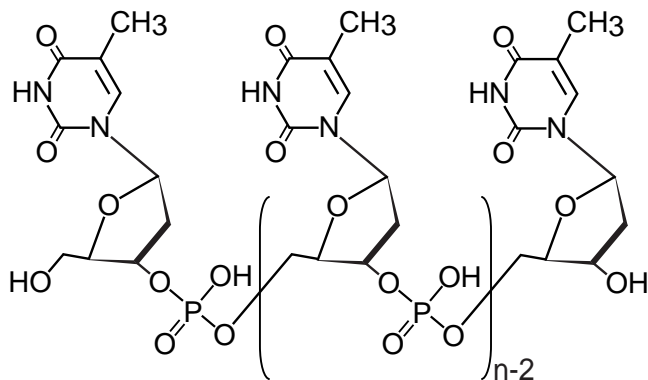

*oligothymidylic acid*  
*d(pT)<sub>12-18</sub>*

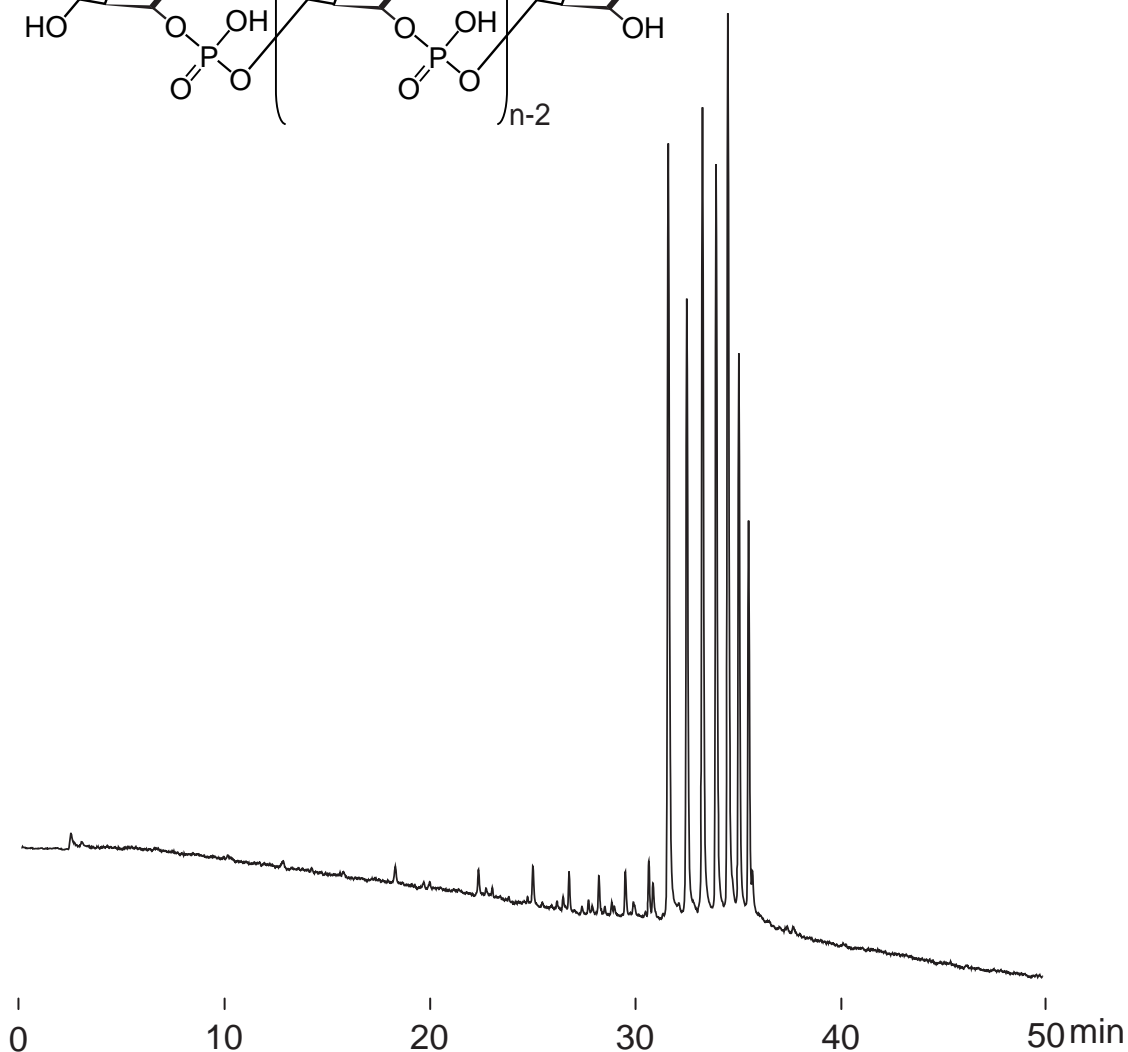

150 x 4.6 mm

A: 100mM ammonium formate

B: methanol / 100mM ammonium formate = 90 / 10

5-35%B (0-50min)

0.4mL/min (15MPa), 60°C, 260nm

4μL (20μg/mL, 100mM ammonium formate)