

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8  
Cadenza CX-C18


Unison UK-C18  
Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

## Fluvastatin ( Hypolipidemic agent )

フルバスタチン (高脂血症治療薬)

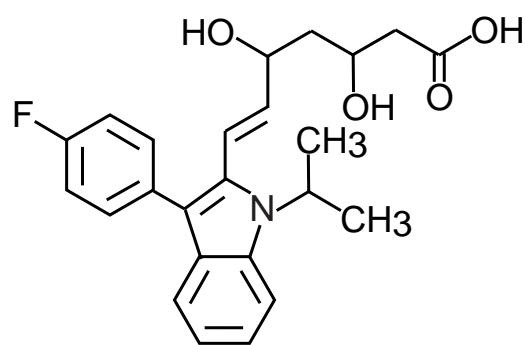

fluvastatin

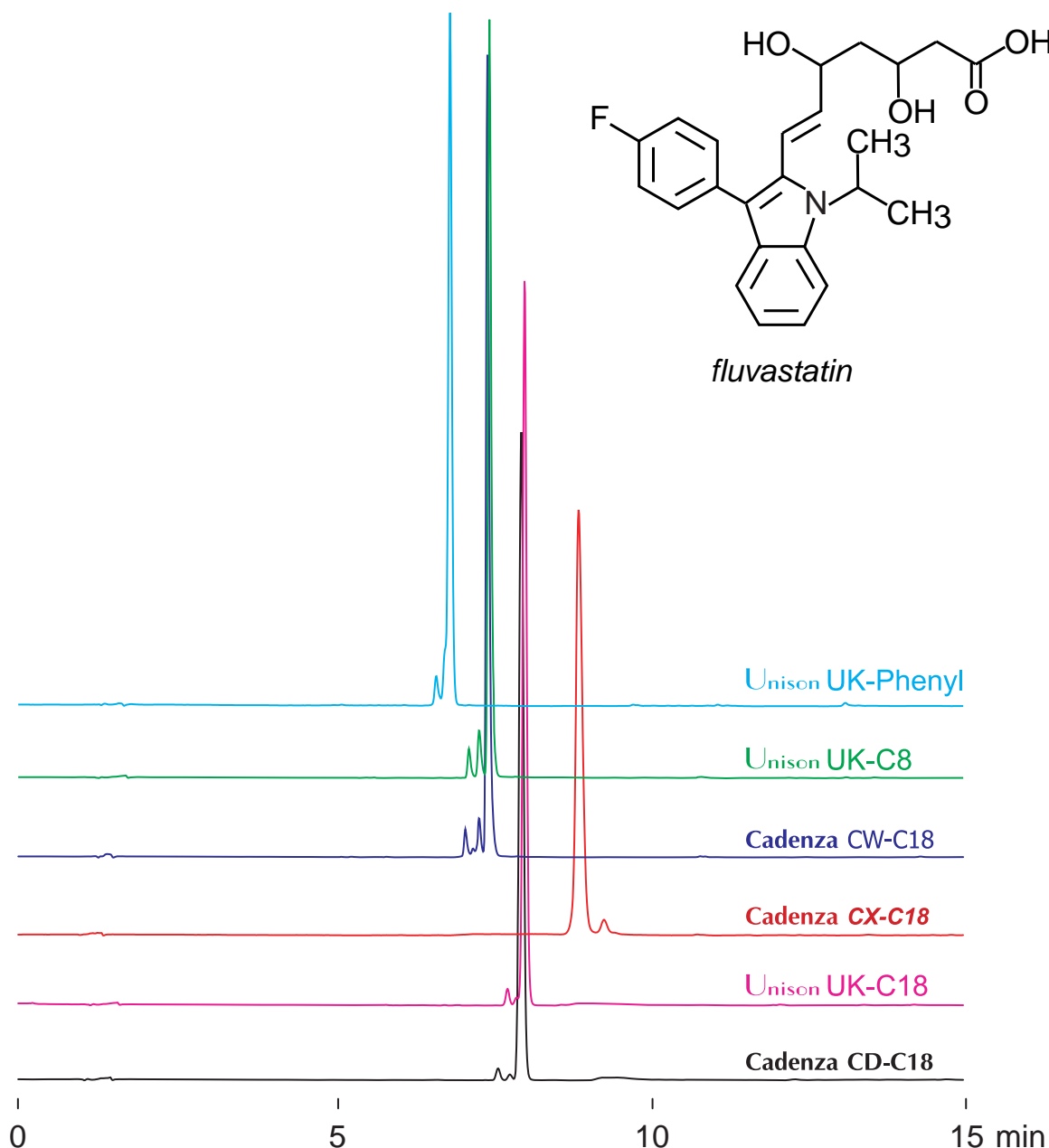

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate= 90 / 10

20-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.25mg/mL, methanol)