


Unison UK-C8

150 x 3 mm

HPLC COLUMN APPLICATION

## LC-MS: Organophosphorus pesticides ( 1 )

有機リン系殺虫剤 ( 1 ) の LC-MS 分析

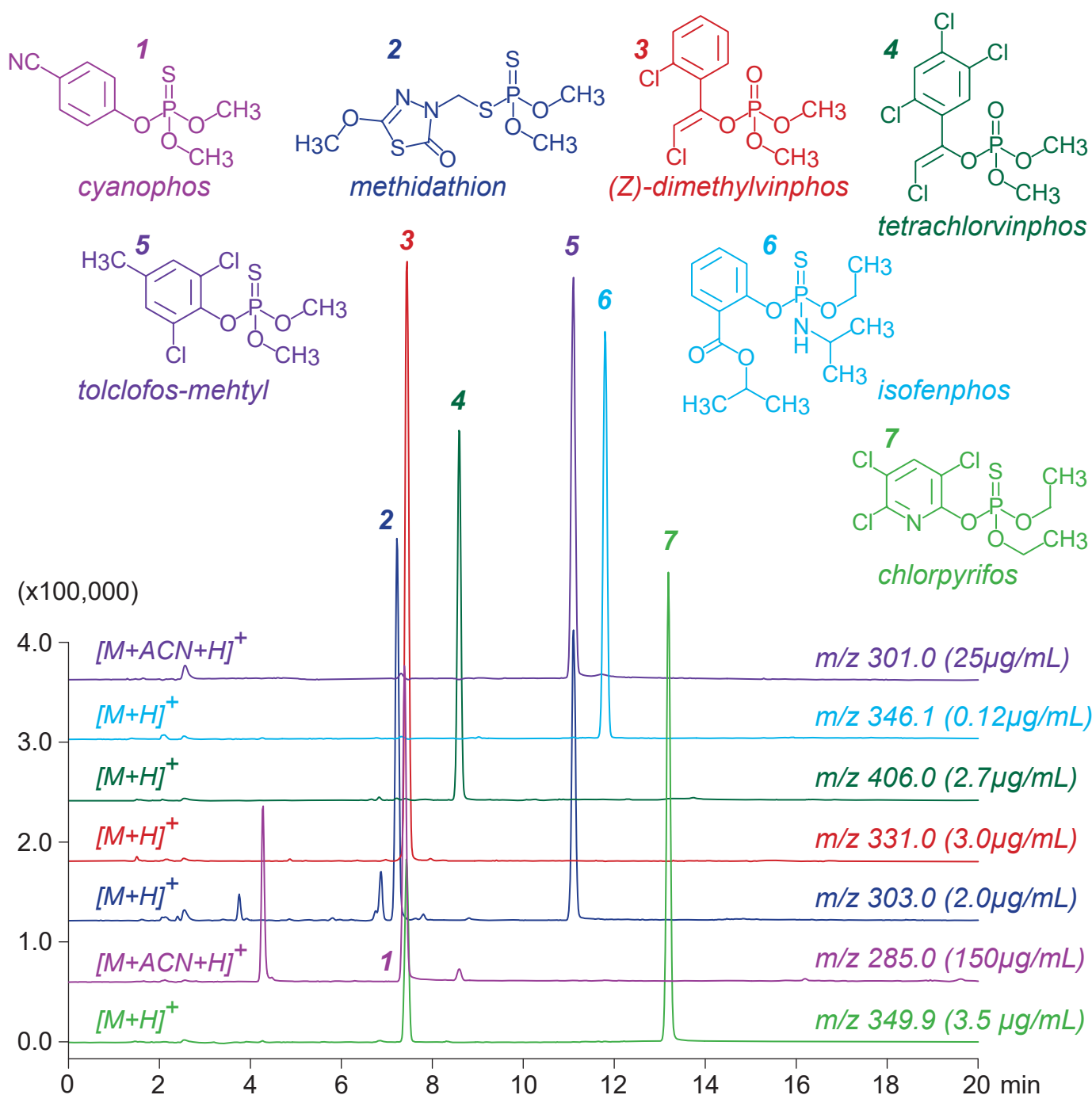

Unison UK-C8, 150 x 3 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

45-85%B (0-15min), 85%B (15-20min)

0.4mL/min (8MPa), 37°C, Single Quad. MS (ESI)

5µL (0.12-150µg/mL, ACN / water / HCOOH = 50 / 50 / 0.1)