


Unison UK-C8

150 x 3 mm

HPLC COLUMN APPLICATION

## LC-MS: Organophosphorus pesticides ( 3 )

有機リン系殺虫剤 ( 3 ) の LC-MS 分析

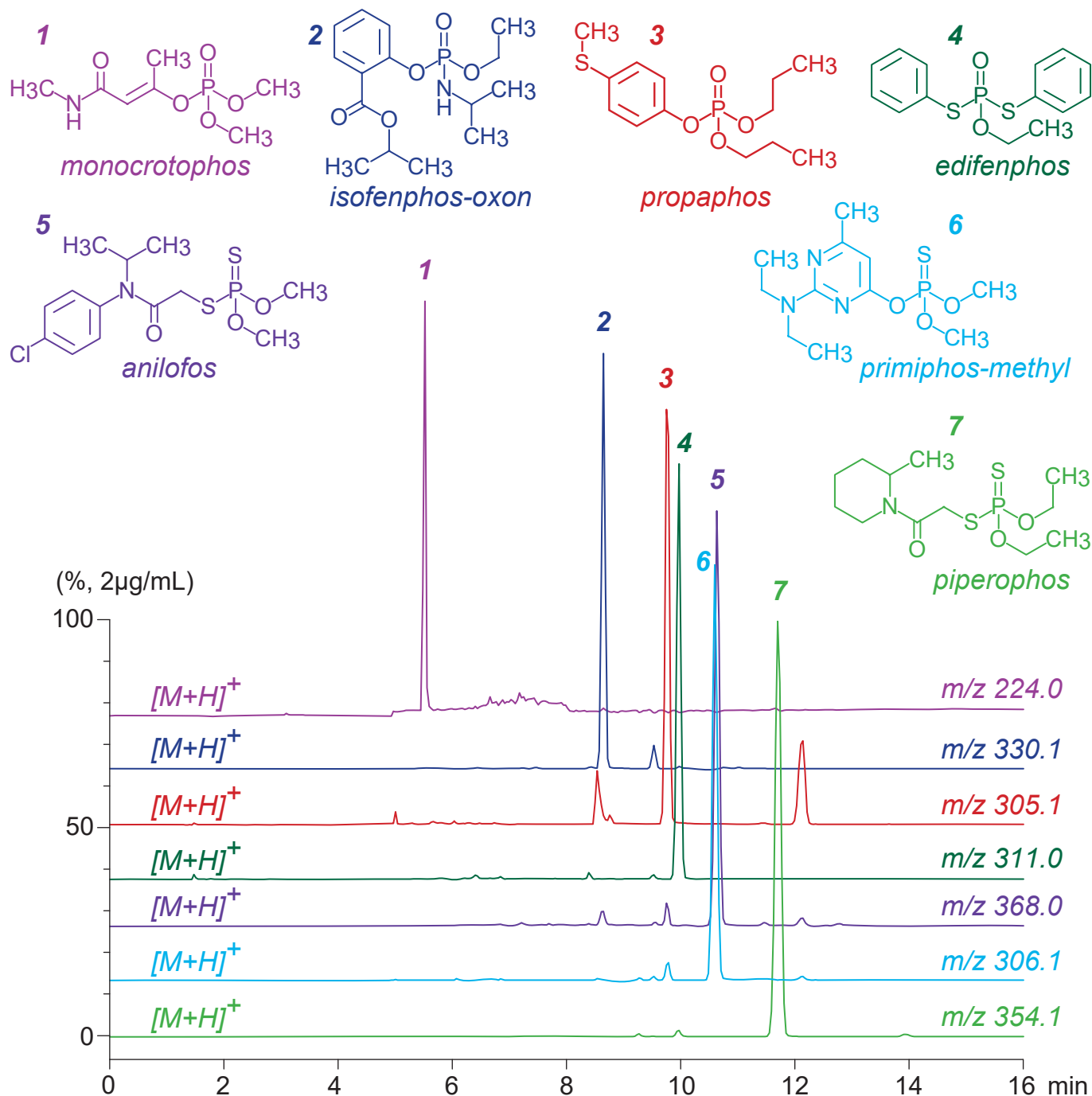

Unison UK-C8, 150 x 3 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-65%B (0-5min), 65%B (5-16min)

0.4mL/min (8MPa), 37°C, Single Quad. MS (ESI)

5μL (2μg/mL, ACN / water / HCOOH = 50 / 50 / 0.1)