


Scherzo SW-C18
 Scherzo SM-C18
 Scherzo SS-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Diphenadione (Indandione anticoagulant)

ジフェナジオン (インダンジオン系抗凝血薬)

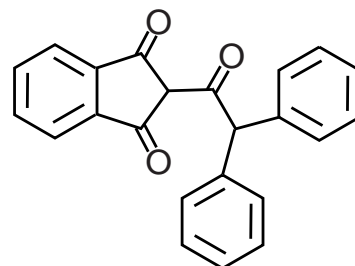

diphenadione
 (diphacinone)

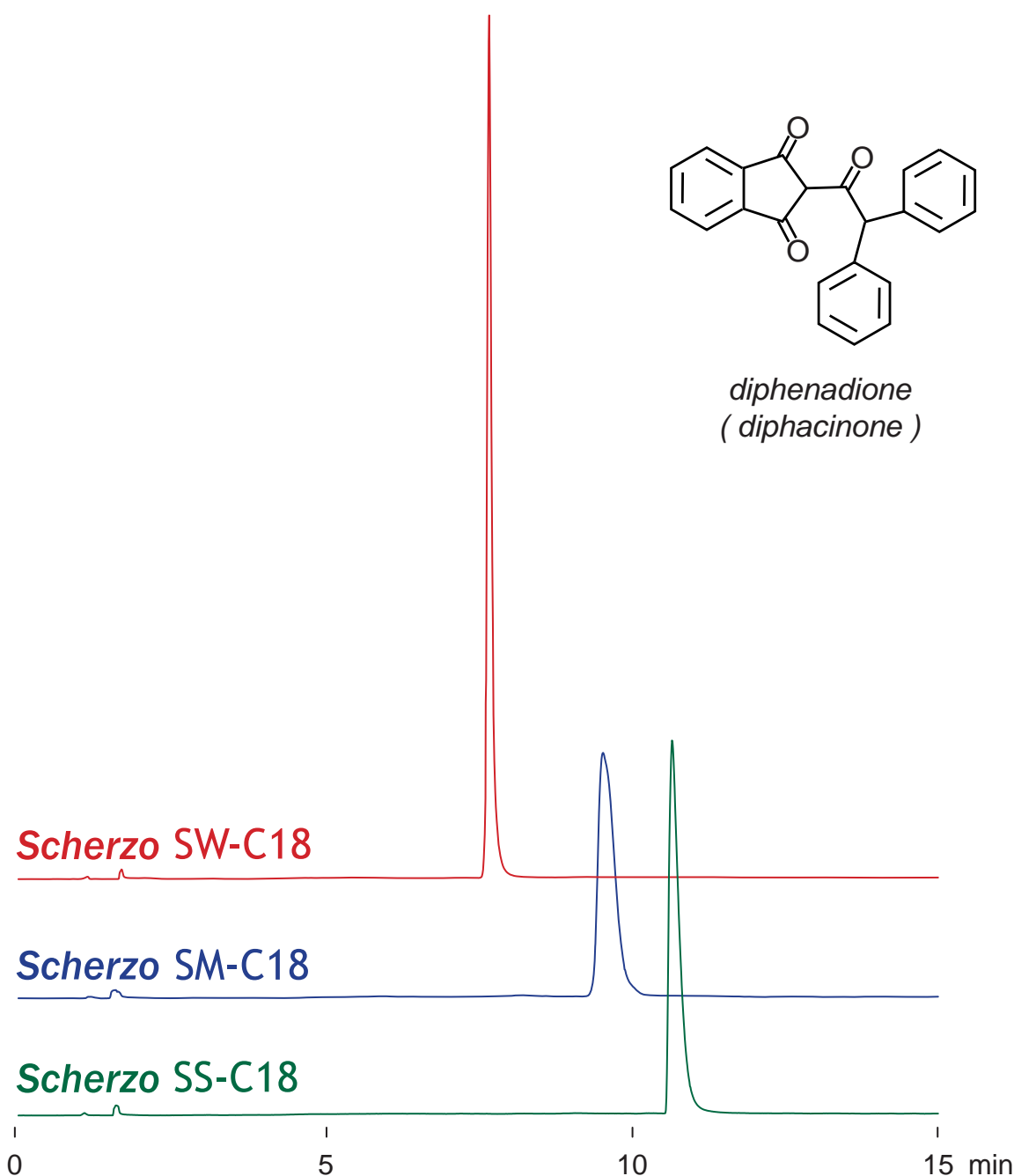

100 x 4.6 mm
 A: 5mM ammonium acetate
 B: acetonitrile / 100mM ammonium acetate= 70 / 30
 0-100%B (0-7min), 100%(7-15min)
 1.0 mL/min (10MPa), 37°C, 280nm
 10µL (0.2mg/mL, acetonitrile / 100mM ammonium acetate = 50 / 50)