

Cadenza CW-C18

Unison UK-C8

Unison UK-C18


100 x 4.6 mm

HPLC COLUMN APPLICATION


Carbazochrome ( Hemostatic )

Troxerutin ( Vascular protectant )


カルバゾクロム ( 止血薬 ), トロキセルチン ( 血管保護薬 )

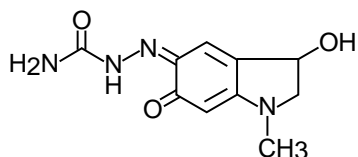

1 carbazochrome (0.05mg/mL)

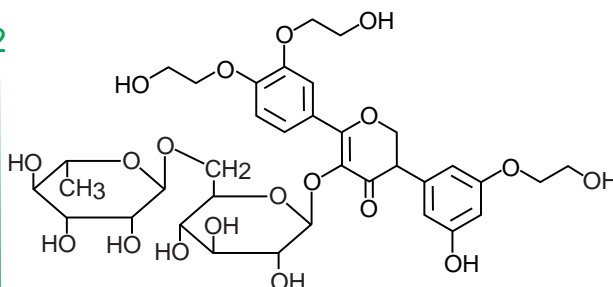

2 troxerutin (0.2mg/mL)

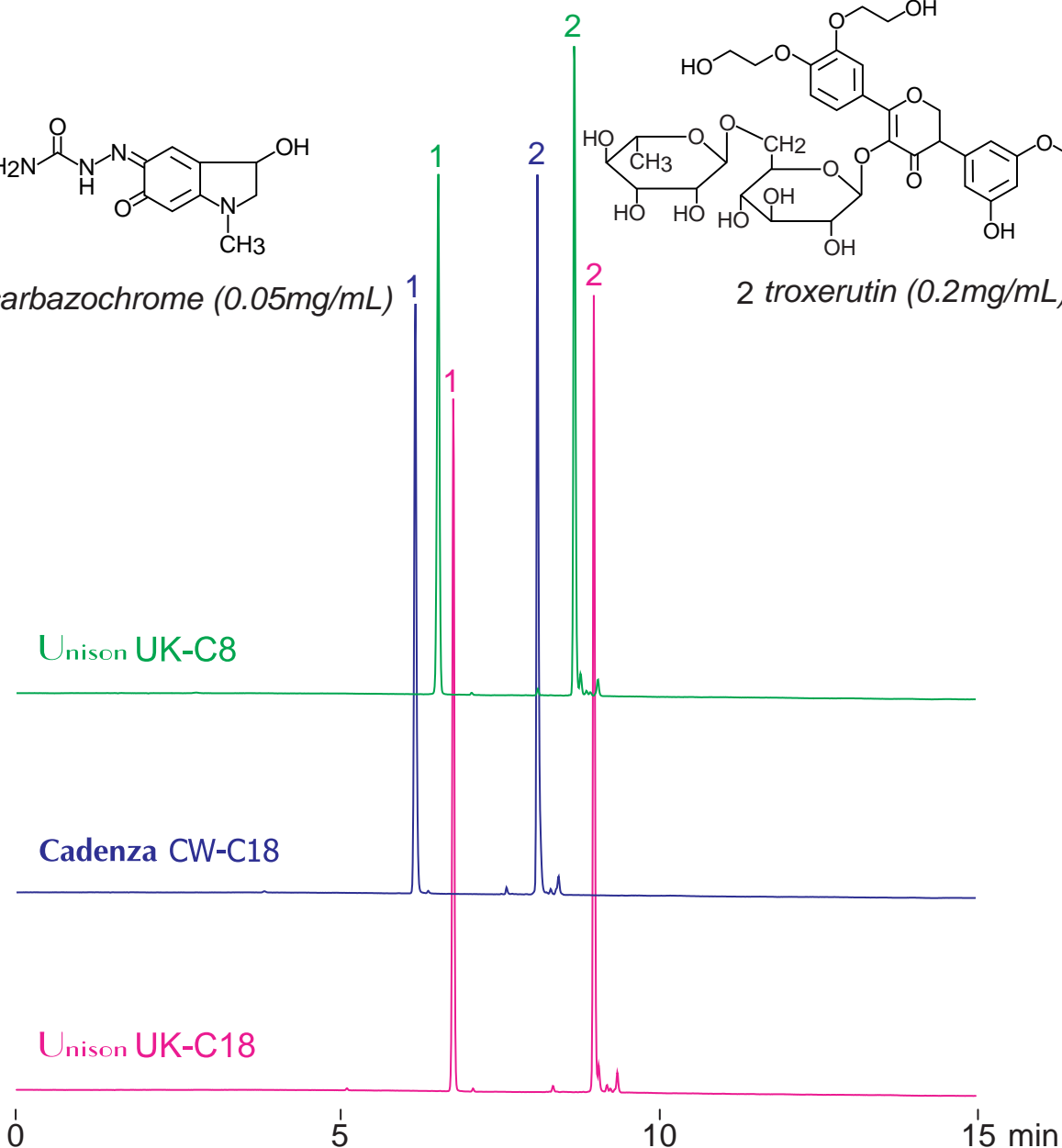

100 x 4.6 mm

A: 100mM ammonium acetate

B: 100mM ammonium acetate / acetonitrile = 10 / 90

0-75%B (0-15min)

1.0mL/min (12MPa), 37°C, 360nm

10µL (0.05-0.2mg/mL, 100mM ammonium acetate )