
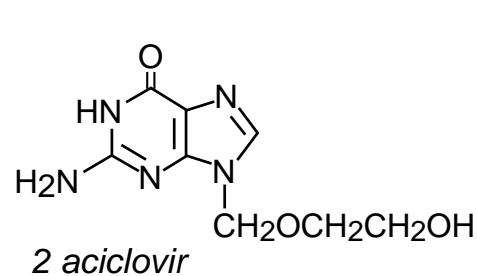


Unison UK-Amino

150 x 4.6 mm

HPLC COLUMN APPLICATION

Ganciclovir, Aciclovir ( Antiviral )  
 ガンシクロビル, アシクロビル ( 抗ウイルス薬 )

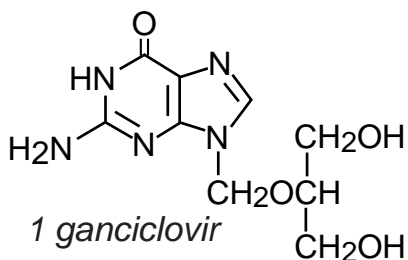

A: 10mM AcONH<sub>4</sub> / ACN = 10 / 90  
 B: 10mM AcONH<sub>4</sub>

A: 100mM AcONH<sub>4</sub> / ACN = 10 / 90  
 B: 100mM AcONH<sub>4</sub>

0 5 min 10 15

Unison UK-Amino, 150 x 4.6 mm  
 0-50%B (0-15min)  
 1.0mL/min (5MPa), 37°C, 260nm  
 10μL (0.1mg/mL, 100mM ammonium acetate )