



Scherzo SS-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Mercaptopurine (Antineoplastic)

メルカプトプリン (抗悪性腫瘍薬)

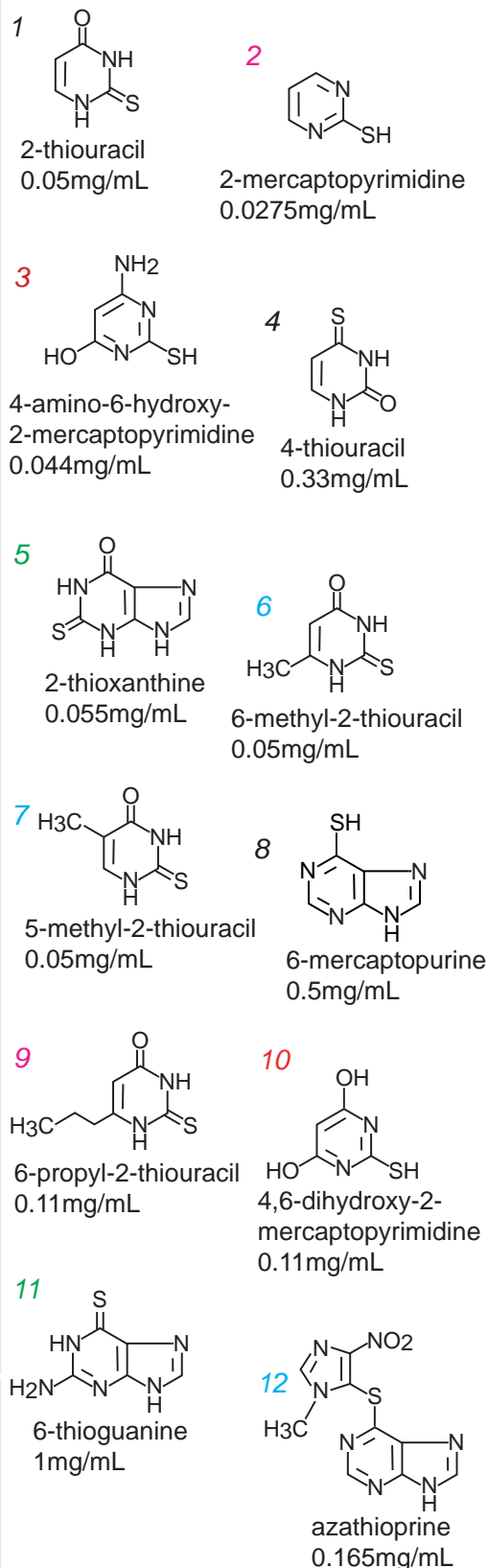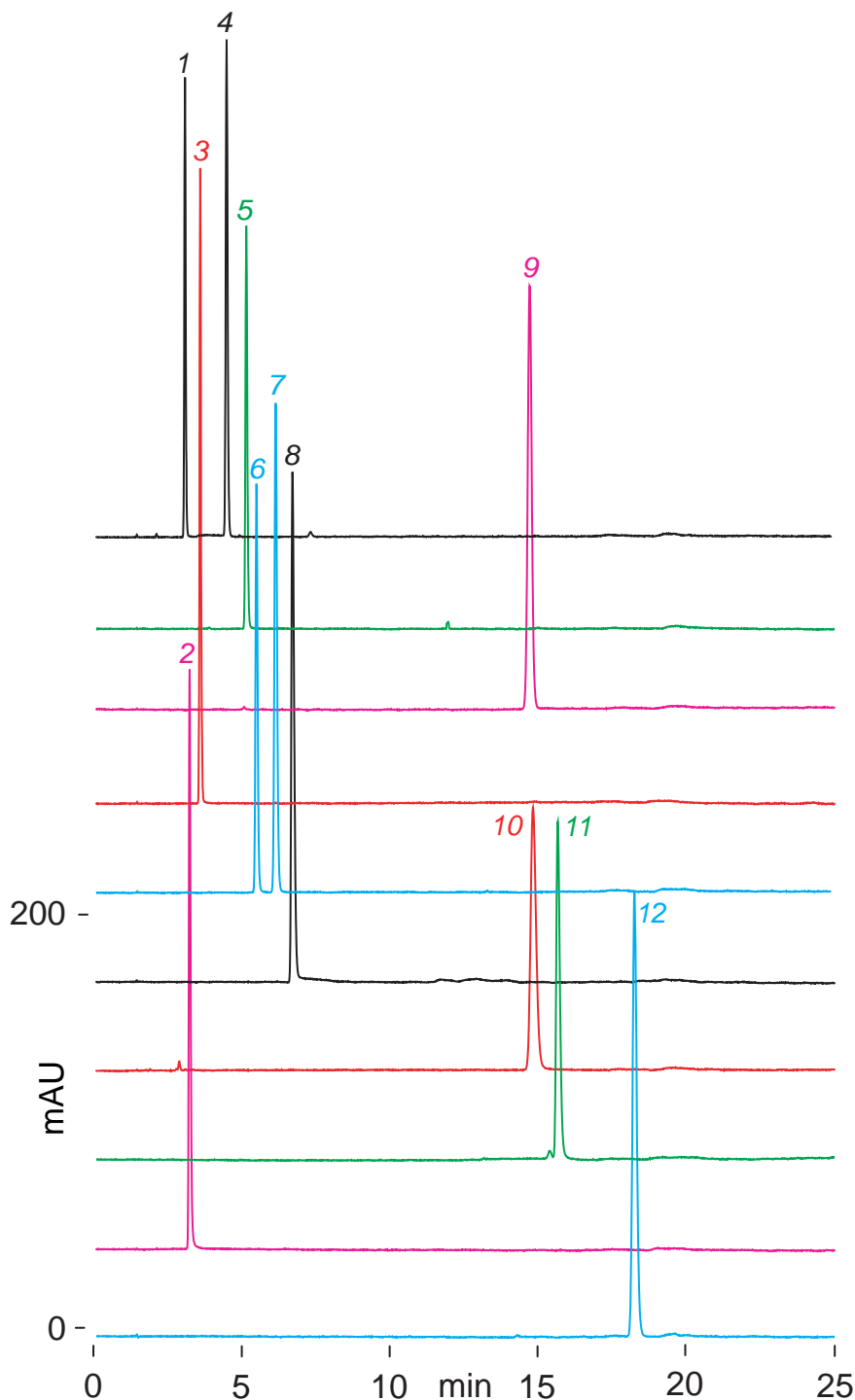

Scherzo SS-C18, 100 x 4.6 mm
 A: water / formic acid = 100 / 0.1
 B: 100mM ammonium formate / acetonitrile = 50 / 50
 0-25%B (0-25min)
 1.0mL/min (9MPa), 37°C, 280nm
 5µL, (33mM HCl (1-2, 4-12), 100mM HCl (3))