



Scherzo SM-C18

100 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Dobesilate (Vasoprotectives)

ドベシル酸 (血管保護薬) の LC-MS 分析

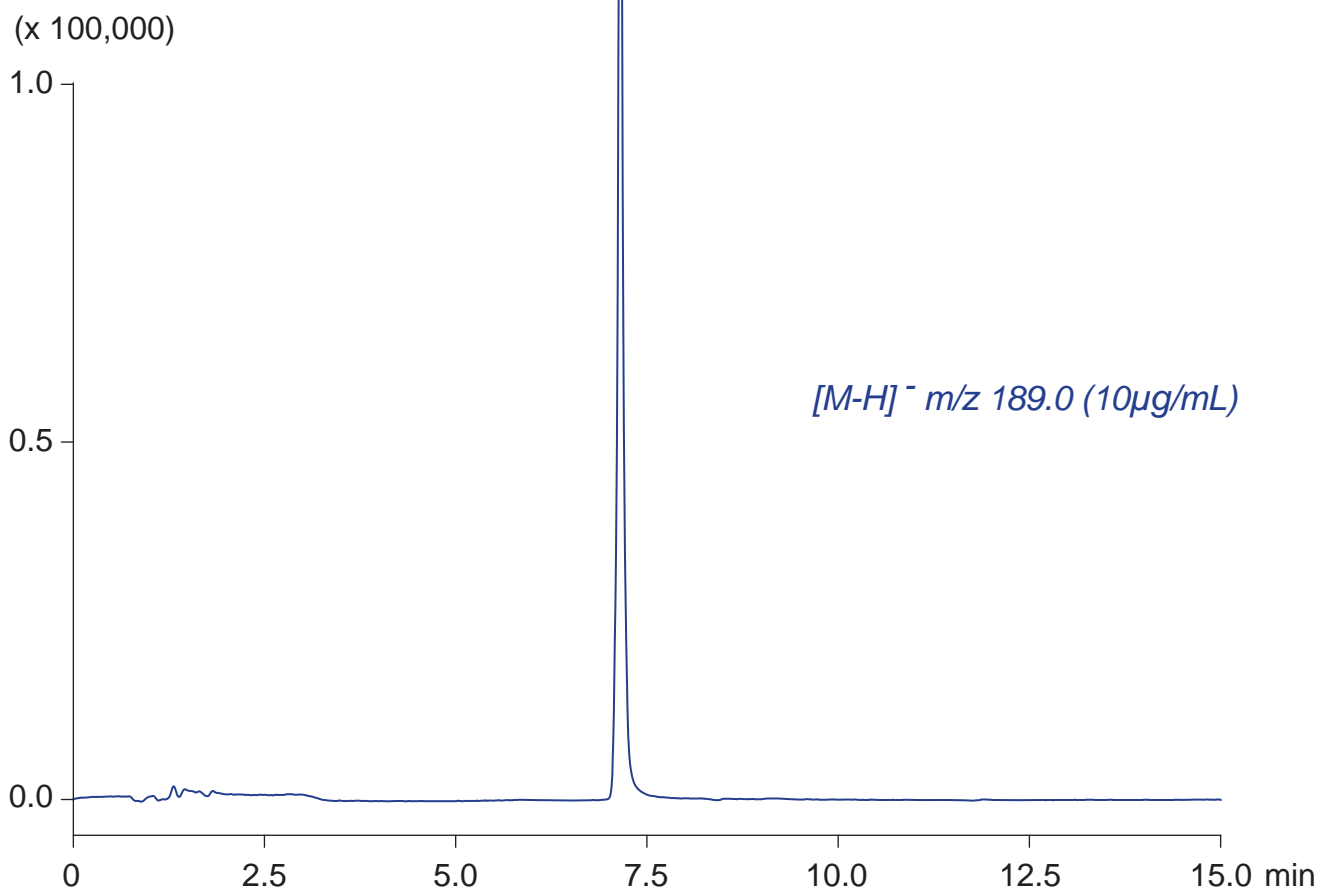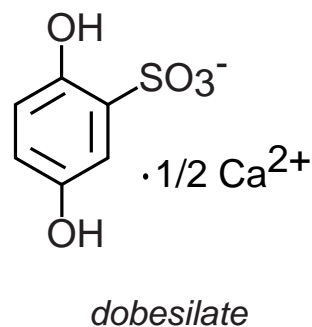

Scherzo SM-C18, 100 x 3 mm

A: acetonitrile / water / formic acid = 10 / 90 / 0.1

B: acetonitrile / 100mM ammonium formate = 50 / 50

0%B (0-1min), 0-100%B (1-10min), 100%B (10-15min)

0.4mL/min (7MPa), 37°C, Single Quad. MS (ESI)

10µL (10µg/mL, water / formic acid = 100 / 1)