

Nardis ND-NX MF


Nardis ND-NX

150 x 2 mm


HPLC COLUMN APPLICATION

Mononucleotides

モノヌクレオチド

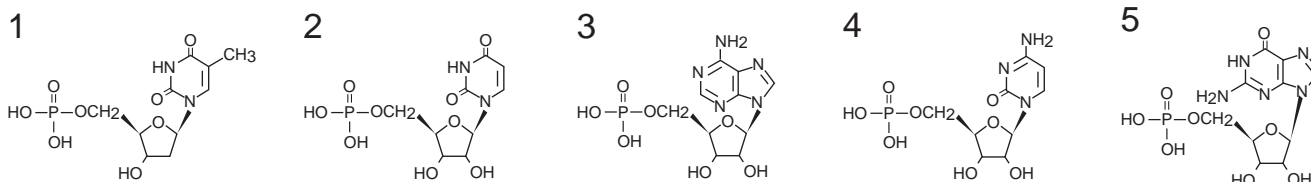

dTMP(0.2mg/mL) UMP(0.2mg/mL) AMP(0.2mg/mL) CMP(0.2mg/mL) GMP(0.2mg/mL)

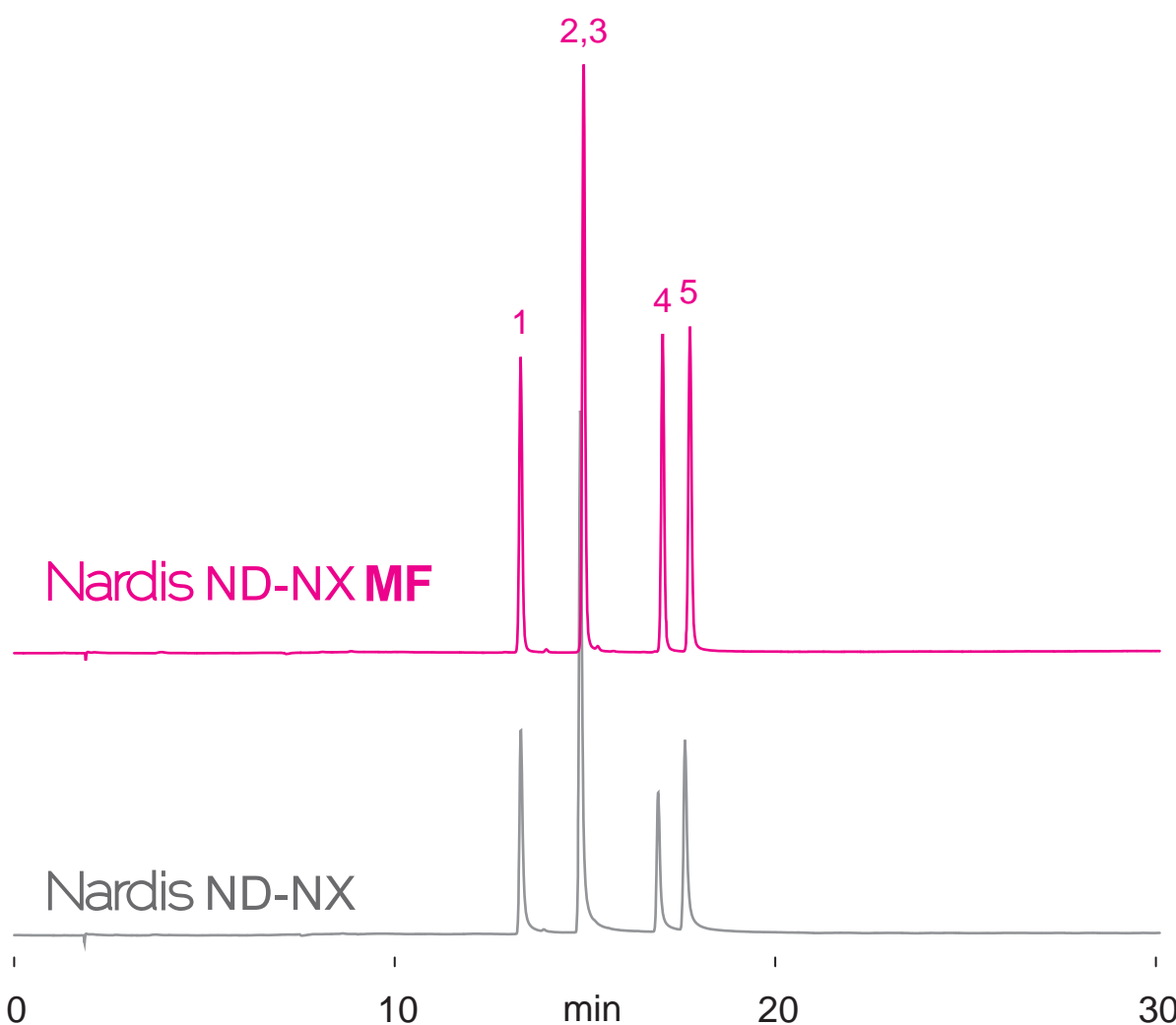

150 x 2 mm

A: acetonitrile

B: acetonitrile / 100mM ammonium acetate = 20 / 80

20-90%B (0-30min)

0.2mL/min (5MPa), 37°C, 260nm

2.0μL (200μg/mL, 100mM ammonium acetate)