



Unison UK-Silica  
Nardis ND-NX

150 x 4.6 mm

HPLC COLUMN APPLICATION

Benfotiamine (Vitamin B1 derivatives)

ベンフォチアミン (ビタミンB誘導体)

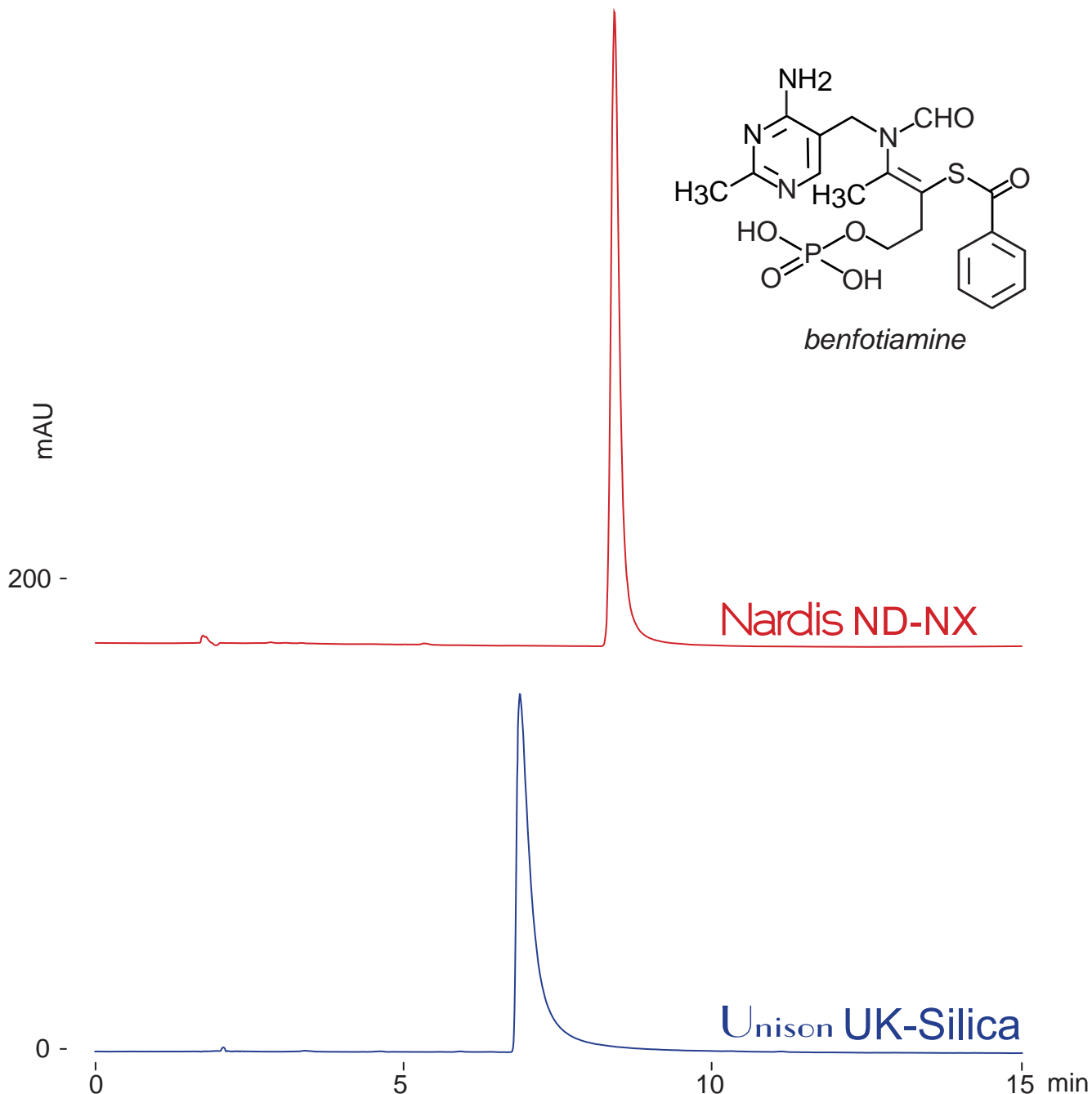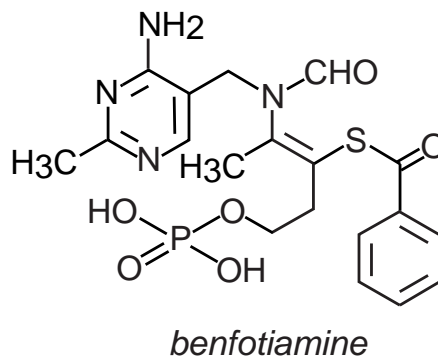

150 x 4.6 mm

A : acetonitrile / acetic acid = 100 / 0.5

B : acetonitrile / water / acetic acid = 50 / 50 / 0.5

20-100%B (0-15min)

1mL/min (7MPa), 37°C, 260nm

10µL (0.1mg/mL, water / acetic acid = 100 / 1)