


Unison UK-Phenyl Cadenza CL-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

250 x 4.6 mm

HPLC COLUMN APPLICATION

tert-Butylphenol isomers

tert-ブチルフェノール異性体

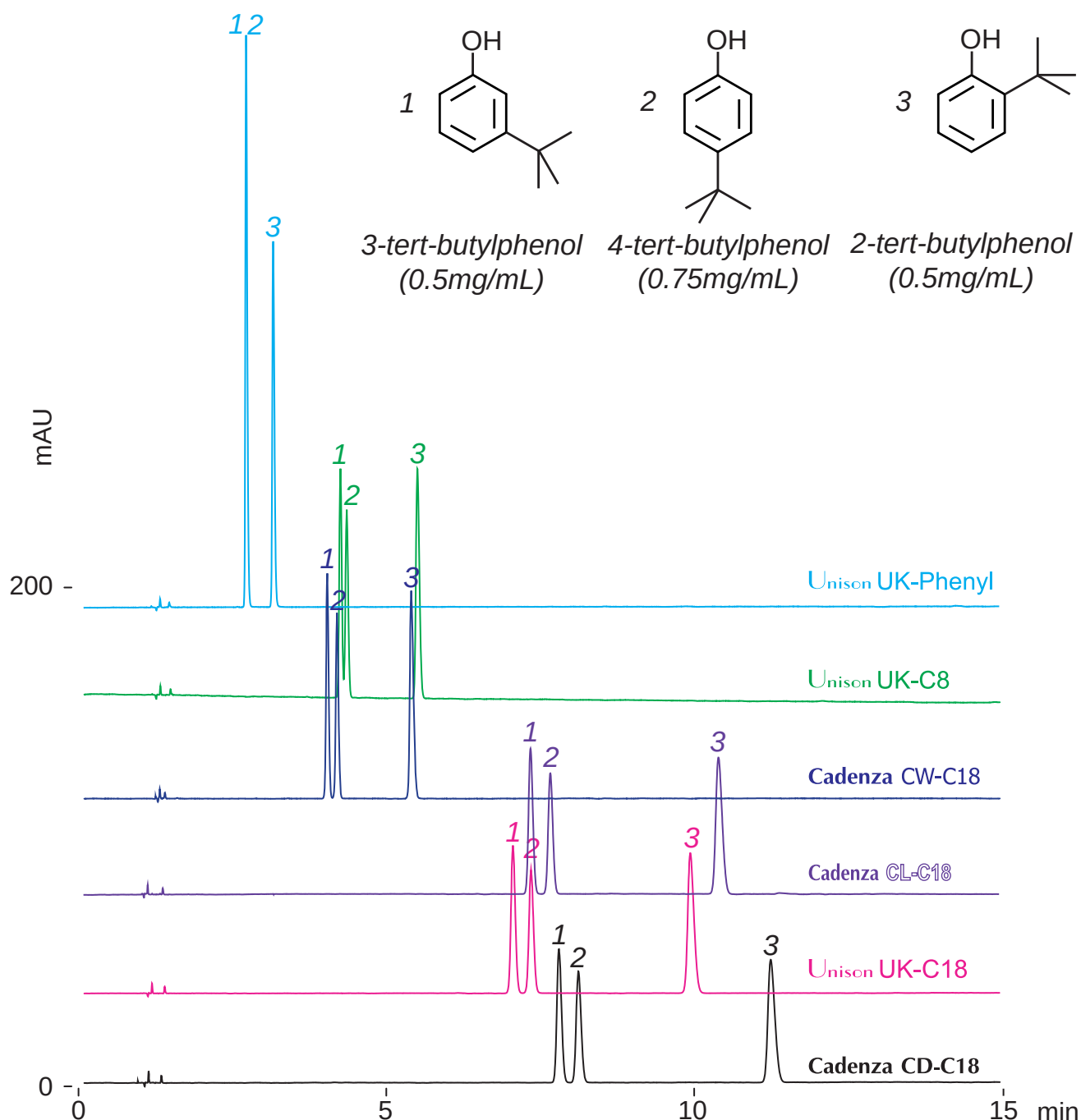

250 x 4.6 mm

A: water / methanol / formic acid = 90 / 10 / 0.1

B: water / methanol / formic acid = 10 / 90 / 0.1

57-62%B (0-40min)

1.0mL/min(20MPa), 50°C, 260nm

10μL (0.5-0.75mg/mL, water / methanol / formic acid = 70 / 30 / 1)