




Scherzo SM-C18
Scherzo SW-C18
Scherzo SS-C18

100 x 3 mm


HPLC COLUMN APPLICATION

LC-MS: Organic disulfides


有機ジスルフィド類の LC-MS 分析

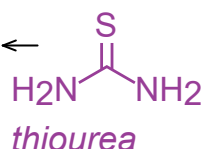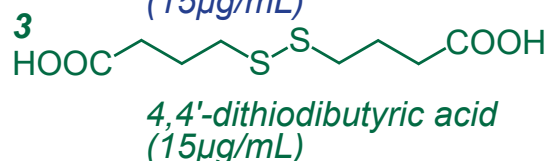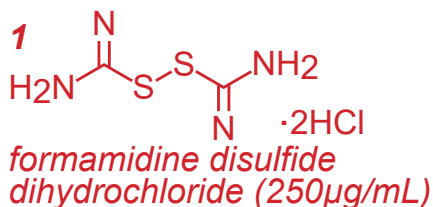


Scherzo SM-C18

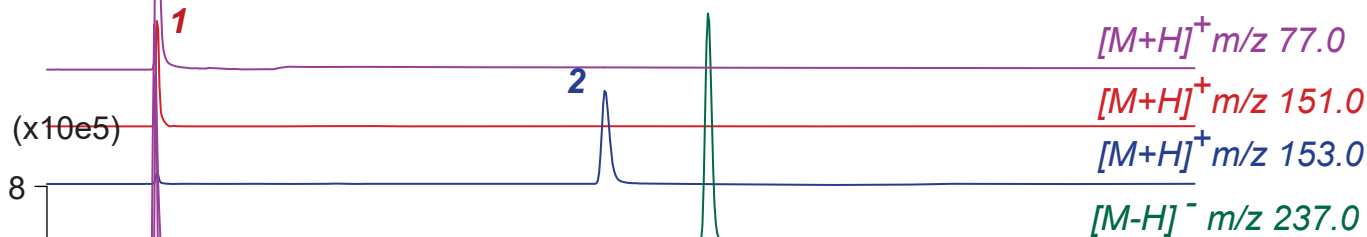

Scherzo SW-C18

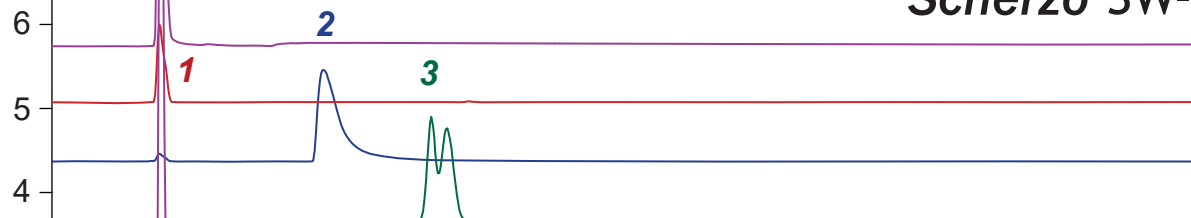

Scherzo SS-C18

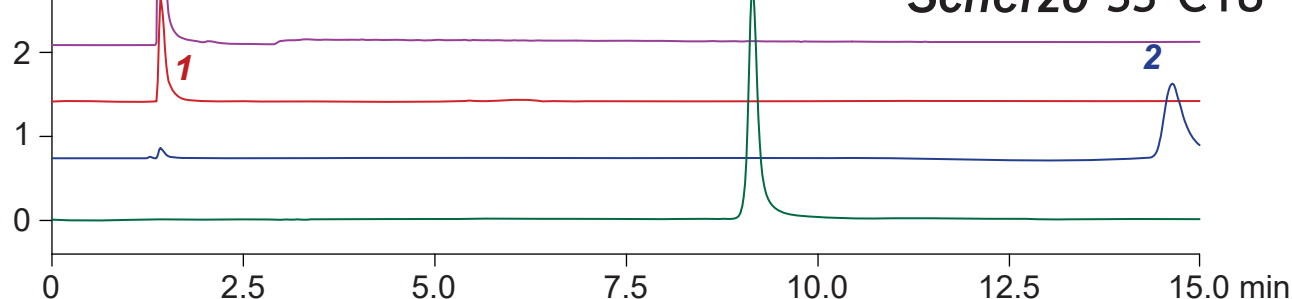

100 x 3 mm

A: 10mM ammonium formate

B: 100mM ammonium acetate / acetonitrile = 80 / 20

0-100%B (0-10min), 100%B (10-15min)

0.4mL/min (7-8MPa), 37°C, Single Quad. MS (ESI)

2µL (15-250µg/mL, MeOH / water = 50 / 50)