

Cadenza CD-C18


Unison UK-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Alizarin, Purpurin (Dye)


アリザリン , プルプリン (染料)

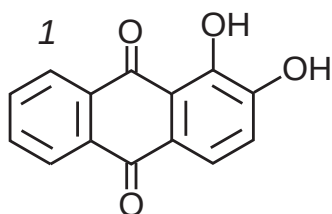

alizarin

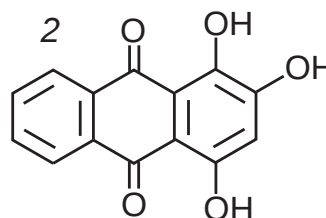

purpurin

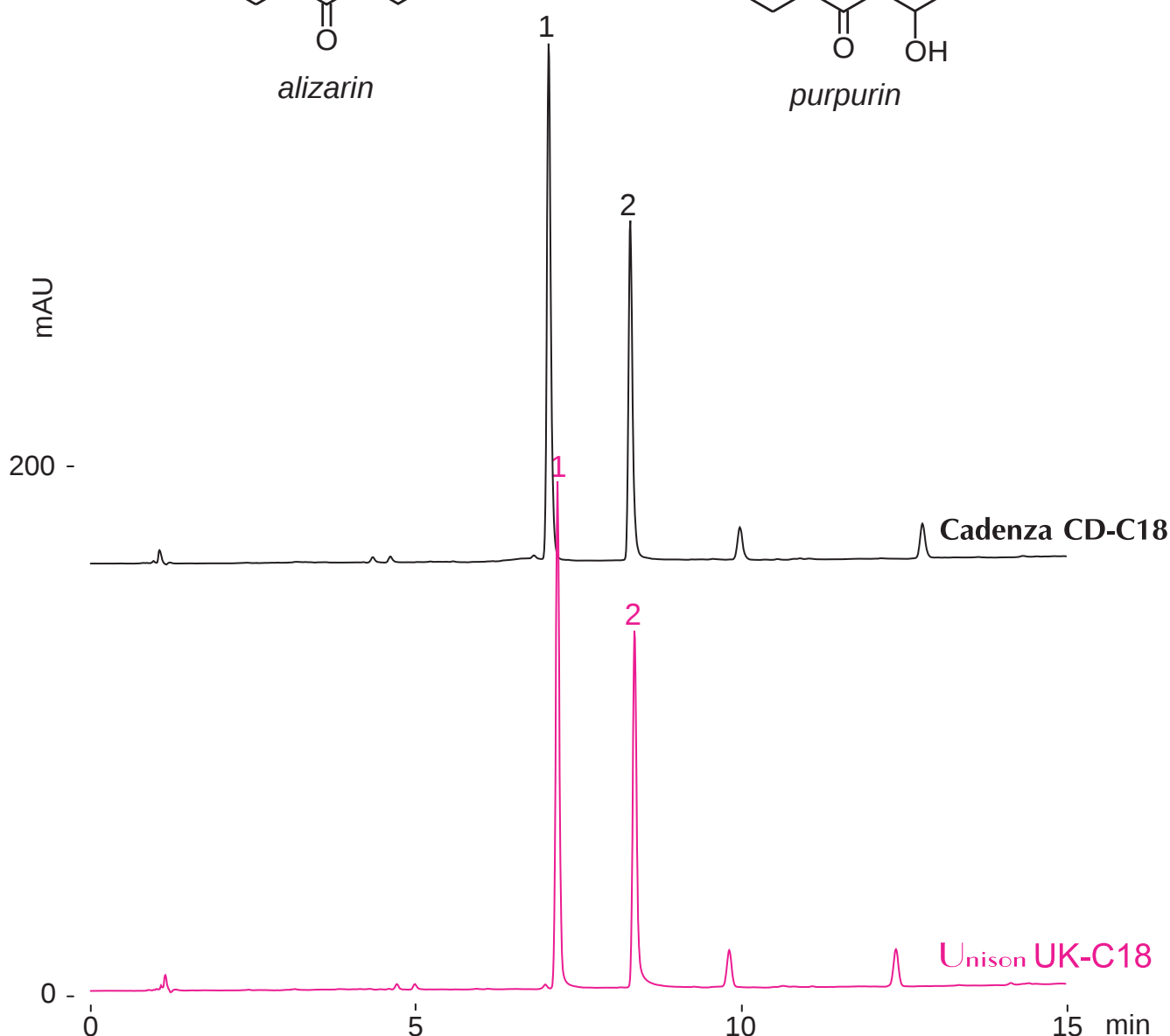

100 x 4.6 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

20-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10μL (0.025mg/mL, water / acetonitrile / formic acid = 33 / 67 / 1)