


Intrada Organic Acid

150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS: Sodium selenite (Therapeutic agent for hyposelenemia)

亜セレン酸ナトリウム (低セレン血症治療薬) の LC-MS 分析

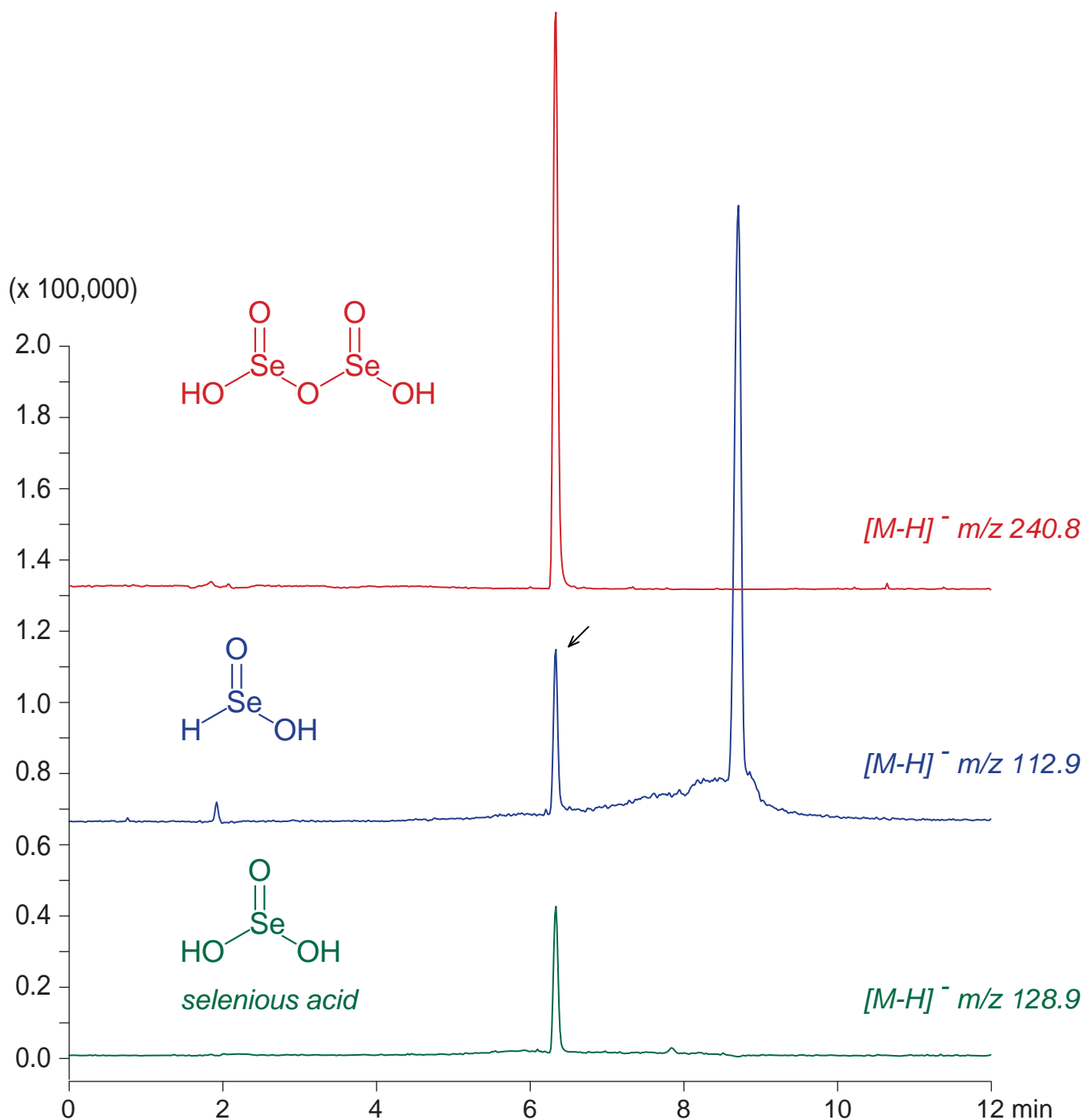

150 x 2 mm

A: acetonitrile / 100mM ammonium formate = 90 / 10

B: acetonitrile / 100mM ammonium formate = 10 / 90

0-100%B (0-7min), 100%B (7-12min)

0.2mL/min (6MPa), 37°C, Single Quad. MS (ESI)

5µL (500µg/mL, acetonitrile / 100mM ammonium formate = 50 / 50)