


Unison UK-Phenyl
Unison UK-C8
Nardis ND-RX
Cadenza CD-C18
Unison UK-C18

150 x 4.6 mm

HPLC COLUMN APPLICATION

Adefovir dipivoxil (Antiviral)

アデホビルジピボキシル (抗ウイルス薬)

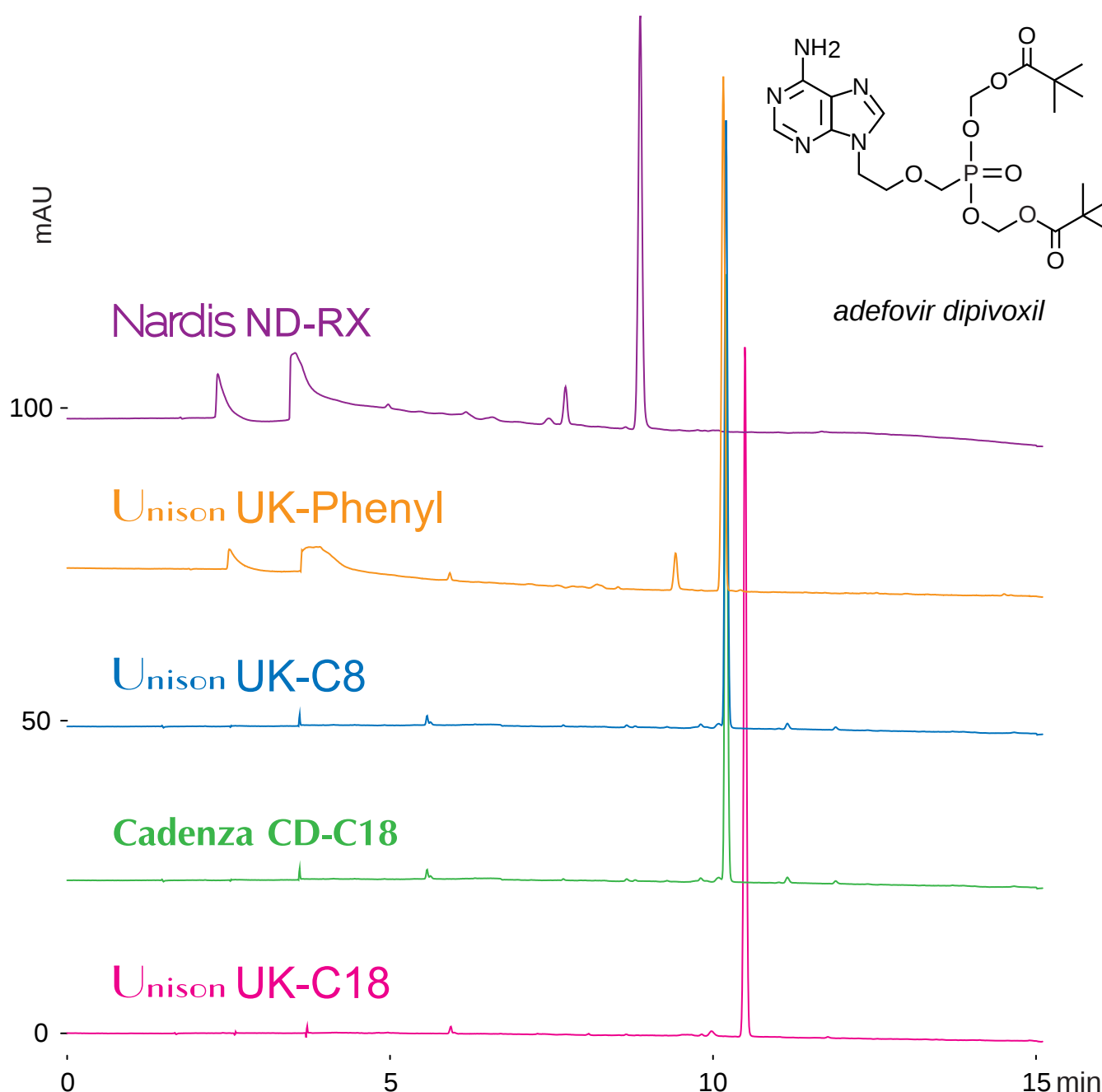

150 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile

0-90%B (0-15min)

1.0mL/min (9-11MPa), 37°C, 260nm

10μL (50μg/mL, water/acetonitrile=50/50)