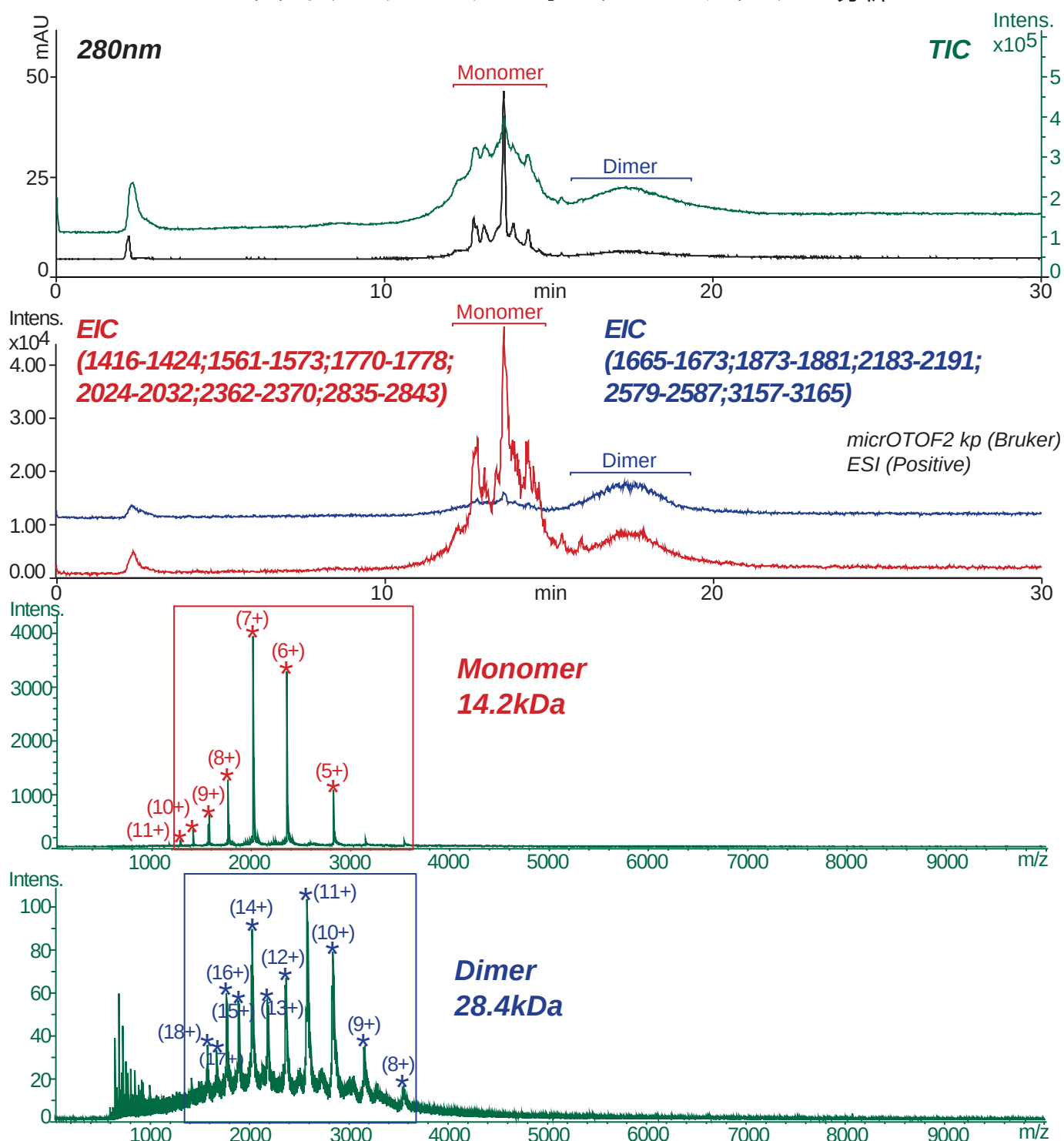


Presto FF-C18

150 x 3 mm

HPLC COLUMN APPLICATION

Non-Porous ODS-MS: Alpha-Lactalbumin

 $\alpha$ -ラクトアルブミンのノンポーラスODSカラム-MS分析


Presto FF-C18, 150 x 3mm

A: 100mM ammonium formate / acetonitrile = 95 / 5

B: 100mM ammonium formate / acetonitrile = 20 / 80

35-60%B (0-30min)

0.2mL/min (26MPa), 37°C

5 $\mu$ L (3mg/mL, 100mM ammonium formate)