


Unison UK-C18

250 x 4.6 mm


Application

Amino Acids - 1


アミノ酸 - 1

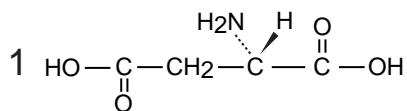


L-aspartic acid (Asp / D)

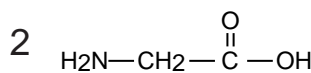


glycine (Gly / G)

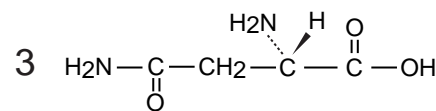


L-asparagine (Asn / N)

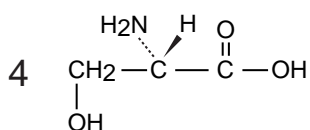


L-serine (Ser / S)

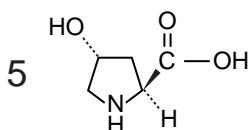

L-hydroxyproline

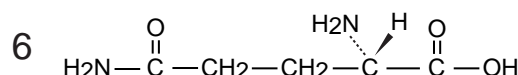

L-glutamine (Gln / Q)

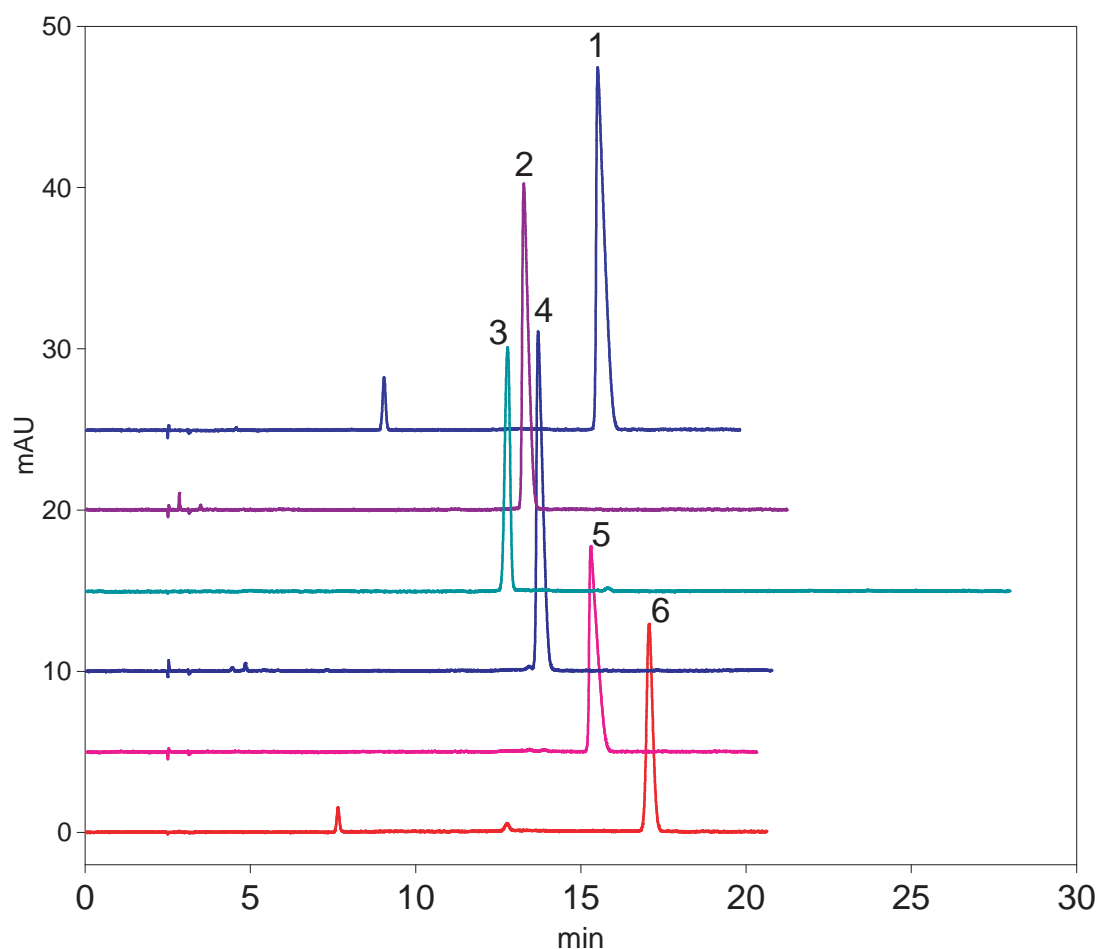

Unison UK-C18, 250 x 4.6 mm  
 50mM H<sub>3</sub>PO<sub>4</sub> + 10mM sodium heptanesulfonate  
 0.8 mL/min, 37 °C, 12.8 MPa, 210 nm, 2-10uL (4-20ug)